

XEMIS: Liquid Xenon Compton Camera for 3γ Imaging

Y. Xing, M. Abaline, S. Acounis, N. Beaupère, J.L. Beney, J. Bert, S. Bouvier, P. Briend, J. Butterworth, T. Carlier, et al.

► To cite this version:

Y. Xing, M. Abaline, S. Acounis, N. Beaupère, J.L. Beney, et al.. XEMIS: Liquid Xenon Compton Camera for 3γ Imaging. 4th International conference on Technology and Instrumentation in Particle Physics, May 2017, Beijing, China. pp.154-158, $10.1007/978-981-13-1316-5_29$. hal-01863097

HAL Id: hal-01863097

https://hal.science/hal-01863097

Submitted on 25 Nov 2021

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

XEMIS: Liquid Xenon Compton Camera for 3γ Imaging

```
Y. Xing<sup>1(⊠)</sup>, M. Abaline<sup>2</sup>, S. Acounis<sup>1</sup>, N. Beaupère<sup>1</sup>, J. L. Beney<sup>1</sup>, J. Bert<sup>3</sup>, S. Bouvier<sup>1</sup>, P. Briend<sup>4</sup>, J. Butterworth<sup>4</sup>, T. Carlier<sup>5</sup>, H. Chanal<sup>2</sup>, M. Cherel<sup>6</sup>, J. P. Cussonneau<sup>1</sup>, M. Dahoumane<sup>7</sup>, L. Gallego-Manzano<sup>1</sup>, D. Giovagnoli<sup>3</sup>, J. Idier<sup>8</sup>, F. Kraeber-Bodéré<sup>5</sup>, P. Le Ray<sup>1</sup>, F. Lefèvre<sup>1</sup>, O. Lemaire<sup>1</sup>, S. Manen<sup>2</sup>, J. Masbou<sup>1</sup>, H. Mathez<sup>7</sup>, E. Morteau<sup>1</sup>, N. Pillet<sup>2</sup>, D. Roy<sup>1</sup>, L. Royer<sup>2</sup>, M. Staempflin<sup>4</sup>, J. S. Stutzmann<sup>1</sup>, R. Vandaele<sup>2</sup>, L. Virone<sup>1</sup>, D. Visvikis<sup>3</sup>, Y. Zhu<sup>1</sup>, and D. Thers<sup>1</sup>
```

¹ SUBATECH, IMT Atlantique, CNRS/IN2P3, Université de Nantes, 44307 Nantes, France Yajing. Xing@subatech. in2p3. fr

² LPC Clermont-Ferrand, 24 Avenue des Landais, Clermont-Ferrand, France
³ INSERM, UMR1101, LaTIM, CHRU Morvan,

2 Avenue Foch, 29600 Brest, France

AIR LIQUIDE Advanced Technologies Division, 38360 Sassenage, France
 Centre Hospitalier Universitaire de Nantes,

1 place Alexis-Ricordeau, 44093 Nantes, France

⁶ INSERM U892 équipe 13, 8 quai Moncousu, 44000 Nantes, France ⁷ CNRS/IN2P3 UMR5822, IPNL Université de Lyon, Lyon, France

⁸ LS2N, Ecole Centrale de Nantes, CNRS/Inp, Université de Nantes, 44307 Nantes. France

Abstract. An innovative liquid xenon Compton camera project, XEMIS (XEnon Medical Imaging System) has been proposed by SUBATECH laboratory, for a new functional medical 3y imaging technique based on the detection in coincidence of 3 γ -rays. The purpose of this 3γ imaging modality is to obtain a 3D image using 100 times less activity than in current PET systems. The combination of a liquid xenon time projection chamber (LXe TPC) and a specific (β^+, γ) radionuclide emitter ⁴⁴Sc is investigated in this concept. In order to provide an experimental demonstration for the use of a LXe Compton camera for 3γ imaging, a succession of R&D programs, XEMIS1 and XEMIS2, have been carried out using innovative technologies. The first prototype XEMIS1 has been successfully validated showing very promising results for energy, spatial and angular resolutions with an ultra-low noise front-end electronics. The second phase dedicated to a 3D imaging of small animals, XEMIS2, is now under installation and qualification, while the characterizations of ionization signal using Monte Carlo simulation has shown preliminary good performances for energy measurement.

Keywords: Medical imaging \cdot 3 γ imaging \cdot Compton camera Liquid xenon

1 Introduction

XEMIS is a XEnon Medical Imaging System directed to the development of personalized medicine, which is based on the 3γ imaging technique and a new liquid xenon Compton camera to realize a low activity medical imaging. This project contains three main phases: the first R&D prototype XEMIS1 has been successfully assessed with good experimental results proving the feasibility of the 3γ imaging concept. The XEMIS2 which is designed for small animal imaging is now under installation and qualification. While for the future research, XEMIS3, a LXe clinical camera for total-body imaging will be investigated.

2 3γ Imaging Technique

The 3γ imaging is a new medical imaging modality proposed by our group, whose principle is to get a 3D location of the radioactive source by detecting simultaneously three gamma-rays [1]. A specific radionuclide that emits quasi-simultaneously a gamma ray and a positron is necessary in this technique. Our best candidate is ⁴⁴Sc, produced by ARRONAX cyclotron while the related radio-pharmaceutical is under researched by Inserm and CHU of Nantes [2].

The position of the radionuclide is obtained by the intersection between the reconstructed line of response (LOR) and Compton cone. Where the LOR is given by the two back-to-back 511 keV gamma-rays resulting from the positron annihilation, and the Compton cone defined from the interaction of the third gamma-ray with a Compton telescope. The opening angle of the cone θ is given by the Compton scattering kinematics, whereas the axis of the cone Δ is determined by the first two interaction positions of the incoming photon inside the detector.

Thanks to the information given by the third gamma ray, the administered activity could be reduced a hundred-time compared to the standard PET. However, in order to get a high-quality image, all interaction vertexes should be detected with a good spatial and energy resolution, which is a big challenge for the experimental feasability. For this purpose, a succession of R&D programs, XEMIS1 and XEMIS2, has been conducted whose recent status will be presented in the following sections.

3 XEMIS1: R&D Prototype

XEMIS1 a small dimension single-phase LXe TPC, designed for testing the possibility of using a LXe Compton telescope for 3γ imaging. XEMIS1 camera consists a cylindrical TPC of 12 cm drift length full of LXe (Fig. 1).

As other liquid noble gases, both scintillation photons and ionization charges are emitted after the interaction of an ionizing particle with LXe. Therefore, at the top of the TPC, it equipped a 1" Hamamatsu photomultiplier tube (PMT) operating in LXe for detecting the VUV light (178 nm). It's used as a trigger of the signal acquisition to get the interaction time t₀. Under the homogeneous drift field generated by 23 copper field rings, the ionized electrons transport to the segmented anode by traversing the Frisch

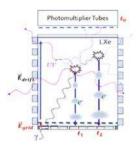




Fig. 1. The schema of the principle of TPC (left) and a view of XEMIS1 TPC (right)

grid that forming the ionization signals represents the deposited energy and position of each interaction vertex. Aims to realize the accurate charges measurement with higher signal/noise ratio, an ultra-low noise (below 100 e⁻ RMS) front-end electronic, IDeF-X operating at LXe temperature have been developed with collaboration.

The performance characterization of XEMIS1 is calibrated by the 511 keV gamma-ray emitted from a low activity ^{22}Na source which is also a $(\beta^+,\,\gamma)$ emitter. The experimental results have shown a 5% energy resolution and a timing resolution of 44.3 \pm 3.0 ns for 511 keV photoelectric events at 1 kV/cm, which is equivalent to a spatial resolution along the z-axis of roughly 100 μm . In addition, a Compton cone measuring test has been done with XEMIS1. The resolution along the LOR is about 8 mm FWHM for 5 cm distant source [3]. This result is limited by the small active area of XEMIS1.

4 XEMIS2: Small Animal Imaging

XEMIS2 is a larger LXe Compton telescope for small animal imaging which is designed for quantifying the potential benefit of the 3γ imaging technique in terms of injected dose. XEMIS2 camera is a monolithic liquid xenon with two back to back cylindrical TPC of 12 cm drift length, where a tube with 10 cm of diameter traverses using to put the little mice. It holds around 200 kg of liquid xenon, totally surrounding the small animal (see Fig. 2). The active volume of detector is covered by 64 Hamamatsu PMTs and two end segmented anodes with a total 20,000 pixels of $3.125 \times 3.125 \text{ mm}^2$, to detect both the VUV scintillation photons and ionization signals produced after interaction of ionizing radiation. A complete Monte Carlo simulation by Geant4 has been done in the first place. The simulation has been realized with a 20 kBq 44 Sc source inside a cylinder full of water (simulating a mouse). Preliminary results shown a quite uniform event detection efficiency in all field of view (FOV). At present, XEMIS2 has been assembled in the laboratory, and prepared to install at the Nantes hospital.

ReStoX Commissioning. As the first LXe camera will be installed in a hospital, it should be compact and safe. For this purpose, a special cryogenic system named ReStoX (Recovery and Storage system of Xenon) has been built with the Air Liquid



Fig. 2. XEMIS2 camera with ReStoX

and SUBATECH collaboration. It could store 210 kg xenon from room temperature to liquid xenon temperature (\sim – 100 °C). In addition, the xenon cooling with liquid nitrogen is very powerful, up to 10 kW. The commissioning of ReStoX has been done in our laboratory. Results have shown that, for 200 kg LXe, without cooling source, the room temperature would be reached from LXe temperature in almost 1 year. And in case of extreme situation, a passive recovery with the gravity could be realized in almost 10 min.

Ionization DAQ Chain. For the low activity imaging, all the interaction signals need to be recorded without dead time during 20 min, that means $\sim 10^4$ charges and time signals/pixel/s. In XEMIS2, there are more than 20000 electronic channels which is a big challenge for the data storage. The IDeF-X [4] and XTRACT [5] have been developed for the ionization DAQ, permitting reading out only the amplitude, time and address of pixel for each ionization waveform in order to control the data storage volume. Finally, there are only 80 process units connected to the 6 FPGA board to reduce the traverse number. There are still lots works to do before the operation of XEMIS2.

5 Conclusion

 3γ imaging would be a new medical modality thanks to LXe technology. The design and construction of XEMIS camera are quite challenging that involves lots of new technologies and innovations. At present, the second prototype XEMIS2 for small animal imaging is under qualification and it will be the first LXe camera installed in a hospital. The expected image quality is very promising with a quite uniform event detecting efficiency in all FOV. In addition, the LXe technology is scalable so that a design of large whole-body camera could be investigated in the future.

References

- 1. Grignon, C., et al.: Nuclear medical imaging using β + γ coincidences from Sc-44 radio-nuclide with liquid xenon as detector medium. Nucl. Instrum. Methods A **571**(1), 142–145 (2007)
- Huclier-Markai, S., et al.: Optimization of reaction conditions for the radiolabeling of DOTA and DOTA-peptide with 44 m/44 Sc and experimental evidence of the feasibility of an in vivo PET generator. Nucl. Med. Biol. 41, e36–e43 (2014)
- Gallego Manzano, L., et al.: XEMIS: a liquid xenon detector for medical imaging. Nucl. Instrum. Methods 787, 89–93 (2015)
- 4. Gevin, O., et al.: IDeF-X V1. 0: A new 16-channel low-noise analog front-end for Cd (Zn) Te detectors. Nucl. Instrum. Methods **567**(1), 140–144 (2006)
- Gallego Manzano, L.: Optimization of a single-phase liquid xenon compton camera for 3γ medical imaging. PhD diss., Nantes, Ecole des Mines (2016)