

Development of a Mass Spectrometry Imaging (MSI) approach to explore Listeria monocytogenes biofilms exposed to a dehumidification stress

Tiago Santos, Laetitia Theron, Christophe Chambon, Didier Viala, Delphine Centeno, Michel Hébraud

▶ To cite this version:

Tiago Santos, Laetitia Theron, Christophe Chambon, Didier Viala, Delphine Centeno, et al.. Development of a Mass Spectrometry Imaging (MSI) approach to explore Listeria monocytogenes biofilms exposed to a dehumidification stress. 8. Colloque du Réseau National Biofilms, Dec 2017, Clermont-Ferrand, France. 1 p., 2017. hal-01901461

HAL Id: hal-01901461

https://hal.science/hal-01901461

Submitted on 22 Oct 2018

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Development of a Mass Spectrometry Imaging (MSI) approach to explore *Listeria monocytogenes* biofilms exposed to a dehumidification stress



<u>Tiago Santos</u>¹, Laëtitia Théron², Christophe Chambon², Didier Viala², Delphine Centeno² and Michel Hébraud^{1,2}







¹ Université Clermont Auvergne, INRA, UMR MEDiS, F-63122 Saint-Genès Champanelle, France ² INRA, Plateforme d'Exploration du Métabolisme (PFEM), F-63122 Saint-Genès Champanelle, France

Introduction

- Human listeriosis cases are due to the ingestion of foods contaminated with *Listeria monocytogenes*. Control of this pathogen is a daily concern for many food industries and represents a major challenge in terms of hygiene and public health. *L. monocytogenes* is capable of forming biofilms and surviving harsh environmental conditions. It can persist in microniches despite the daily procedures of cleaning disinfection. These procedures also cause significant variations in the relative air humidity and the dehumidification of the premises is an important step for the control of microbial contaminations. Despite this, part of the population is able to resist these treatments. The molecular mechanisms by which biofilms adapt to dehumidification, often applied empirically, are not well known and it is important to develop tools to better decipher them.
- Matrix-assisted laser desorption/ionization time-of-flight imaging mass spectrometry (MALDI-TOF IMS) is a surface-sampling technology that can determine spatial information and relative abundance of analytes directly from biological samples. Basically applied to animal or human tissues in biomedical research, applications have diversified in other scientific fields. By assimilating a microbial biofilm to a tissue, this approach has been implemented to investigate its potential to provide relevant data on molecular mechanisms that allow cells to adapt to air dehumidification stress (moderate desiccation).
- 🔆 This study aims to develop an IMS approach and an in situ identification method to explore the protein expression within L. monocytogenes biofilms exposed to air dehumidification.

Workflow


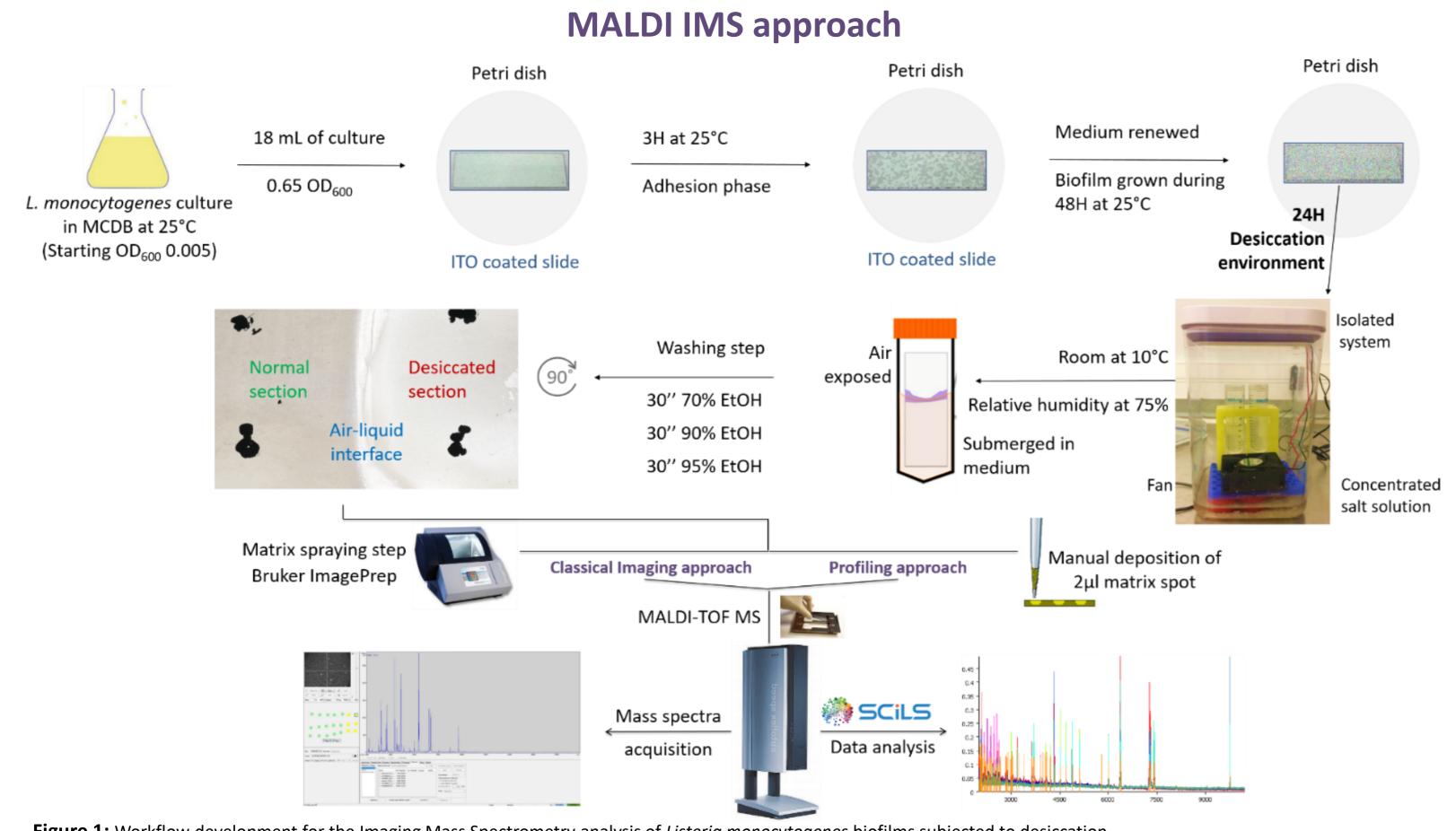


Figure 1: Workflow development for the Imaging Mass Spectrometry analysis of Listeria monocytogenes biofilms subjected to desiccation

Results

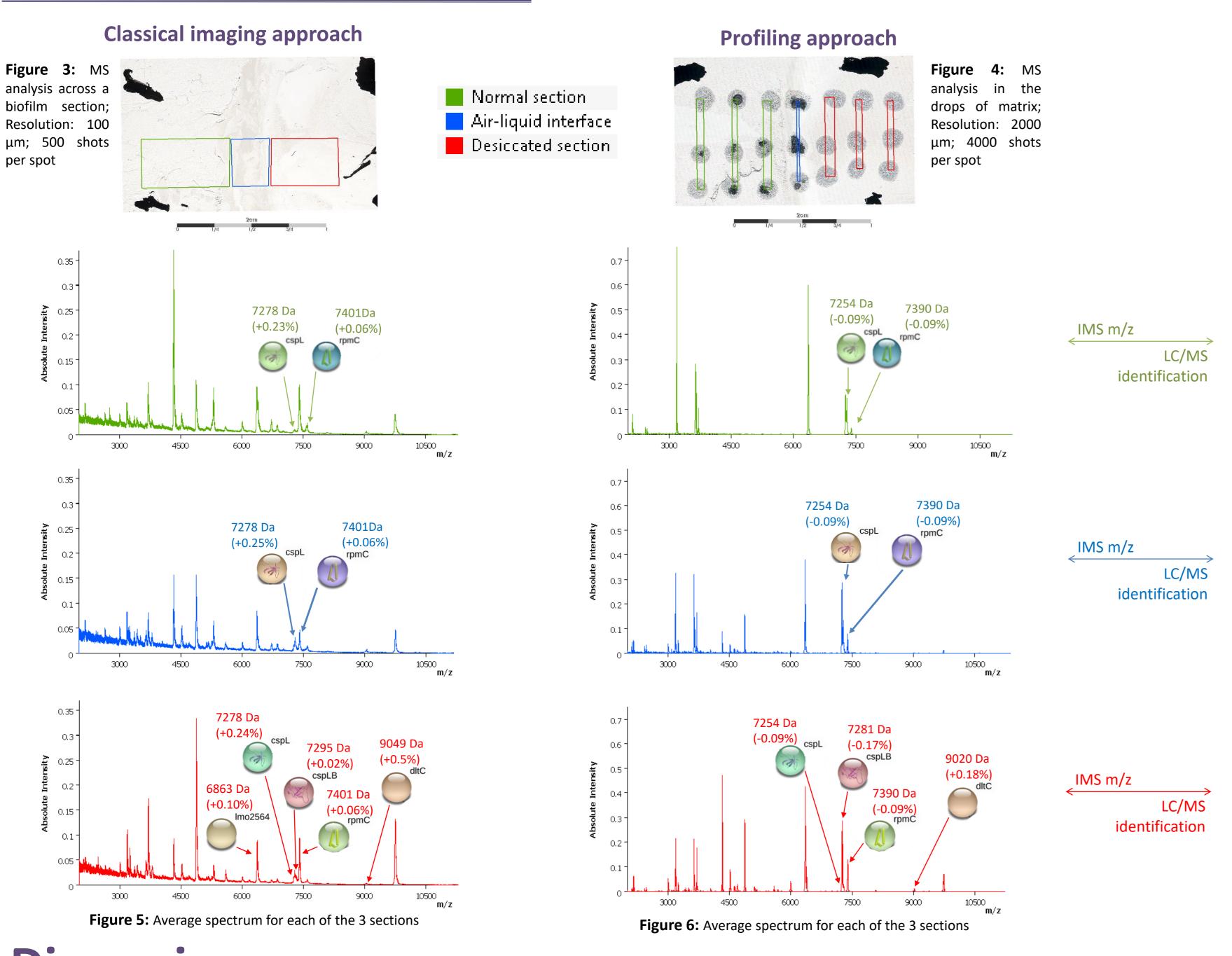

Discussion

- The 3 sections of analysis presented good spectra in both approaches, classical imaging and profiling. The profiling showed lower level of background noise and higher intensity of peaks than the classical approach. This is due to the higher ratio matrix-sample present in the profiling and also the superior signal accumulation obtained in one spot of analysis.
- The *in situ* extraction of proteins resulted in a high and reproducible number of proteins identified. The higher amount of proteins identified in the interface and desiccated section could be related to their exposition to the stress condition, and therefore an increase of cell lysis and proteins available for extraction.
- The back correlation between the significative m/z intervals and the *in situ* identification lead to 5 different proteins being assigned to the mass spectra (mass tolerance of 0.5%): In the normal section and interface, a cold shock protein (CspLA) and a 50S ribosomal protein L29; From the desiccated section, a second CspL (CspLB) which is involved in adaptation to atypical conditions, a tautomerase Imo2564 and the membrane-associated protein D-alanine-poly(phosphoribitol) ligase, which is involved in the pathway lipoteichoic acid biosynthesis. Loss of D-alanylation of lipoteichoic acids alters the cell surface charge and results in reduced biotic attachment and biofilm production.

In situ extraction of proteins and identification by nanoLC-MS/MS


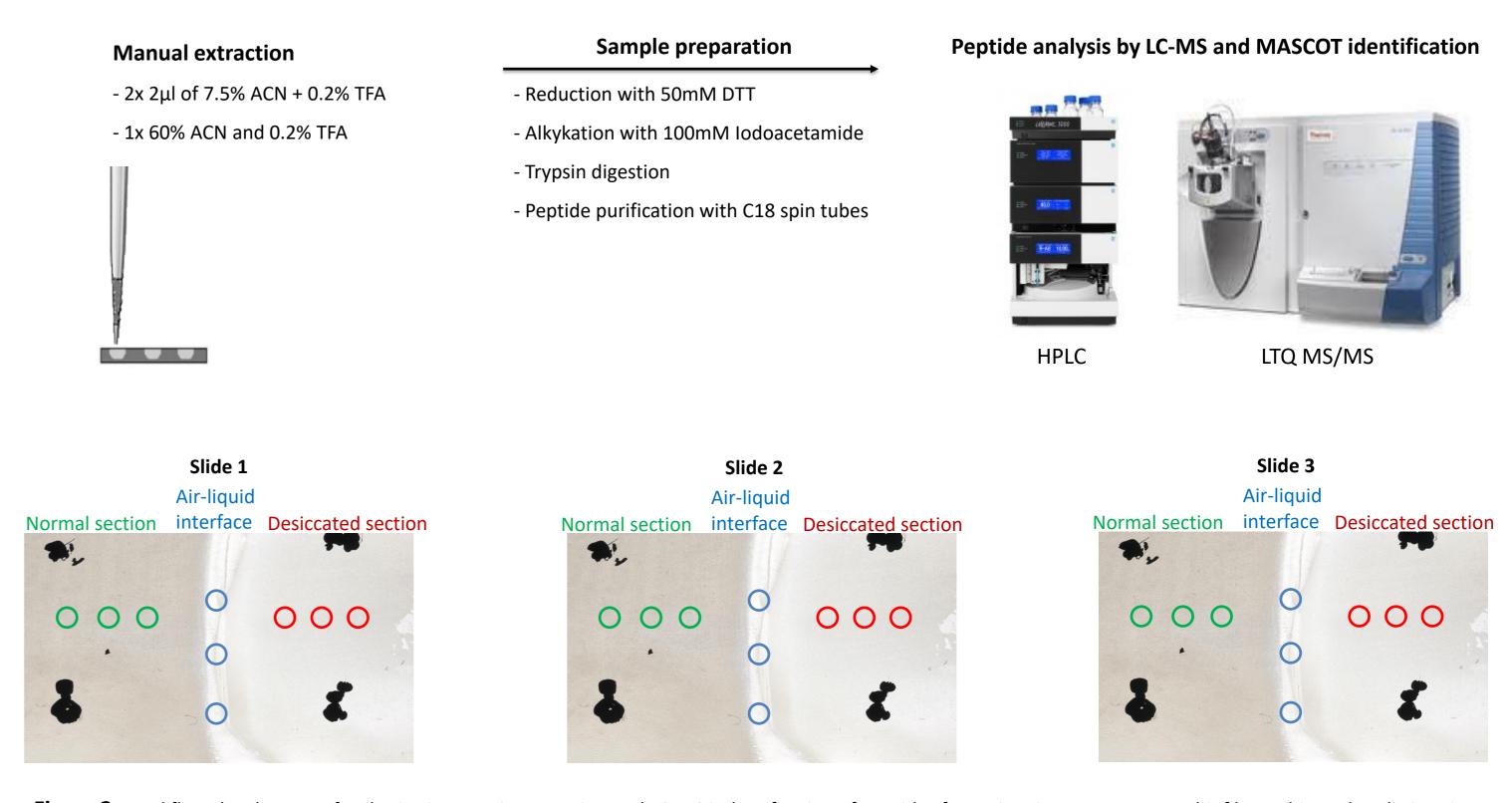


Figure 2: Workflow development for the in situ protein extraction and LC-MS indentification of peptides from Listeria monocytogenes biofilms subjected to desiccation

Desiccated 2

Interface 3

Interface 2

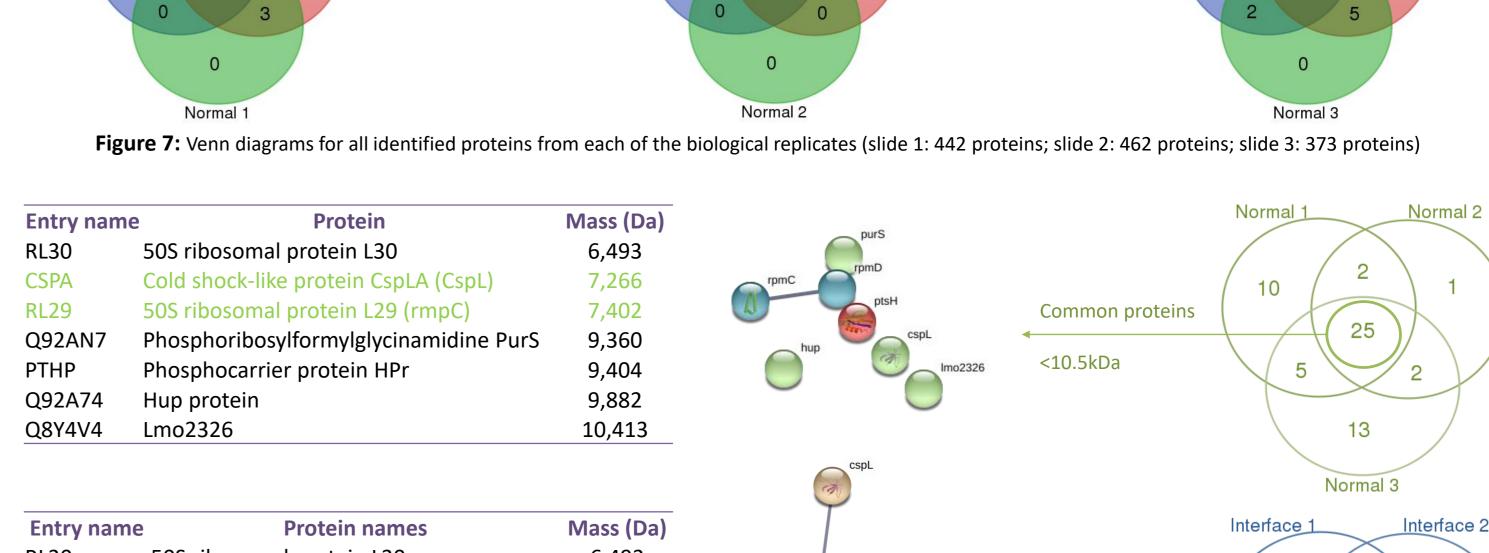

Entry name	Protein names	Mass (Da)
RL30	50S ribosomal protein L30	6,493
CSPA	Cold shock-like protein CspLA (CspL)	7,266
RL29	50S ribosomal protein L29	7,402
Q92AN7	Phosphoribosylformylglycinamidine PurS	9,360
PTHP	Phosphocarrier protein HPr	9,404
Q92A74	Hup protein	9,882
RS17	30S ribosomal protein S17	10,036
CH10	10 kDa chaperonin (GroES protein)	10,064
Entry name	e Protein names	Mass (Da)
RL30	50S ribosomal protein L30	6,493
RS21	30S ribosomal protein S21	6,846
Y2564	Probable tautomerase	6,861
CSPA	Cold shock-like protein CspLA (CspL)	7,266
CSPB	Cold shock-like protein CspLB (CspB)	7,298
RL29	50S ribosomal protein L29	7,402
Y1028	UPF0356 protein lmo1028	8,292
Q8Y3X2	Lmo2707 protein	8,337
Q92D14	Lmo1008 protein	8,956
ムノムレエ す	LINGTOOD DIOLEIN	0.550

RS19 30S ribosomal protein S19 10,475 **Table 1:** Description of the proteins present in Figure 9; Highlighted proteins correspond to the ones correlated with m/z intervals from IMS

D-alanine--poly(phosphoribitol)

Phosphoribosylformylglycinamidine PurS

30S ribosomal protein S20

Phosphocarrier protein HPr

UPF0237 protein Imo0533

30S ribosomal protein S17

10 kDa chaperonin (GroES protein)

Hup protein

RS20

Q92AN

Q92A74

PTHP

Y533

RS17

CH10

Q8Y4V4

Conclusions

and profiling (mass tolerance: 0.5%)

Lmo2326 protein

Interface 1

Desiccated 1

These MALDI IMS and innovative protein identification approaches contributes to dissect the spatial proteome of an intact bacterial biofilm giving a new

identified, with mass below 10.5kDa,

from each of the three section

9,169

9,360

9,882

9,943

10,036

10,064

10,413

spatial proteome of an intact bacterial biofilm giving a new insight into protein expression relating to biofilm adaptation.

