



# Hybrid Electron Microscopy Normal Mode Analysis with Scipion

Mohamad Harastani, Carlos Oscar S. Sorzano, Slavica Jonic

## ► To cite this version:

Mohamad Harastani, Carlos Oscar S. Sorzano, Slavica Jonic. Hybrid Electron Microscopy Normal Mode Analysis with Scipion. *Protein Science*, 2020, 29 (1), pp.223-236. 10.1002/pro.3772 . hal-02388487

**HAL Id: hal-02388487**

**<https://hal.science/hal-02388487>**

Submitted on 1 Dec 2019

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

# **Hybrid Electron Microscopy Normal Mode Analysis (HEMNMA) with Scipion**

Mohamad Harastani <sup>1</sup>, Carlos Oscar S. Sorzano <sup>2</sup>, Slavica Jonić <sup>1\*</sup>

<sup>1</sup> Sorbonne Université, UMR CNRS 7590, Muséum National d'Histoire Naturelle, IRD,  
Institut de Minéralogie, de Physique des Matériaux et de Cosmochimie, IMPMC, 75005 Paris,  
France

<sup>2</sup> Biocomputing Unit, Centro Nacional de Biotecnología – CSIC, 28049 Madrid, Spain.

\* Corresponding author

## **Contact details of the corresponding author:**

Dr. Slavica Jonić  
Sorbonne Université  
CC 115  
4 Place Jussieu  
75005 Paris, France  
Tel: +33 1 44 27 72 05  
Fax: +33 1 44 27 37 85  
E-mail: [slavica.jonic@upmc.fr](mailto:slavica.jonic@upmc.fr)

**Running title:** HEMNMA with Scipion

**Total number of manuscript pages:** 29

**Total number of figures:** 8

**Abstract:** Hybrid Electron Microscopy Normal Mode Analysis (HEMNMA) method was introduced in 2014. HEMNMA computes normal modes of a reference model (an atomic structure or an electron microscopy map) of a molecular complex and uses this model and its normal modes to analyze single-particle images of the complex to obtain information on its continuous conformational changes, by determining the full distribution of conformational variability from the images. An advantage of HEMNMA is a simultaneous determination of all parameters of each image (particle conformation, orientation, and shift) through their iterative optimization, which allows applications of HEMNMA even when the effects of conformational changes dominate those of orientational changes. HEMNMA was first implemented in Xmipp and was using MATLAB for statistical analysis of obtained conformational distributions and for fitting of underlying trajectories of conformational changes. A HEMNMA implementation independent of MATLAB is now available as part of a plugin of Scipion V2.0 (<http://scipion.i2pc.es>). This plugin, named ContinuousFlex, can be installed by following the instructions at <https://pypi.org/project/scipion-em-continuousflex>. In this article, we present this new HEMNMA software, which is user-friendly, totally free, and open-source.

**Keywords:** cryo-electron microscopy; single particle analysis; continuous conformational changes; normal mode analysis; structure; dynamics; software

**50-75 word statement for a broader audience:** This article presents Hybrid Electron Microscopy Normal Mode Analysis (HEMNMA) software that allows analyzing single-particle images of a complex to obtain information on continuous conformational changes of the complex, by determining the full distribution of conformational variability from the images. The HEMNMA software is user-friendly, totally free, open-source, and available as

part of ContinuousFlex plugin (<https://pypi.org/project/scipion-em-continuousflex>) of Scipion V2.0 (<http://scipion.i2pc.es>).

## INTRODUCTION

Cryo-electron microscopy (cryo-EM) has become comparable to X-ray crystallography with regards to the obtainable resolution of structures of biomolecular complexes, which are now increasingly determined at near-atomic resolution.<sup>1-10</sup> No requirement for sample crystallization and the possibility to elucidate multiple conformations of a complex from the same sample are among the main advantages of cryo-EM. Characterizing the different conformations that can coexist is essential for understanding how the complexes function and addressing their dynamics.

To achieve near-atomic resolution of 3D reconstructions, the classical approach is to collect a large number of images of complexes (particles) at random and unknown orientations within a thin layer of vitreous ice, then, perform 2D and 3D classifications into an initially set number of classes and, finally, perform 3D reconstruction using only those particles that have the most consistent views and conformations (those that contribute to the highest-resolution class averages) while removing all other particles.<sup>8-14</sup> Such “selection” of particles may obscure information on a possibly larger conformational variability as some conformational states may be thrown away blindly instead of being elucidated. Therefore, non-classification-based methods are required to extract from images the full distribution of conformational variability (the so-called conformational space or landscape in which the images are mapped) and to assemble 3D reconstructions from identified, more or less dense regions in this space (denser regions contain more frequent conformational states and less dense regions contain less frequent states). Such methods are necessary for studying continuous conformational changes of complexes (the concept including the possibility of unequally distributed conformers, e.g. possible existence of more stable conformers) and are referred to as continuous-state methods. Many classification-based methods (referred to as discrete-state methods) can be found in the

literature.<sup>15-28</sup> While the majority of the discrete-state methods require that the number of classes is set initially, setting this number in some of the methods is less arbitrary as it is based on making a balance between intraclass and interclass variances after a statistical analysis of the data, typically using an eigenvector analysis of the covariance estimated from images or from volumes reconstructed from random subsets of images.<sup>18; 19; 21</sup> The covariance eigenvectors estimation has also been introduced in the context of continuous conformational heterogeneity.<sup>29; 30</sup> The latter methods are closely related to continuous-state methods whose development started recently and is currently an active field of research.<sup>31-37</sup>

We introduced the continuous-state method referred to as Hybrid Electron Microscopy Normal Mode Analysis (HEMNMA) in 2014.<sup>31; 32</sup> HEMNMA computes normal modes of a reference model (an atomic structure or an EM map) of a molecular complex and uses this model and its normal modes to analyze single-particle images of the complex iteratively in order to extract information on the full distribution of conformational variability from these images. In each iteration, HEMNMA aims at simultaneously solving both orientational and conformational heterogeneity of images using the reference model and its normal modes.

Although often considered as “model-free”, the continuous-state methods from other groups do use a model (e.g., the model used by Dashti et al (2014)<sup>33</sup> is EMDB-1067 density map). This model is used to first solve the orientational heterogeneity of images, by determining the particle 3D orientation and 2D shift in each image assuming conformational homogeneity of images (e.g., the particle orientation and shift in each image can be determined by the standard projection matching of the image with the model). Then, the conformational heterogeneity of images is solved (i.e., the conformations are determined) using the image orientations and shifts determined in the previous step. These methods do not refine

orientations and shifts when determining conformations, assuming that the effects of orientational changes dominate those of conformational changes so that each image can be associated with the correct orientation.

An advantage of HEMNMA is a simultaneous determination of all parameters of each image (particle conformation, orientation, and shift) through their iterative optimization, which allows applications of HEMNMA even when the effects of conformational changes dominate those of orientational changes. The simultaneous determination of the parameters is done by iterative matching of the image with projections of the reference model being deformed using a linear combination of normal modes (the conformation is estimated by determining the unknown coefficients of the linear combination, i.e. by determining the displacement amplitudes along normal modes). The reference model does not need to be an atomic-resolution model and can also be an EM map. This EM map can be obtained by 3D reconstruction from combined data as if the data was conformationally homogeneous (e.g., by the so-called random sample consensus approach that identifies the best EM maps from those reconstructed using random subsets of images and random image orientations,<sup>38</sup> maximum likelihood optimization of a number of 3D reconstructions starting from random data subsets,<sup>39</sup> or other methods<sup>40</sup>).

HEMNMA and other continuous-state methods are fundamentally different from classification-based (discrete-state) methods. For instance, the discrete-state method of Haselbach et al (2018)<sup>27</sup> uses standard Relion<sup>39</sup> software for discrete classification of images into a number of 3D maps and standard principal component analysis to cluster the maps. As similar classification-based methods, it is prone to assigning images to lower-noise density maps independently of whether or not the images belong there. For further reading on the

different methods, the reader can refer to some of the recent methods reviews.<sup>41-45</sup>

The first version of HEMNMA was implemented in Xmipp<sup>46</sup> and was using MATLAB for statistical analysis of obtained conformational distributions and for fitting of underlying trajectories of conformational changes.<sup>31; 32</sup> A MATLAB-independent version of HEMNMA was recently implemented in Scipion<sup>47</sup> and is currently available as part of a plugin of Scipion V2.0 (<http://scipion.i2pc.es>). This plugin, named ContinuousFlex, can be installed by following the instructions at <https://pypi.org/project/scipion-em-continuousflex>. The new HEMNMA software is totally free and open-source, and the new graphical interface is even more user-friendly than the previous one. In this article, we present this new HEMNMA software.

## **HEMNMA WITHIN SCIPION: STEPS, PARAMETERS, GRAPHICAL INTERFACE**

The steps to follow are listed in the HEMNMA menu on the left side of the Scipion project window and are numbered from 1 to 6. When the different steps are executed, a tree-like structure of the project appears on the right side of the window, with each block corresponding to a step or a substep (**Figs. 1-2**). In this section, we describe the steps (Steps 1-6), parameters, and graphical interface of HEMNMA. Each parameter is also described in the graphical interface, via a help message that can be displayed by clicking on the question mark next to the corresponding parameter field. The majority of the parameters are set to the values that usually produce good results. These default values are visible in the graphical interface and can be modified by the user. The parameters whose values are expected to be less frequently changed are hidden by default. By selecting “Advanced” as “Expert Level” mode (“Normal” is the default mode), the values of these “advanced” parameters can be visualized (the parameter names showing up on a grey background) and modified (e.g., **Fig 3A**).



**Step 1. Reference model:** The reference model to import can be a PDB file with atomic coordinates or a density volume (e.g., an EM map or a simulated map from a PDB structure). Both import options are provided (steps 1.a and 1.b1 in the project tree, **Figs. 1-2**). **Fig. 1** and **Fig. 2** show the project tree for an input atomic structure and for an input density volume, respectively. Before the next step, the input density volume must be converted into a PDB-format file. We provide a tool for the volume conversion into a set of 3D Gaussian functions (step 1.b2 in the project tree, **Fig. 2**), the so-called pseudoatoms,<sup>48</sup> whose coordinates are written into a PDB file. More precisely, the conversion is performed so that the input density volume is represented with Gaussian functions of a given (input) standard deviation (“Pseudoatom radius” input parameter expressed in voxels), by minimizing the normalized volume approximation error<sup>48</sup> to reach a given (input) value of this error (“Volume approximation error” parameter expressed in percent) (**Fig. 3**). Note that the default value of the “Volume approximation error” parameter is 5 %, which can be modified by first selecting “Advanced” as “Expert Level” mode to visualize the parameter (**Fig. 3A**). Through this conversion procedure, some noise may be removed from the input density volume, as shown elsewhere.<sup>49</sup> However, strong particle background noise should be removed using external tools or the masking procedures available in HEMNMA (“Mask mode” parameter in **Fig. 3A** allows entering the name of a binary mask file or the threshold value below which the densities should be removed). For faster processing, the volume-to-pseudoatoms conversion was parallelized and should be run using several CPU threads (sharing the same memory). The number of threads can be specified via the “Threads” parameter (**Fig. 3A**). The results of the volume-to-pseudoatoms conversion can be visualized with Chimera<sup>50</sup> (**Fig. 3B**) by clicking on “Analyze Results” button in the project window (**Fig. 2**).

After a displacement with normal modes, both an atomic reference model and a pseudoatomic reference model (from a reference density volume) are converted into a density volume whose projections are compared with images to determine the conformations in these images (Step 5). The higher the resolution of the reference model, the more accurate the normal modes and the conformational determination will be. Thus, the most accurate conformational determination is expected with input atomic models and high-resolution EM maps.

The reference atomic model can be a full-atomic model or a coarse-grained model (e.g., containing C $\alpha$  atoms only). The computation of normal modes (Step 2) requires information on the shape of the entire complex. An atomic model lacking a few small regions (e.g., disordered regions) may still be used as the reference model. If the atomic model lacks large portions of a complex and an EM map can be obtained for the entire complex, this EM map should be used as the reference model. An alternative reference model could be a hybrid model of the entire complex obtained by modeling the missing parts in the atomic model or the density volume computed from this hybrid model.

Classical, classification-based approaches can be used to obtain a reference EM map for HEMNMA to analyze the entire heterogeneous set of images or an image subset with heterogeneity reduced through classification, as shown elsewhere<sup>31</sup>. The use of classification-based approaches before HEMNMA may be particularly useful in the case of combined discrete and continuous conformational heterogeneity (e.g., a mixture of flexible complexes bound and unbound with ligands).

There are no restrictions regarding the number of atoms in input PDB files or the size of input EM maps, except that the maps should have a cubic shape. There are no restrictions regarding

the content of the atomic model, such as the presence of DNA, ligands, etc. However, in the case of a heterogeneous set of particles bound and unbound with ligands, it may be interesting to remove the ligand from the reference model, as shown elsewhere.<sup>31; 51</sup>

**Step 2. Normal mode analysis (NMA):** Normal modes of the atomic or pseudoatomic representation of the reference model (“Input structure” field in **Fig. 4A**) are computed using Tirion’s elastic network model<sup>52</sup>. In this model, atoms or pseudoatoms interact if they are connected with elastic springs. The interaction cut-off distance parameter determines how many atoms or pseudoatoms will be connected with elastic springs. The following two options exist for setting this value, via “Cut-off mode” parameter (**Fig. 4A**): “absolute” and “relative”. If “absolute” is chosen for “Cut-off mode” (unshown in **Fig. 4A**), the interaction cut-off distance value is an input parameter that should be set directly by the user (by setting “Cut-off distance” parameter, expressed in angstroms, which shows up only if “absolute” is chosen for “Cut-off mode”). The default value of “Cut-off distance” is 8 Å, which is the empirical recommended value for atomic structures. The “Cut-off distance” value of 8 Å implies building an elastic network in which two atoms are connected with a spring (interact with each other) only if their distance is smaller than 8 Å. If “relative” is chosen for “Cut-off mode”, the interaction cut-off distance value is computed automatically based on a given (input) percentage of distances (“Cut-off percentage” parameter expressed in percent, **Fig. 4A**) that should be below this value. The default value of “Cut-off percentage” is 95 %, which is the empirical recommended value for pseudoatomic structures. The “Cut-off percentage” value of 95 % implies building an elastic network in which two pseudoatoms are connected with a spring only if their distance is smaller than the one below which 95 % of distances are. The larger the interaction cut-off distance for a given complex, the more rigid the elastic network will be. A too large interaction cut-off distance produces a too rigid elastic network,

which can be detected and corrected by analyzing the normal-mode animations with the provided graphical interface (NMA results viewer in **Fig. 4B**).

Normal modes are computed by diagonalizing a square matrix of second derivatives of the potential energy function of the structure (Hessian matrix). One dimension of the Hessian matrix is 3 times the number of atoms or pseudoatoms. For faster Hessian diagonalization with input atomic structures (large Hessians due to large numbers of atoms), the structure is split into blocks of several residues, each block having 6 degrees of freedom (3 rotations and 3 translations), which reduces the Hessian dimension to 6 times the number of blocks and is known as Rotation Translation Block (RTB) method.<sup>53</sup> The number of residues per RTB block should be set by the user, and its default value in the graphical interface is 10 (**Fig. 4A**). The larger the number of residues per block (i.e., the smaller the number of blocks), the faster the computation of normal modes and the more rigid the elastic network will be. A too large number of residues per block will produce a too rigid elastic network, which can be detected and corrected by analyzing the normal-mode animations with the provided graphical interface (NMA results viewer, **Fig. 4B**).

The number of computed normal modes is 3 times the number of pseudoatoms (for a pseudoatomic structure) or 6 times the number of RTB blocks (for an atomic structure). However, the number of modes written on the disk and animated should be set by the user (“Number of modes” parameter, **Fig. 4A**). It is usually enough to save animations for 20-100 modes (e.g., the most relevant modes to experimentally observed conformational changes of low-symmetry structures are generally among the first 10-20 lowest-frequency non-rigid-body modes). The computed normal modes can be inspected by visualizing their animations, collectivities<sup>54</sup> and scores<sup>32</sup>, via the NMA results viewer (**Fig. 4B**) that can be obtained by

clicking on “Analyze Results” button in the project window. The collectivity<sup>54</sup> indicates how much the atoms or pseudoatoms move together with a given mode. It is computed for the entire complex without differentiating between its domains and it is normalized between  $1/N$  (maximally localized motion) and 1 (maximally collective motion), where  $N$  is the number of atoms or pseudoatoms. Highly collective low-frequency modes have been shown to be functionally relevant<sup>55-58</sup>. Therefore, we also provide the score<sup>32</sup> that combines the collectivity and frequency criteria, by penalizing the modes with higher frequencies and lower collectivities. The score is normalized between  $1/M$  (highest collectivity and lowest frequency motion) and 1 (lowest collectivity and highest frequency motion), where  $M$  is the number of modes.

The list of modes, with their collectivities and scores, can be obtained using “Display output Normal Modes?” in the NMA results viewer (**Fig. 4B**). The interface allows ordering the modes according to the increasing or decreasing values of the collectivity or score measures. The modes most relevant to the actual conformational change are expected to be among those with the lowest scores. The number of the lowest-score modes to inspect should be larger in the case of highly symmetrical structures (e.g., icosahedral-symmetry viruses) as highly collective modes may also exist at higher frequencies in such cases. To help the user decide which modes should certainly not be selected for the image analysis (Step 5), the modes with collectivities below 0.15 (very localized motions) are unchecked in the list of modes obtained via “Display output Normal Modes?” (**Fig. 4B**). The value of 0.15 for this collectivity threshold can be modified by the user (“Threshold on collectivity” parameter, **Fig. 4A**).

The motion trajectory along a normal mode is saved in a text file by concatenating the frames (in PDB format) of the coordinate displacement along this mode. The default number of

frames is 10 and the default amplitude of the displacement along normal mode is 50 (these parameters can be modified via the “Animation” tab, unopened in **Fig. 4A**). This trajectory file is animated with VMD,<sup>59</sup> by selecting “Display mode animation with VMD?” and specifying the mode to be animated in the “Mode number” field (**Fig. 4B**). VMD provides tools to save the motion trajectories in Animated GIF or MPEG movie formats that can be played with more standard movie players.

The modes can also be inspected by plotting the shifts of atoms (or pseudoatoms) along the specified mode and by plotting their maximum shifts over all modes via “Plot mode distance profile?” and “Plot max distance profile?”, respectively (**Fig. 4B**). For instance, a too large shift of one or a very few atoms (or pseudoatoms) with respect to the shift of other atoms (or other pseudoatoms) is typically a sign that the value of the atomic (or pseudoatomic) interaction cut-off distance used for computing normal modes is not optimal. In this case, normal modes should be recomputed using a modified value of this parameter.

**Step 3. Information:** At this step, entitled “Stop here or continue” in the HEMNMA menu (**Figs. 1-2**), we remind the user that HEMNMA software may also be used for NMA only (the processing can be stopped after computing and analyzing normal modes at Step 2). The following steps should be performed if aiming at analyzing conformational heterogeneity in images using the normal modes computed in Step 2.

**Step 4. Images:** This step allows importing images that will be analyzed with normal modes. The images should have a square shape. They should have a power of 2 pixels in each dimension when using the so-called “wavelet & splines” method for rigid-body alignment in combination with the elastic alignment (see Step 5 for more information about the alignment

methods). The larger the size of images, the longer the image analysis time will be. The image size can be reduced if the speed is an issue. The image size of 128x128 pixels is usually a good compromise.

**Step 5. Conformational distribution:** In this step, a combined elastic and rigid-body alignment is performed between particle images and the reference model to calculate the parameters of orientation, translation, and elastic deformation (amplitudes of normal modes) of the reference model that best describe the given particle image.

**Fig. 5A** shows the interface (“Input” tab) for selecting the images to be analyzed (“Input particles” parameter) and the normal modes to use for the image analysis (“Modes selection” parameter) from an entire list of computed normal modes (“Normal modes” parameter). The modes to use for the image analysis can be selected as the modes with the highest collectivities or least scores, with or without using previous knowledge about possible movements. The 6 lowest-frequency modes (modes 1-6, where the index of the mode corresponds to the mode order according to the increasing frequency) describe rigid-body movements of the structure and are often referred to as rigid-body normal modes. These modes should not be selected because rigid-body movements are taken into account through the iterative combined rigid-body and elastic alignment, as explained in the next paragraph. Usually, this alignment is performed with 1-10 selected normal modes excluding modes 1-6.

The interface for setting the parameters of the combined elastic and rigid-body alignment is shown in **Fig. 5B** (“Combined elastic and rigid-body alignment” tab). We use Powell’s optimization (more precisely, Powell’s trust region method<sup>60</sup>) to iteratively estimate the elastic deformation (normal-mode amplitudes) of the reference model by maximizing the similarity

between the particle image and the best matching projection of the deformed model (objective function, in optimization terminology). In each iteration of Powell’s method, an estimate of the normal-mode amplitudes is used to displace the atoms or the pseudoatoms along normal modes and the obtained, elastically deformed structure is converted into a density volume that is then rigid-body aligned with the particle image. The rigid-body alignment allows the determination of the orientation and translation parameters (rigid-body parameters) and the objective-function evaluation that is then used to better estimate the normal-mode amplitudes (elastic deformation parameters) for the next iteration of Powell’s method.

The values of the parameters of Powell’s method are set internally in the code and are not modifiable via the interface except for the scaling factor of the initial trust-region radius (“Elastic-alignment trust region scale” parameter in **Fig. 5B**). The default value of the scaling factor is 1, which generally works well. When expecting larger conformational changes, it may be interesting to increase the scaling factor (typically to a value between 1 and 2). The rigid-body alignment can be performed using one of the following two methods: 1) “projection matching” that stands for the standard reference-library-based (discrete) projection matching in real space (faster but less accurate method); and 2) “wavelets & splines” that stands for a discrete projection matching in wavelet space<sup>61</sup> followed by a continuous matching in 3D Fourier space based on spline interpolation<sup>62</sup> (slower but more accurate method; in particular, more robust to noise). The default choice is “wavelets & splines” (“Rigid-body alignment method” in **Fig. 5B**). The angular sampling step<sup>61</sup> for the reference projection library computation is 10° by default and can be modified via the interface (“Discrete angular sampling” parameter expressed in degrees, **Fig. 5B**).

This combined rigid-body and elastic alignment is the task that requires the most



computation. Therefore, it was MPI parallelized to simultaneously process  $N$  particle images using  $N$  MPI cores on CPU computers, clusters or supercomputers. The user specifies the number of MPI cores to be used (“MPI” parameter in **Fig. 5A**). On an Intel Xeon CPU at 2.9 GHz, the analysis of one particle image of size 128×128 pixels with 3, 6, and 9 normal modes may respectively take 2, 5, and 10 min when using “projection matching” rigid-body alignment and 5, 13, and 25 min when using “wavelets & splines” rigid-body alignment.

Normal-mode amplitudes resulting from the image analysis with normal modes can be visualized in 1 dimension (histogram) as well as in 2 or 3 dimensions (**Fig. 5C**). Points in a 2D and 3D space correspond to images with assigned orientations, translations, and normal-mode displacement amplitudes with respect to the reference model. Close points in this space correspond to similar conformations and vice versa.

**Step 6. Dimension reduction, clusters, and trajectories:** The normal-mode displacement amplitudes obtained by image analysis in Step 5 can be projected onto a space of lower dimension (1D, 2D, or 3D) (**Fig. 6**), using the technique of principal component analysis (PCA) or one of several other dimension reduction techniques.<sup>63</sup> To perform the dimension reduction, the user specifies the conformational distribution to be analyzed (the image analysis results obtained in Step 5), the desired lower dimension, and one of the available linear or non-linear dimensionality reduction techniques via “Conformational distribution”, “Reduced dimension”, and “Dimensionality reduction method” fields, respectively (**Fig. 6A**). The “Extra params” field (**Fig. 6A**) allows modifying the parameters of the dimensionality reduction techniques (the parameters and their default values are listed in the help message linked to the “Dimensionality reduction method” field). The available dimensionality reduction techniques are PCA, Kernel PCA, Probabilistic PCA, Local Tangent Space

Alignment (LTSA), Linear LTSA, Diffusion Map, Linearity Preserving Projection, Laplacian Eigenmap, Hessian Locally Linear Embedding, Stochastic Proximity Embedding, and Neighborhood Preserving Embedding.<sup>63</sup> The dimension reduction viewer (**Fig. 6B**) can be obtained by clicking on “Analyze Results” button in the project window. It allows displaying the axes of the low-dimensional space (**Fig. 6C**) specified via the “Display normal-mode amplitudes in the low-dimensional space” field (**Fig. 6B**). Also, the dimension reduction viewer (**Fig. 6B**) allows opening the tool for making animations (“Trajectories Tool”, **Fig. 7A**) and the tool for computing 3D reconstructions (“Clustering Tool”, **Fig. 8A**).

While the advantage of non-linear dimensionality reduction techniques is better suitability to non-linear manifold data representations, the advantage of linear dimensionality reduction techniques (PCA, Linear LTSA, Linearity Preserving Projection, Probabilistic PCA, and Neighborhood Preserving Embedding) is that they allow mapping from the low-dimensional space back to the original space. HEMNMA uses this property to generate animated trajectories of conformational changes. More precisely, “Trajectories Tool” (**Fig. 7A**) allows recording and visualizing the displacement of the reference model in the low-dimensional space. The specification of the displacement trajectory to be animated requires the user’s interaction. The outlier data points can be removed by providing logical (Boolean) expressions (“Expression” field in **Fig. 7A**). The trajectory is specified by the coordinates of 10 points in the low-dimensional space (red points in **Fig. 7B**). The user may select all of these 10 points (by clicking on the plot to select each point) or 8 points may be automatically placed on a line between 2 points selected by the user (the first and last points of the trajectory). The position of the initially placed points may be changed by dragging the points. As in the case of animations of normal modes in Step 2, the trajectory is saved in a text file by concatenating the PDB-format frames of the coordinate displacement. The trajectory can be

animated with VMD (**Fig. 7C**) or saved with VMD in Animated GIF or MPEG movie formats for playing with other movie players.

Additionally, conformational changes can be analyzed in terms of 3D reconstructions from images in the low-dimensional space. “Clustering Tool” (**Fig. 8A**) allows making groups of close points (**Fig. 8B**), corresponding to images with similar conformational states, and computing 3D reconstructions from these groups (**Fig. 8C**). A group of points is specified by providing logical (Boolean) expressions (“Expression” field in **Fig. 8A**) or by clicking on the plot and dragging to add points to the group. Each point in the selected group of points is denoted by a yellow circle (**Fig. 8B**). Clicking on “Create Cluster” button (**Fig. 8A**) saves the selected group of points and performs 3D reconstruction from this group. The saved group can be inspected by displaying the 3D reconstruction results (slices and isosurface of the reconstructed volume, **Fig. 8C**) and the images present in the group (the images are not shown in **Fig. 8**). HEMNMA uses a fast Fourier-space method for 3D reconstruction from the selected group of points. The 3D reconstruction can be performed with other reconstruction methods (available in Scipion or other software packages), using the output metadata (text) file with the rigid-body and elastic alignment parameters corresponding to the selected group of points. The 3D reconstructions represent the average states from the corresponding groups of images. One should make as homogeneous groups of points as possible with a sufficient number of points per group in order to obtain high-resolution 3D reconstructions (note that the example in **Fig. 8B** does not show the most optimal grouping of points).

As explained, HEMNMA provides a full conformational variability landscape before grouping images into 3D reconstructions (images with similar conformations on the landscape). As such, HEMNMA can detect the conformations and motions that are

undetectable with traditional classification-based approaches, which was extensively studied elsewhere.<sup>31</sup>

## **PERSPECTIVES**

A full description of conformational heterogeneity is important for both biology and drug design. Cryo-EM has been under continuous development since its beginning (early '80s), which was even accelerated with the latest instrumental developments including direct electron detector devices (DDD cameras). Allowing a routine near-atomic resolution of structures and a routine full description of conformational variability are currently two main challenges. Indeed, the recent cryo-EM advances such as DDD cameras, phase plates, and sample motion correction have reduced noise and improved contrast in images, which makes elucidation of conformations from images more accurate. Methods for a full description of continuous conformational variability are being developed. However, these methods will need to be more efficient and user-friendly to allow routine use of cryo-EM for such studies. HEMNMA is user-friendly software that has been developed for determining the full distribution of continuous conformational variability from cryo-EM images. Still, HEMNMA needs to be made faster in order to allow an efficient high-resolution description of this variability. For instance, the analysis of  $10^4$  particle images based on the elastic alignment with 6 normal modes and the “wavelets & splines” rigid-body alignment may take 17 h using 128 MPI cores on 2.9 GHz Intel Xeon CPUs, which means that the analysis of  $10^6$  particle images would require around 70 days of use of the same 128 MPI cores. To get a high-resolution description of the full conformational landscape of a complex, several millions of particle images would need to be analyzed. A combination of HEMNMA with a deep learning approach will be implemented in the future to speed up processing of such large numbers of particle images.

HEMNMA is well suited to compact-support particle shapes and may be adapted to deal with helical symmetry particles in the future.

## ACKNOWLEDGMENTS

We acknowledge the help of José Miguel de la Rosa-Trevín, CNB-CSIC (currently at SciLifeLab, Stockholm University) with porting the HEMNMA code from Xmipp to Scipion V1.0; the support of the French National Research Agency – ANR (ANR-19-CE11-0008-01 to SJ); the access to the HPC resources of CINES and IDRIS granted by GENCI (2019-A0070710998 to SJ); and the support of French Microscopy Society – Sfμ and Sorbonne University – SU (2019 Sfμ Master internship grant and 2019 SU PhD scholarship grant to MH).

## REFERENCES

1. Zhang X, Settembre E, Xu C, Dormitzer PR, Bellamy R, Harrison SC, Grigorieff N (2008) Near-atomic resolution using electron cryomicroscopy and single-particle reconstruction. *Proc Natl Acad Sci U S A* 105:1867-1872.
2. Yu X, Jin L, Zhou ZH (2008) 3.88 Å structure of cytoplasmic polyhedrosis virus by cryo-electron microscopy. *Nature* 453:415-419.
3. Liao M, Cao E, Julius D, Cheng Y (2013) Structure of the trpv1 ion channel determined by electron cryo-microscopy. *Nature* 504:107-112.
4. Bai XC, Yan C, Yang G, Lu P, Ma D, Sun L, Zhou R, Scheres SH, Shi Y (2015) An atomic structure of human gamma-secretase. *Nature* 525:212-217.
5. Bartesaghi A, Merk A, Banerjee S, Matthies D, Wu X, Milne JL, Subramaniam S (2015) 2.2 Å resolution cryo-em structure of beta-galactosidase in complex with a cell-

- permeant inhibitor. *Science* 348:1147-1151.
6. Nguyen TH, Galej WP, Bai XC, Oubridge C, Newman AJ, Scheres SH, Nagai K (2016) Cryo-em structure of the yeast u4/u6.U5 tri-snrnp at 3.7 a resolution. *Nature* 530:298-302.
  7. Natchiar SK, Myasnikov AG, Kratzat H, Hazemann I, Klaholz BP (2017) Visualization of chemical modifications in the human 80s ribosome structure. *Nature* 551:472-477.
  8. Zhang J, Ma J, Liu D, Qin S, Sun S, Zhao J, Sui SF (2017) Structure of phycobilisome from the red alga *griffithsia pacifica*. *Nature* 551:57-63.
  9. Liu Y, Gonen S, Gonen T, Yeates TO (2018) Near-atomic cryo-em imaging of a small protein displayed on a designed scaffolding system. *Proceedings of the National Academy of Sciences* 115:3362.
  10. Fan X, Wang J, Zhang X, Yang Z, Zhang J-C, Zhao L, Peng H-L, Lei J, Wang H-W (2019) Single particle cryo-em reconstruction of 52 kda streptavidin at 3.2 angstrom resolution. *Nature Communications* 10:2386.
  11. Khatter H, Myasnikov AG, Natchiar SK, Klaholz BP (2015) Structure of the human 80s ribosome. *Nature* 520:640-645.
  12. He Y, Yan C, Fang J, Inouye C, Tjian R, Ivanov I, Nogales E (2016) Near-atomic resolution visualization of human transcription promoter opening. *Nature* 533:359-365.
  13. Banerjee S, Bartesaghi A, Merk A, Rao P, Bulfer SL, Yan Y, Green N, Mroczkowski B, Neitz RJ, Wipf P et al. (2016) 2.3 a resolution cryo-em structure of human p97 and mechanism of allosteric inhibition. *Science* 351:871-875.
  14. Abeyrathne PD, Koh CS, Grant T, Grigorieff N, Korostelev AA (2016) Ensemble cryo-em uncovers inchworm-like translocation of a viral ires through the ribosome. *Elife* 5.
  15. van Heel M, Stoffer-Meilicke M (1985) Characteristic views of e. Coli and b.

- Stearothermophilus 30s ribosomal subunits in the electron microscope. *Embo J* 4:2389-2395.
16. Penczek PA, Frank J, Spahn CM (2006) A method of focused classification, based on the bootstrap 3d variance analysis, and its application to ef-g-dependent translocation. *J Struct Biol* 154:184-194.
  17. Scheres SH, Gao H, Valle M, Herman GT, Eggermont PP, Frank J, Carazo JM (2007) Disentangling conformational states of macromolecules in 3d-em through likelihood optimization. *Nat Methods* 4:27-29.
  18. Elad N, Clare DK, Saibil HR, Orlova EV (2008) Detection and separation of heterogeneity in molecular complexes by statistical analysis of their two-dimensional projections. *J Struct Biol* 162:108-120.
  19. Simonetti A, Marzi S, Myasnikov AG, Fabbretti A, Yusupov M, Gualerzi CO, Klaholz BP (2008) Structure of the 30s translation initiation complex. *Nature* 455:416-420.
  20. Sorzano COS, Bilbao-Castro JR, Shkolnisky Y, Alcorlo M, Melero R, Caffarena-Fernandez G, Li M, Xu G, Marabini R, Carazo JM (2010) A clustering approach to multireference alignment of single-particle projections in electron microscopy. *Journal of Structural Biology* 171:197-206.
  21. Penczek PA, Kimmel M, Spahn CM (2011) Identifying conformational states of macromolecules by eigen-analysis of resampled cryo-em images. *Structure* 19:1582-1590.
  22. Scheres SH (2012) A bayesian view on cryo-em structure determination. *J Mol Biol* 415:406-418.
  23. Lyumkis D, Brilot AF, Theobald DL, Grigorieff N (2013) Likelihood-based classification of cryo-em images using frealign. *J Struct Biol* 183:377-388.
  24. Fischer N, Neumann P, Konevega AL, Bock LV, Ficner R, Rodnina MV, Stark H (2015)

- Structure of the e. Coli ribosome-ef-tu complex at  $<3$  Å resolution by cs-corrected cryo-em. *Nature* 520:567-570.
25. Scheres SHW, Crowther RA. 2016. Chapter six - processing of structurally heterogeneous cryo-em data in relion. *Methods in enzymology*. Academic Press. p. 125-157.
  26. Punjani A, Rubinstein JL, Fleet DJ, Brubaker MA (2017) Cryosparc: Algorithms for rapid unsupervised cryo-em structure determination. *Nature Methods* 14:290.
  27. Haselbach D, Komarov I, Agafonov DE, Hartmuth K, Graf B, Dybkov O, Urlaub H, Kastner B, Luhrmann R, Stark H (2018) Structure and conformational dynamics of the human spliceosomal b(act) complex. *Cell* 172:454-464 e411.
  28. Grant T, Rohou A, Grigorieff N (2018) Cistem, user-friendly software for single-particle image processing. *eLife* 7:e35383.
  29. Katsevich E, Katsevich A, Singer A (2015) Covariance matrix estimation for the cryo-em heterogeneity problem. *SIAM J Imag Sci* 8:126-185.
  30. Tagare HD, Kucukelbir A, Sigworth FJ, Wang H, Rao M (2015) Directly reconstructing principal components of heterogeneous particles from cryo-em images. *J Struct Biol* 191:245-262.
  31. Jin Q, Sorzano CO, de la Rosa-Trevin JM, Bilbao-Castro JR, Nunez-Ramirez R, Llorca O, Tama F, Jonic S (2014) Iterative elastic 3d-to-2d alignment method using normal modes for studying structural dynamics of large macromolecular complexes. *Structure* 22:496-506.
  32. Sorzano CO, de la Rosa-Trevin JM, Tama F, Jonic S (2014) Hybrid electron microscopy normal mode analysis graphical interface and protocol. *J Struct Biol* 188:134-141.
  33. Dashti A, Schwander P, Langlois R, Fung R, Li W, Hosseinizadeh A, Liao HY, Pallesen J, Sharma G, Stupina VA et al. (2014) Trajectories of the ribosome as a brownian nanomachine. *Proc Natl Acad Sci U S A* 111:17492-17497.



34. Frank J, Ourmazd A (2016) Continuous changes in structure mapped by manifold embedding of single-particle data in cryo-em. *Methods* 100:61-67.
35. Anden J, Singer A (2018) Structural variability from noisy tomographic projections. *SIAM J Imaging Sci* 11:1441-1492.
36. Dashti A, Hail DB, Mashayekhi G, Schwander P, des Georges A, Frank J, Ourmazd A  
Conformational dynamics and energy landscapes of ligand binding in ryr1. *bioRxiv* 167080.
37. Lederman RR, Singer A Continuously heterogeneous hyper-objects in cryo-em and 3-d movies of many temporal dimensions. *arXiv* 1704.02899.
38. Vargas J, Alvarez-Cabrera AL, Marabini R, Carazo JM, Sorzano CO (2014) Efficient initial volume determination from electron microscopy images of single particles. *Bioinformatics* 30:2891-2898.
39. Scheres SH (2012) Relion: Implementation of a bayesian approach to cryo-em structure determination. *J Struct Biol* 180:519-530.
40. Reboul CF, Eager M, Elmlund D, Elmlund H (2018) Single-particle cryo-em-improved ab initio 3d reconstruction with simple/prime. *Protein Sci* 27:51-61.
41. Klaholz BP (2015) Structure sorting of multiple macromolecular states in heterogeneous cryo-em samples by 3d multivariate statistical analysis. *Open J Stat* 5:820-836.
42. Jonic S (2017) Computational methods for analyzing conformational variability of macromolecular complexes from cryo-electron microscopy images. *Curr Opin Struct Biol* 43:114-121.
43. Cossio P, Hummer G (2018) Likelihood-based structural analysis of electron microscopy images. *Curr Opin Struct Biol* 49:162-168.
44. Sorzano COS, Jimenez A, Mota J, Vilas JL, Maluenda D, Martinez M, Ramirez-Aportela E, Majtner T, Segura J, Sanchez-Garcia R et al. (2019) Survey of the analysis of

- continuous conformational variability of biological macromolecules by electron microscopy. *Acta Crystallogr F Struct Biol Commun* 75:19-32.
45. Serna M (2019) Hands on methods for high resolution cryo-electron microscopy structures of heterogeneous macromolecular complexes. *Frontiers in molecular biosciences* 6:33-33.
  46. Sorzano CO, Marabini R, Velazquez-Muriel J, Bilbao-Castro JR, Scheres SH, Carazo JM, Pascual-Montano A (2004) Xmipp: A new generation of an open-source image processing package for electron microscopy. *J Struct Biol* 148:194-204.
  47. de la Rosa-Trevin JM, Quintana A, del Cano L, Zaldivar A, Foche I, Gutierrez J, Gomez-Blanco J, Burguet-Castell J, Cuenca-Alba J, Abrishami V et al. (2016) Scipion: A software framework toward integration, reproducibility and validation in 3d electron microscopy. *Journal of Structural Biology* 195:93-99.
  48. Jonic S, Sorzano COS (2016) Coarse-graining of volumes for modeling of structure and dynamics in electron microscopy: Algorithm to automatically control accuracy of approximation. *IEEE J Sel Topics Signal Process* 10:161 - 173.
  49. Jonic S, Vargas J, Melero R, Gomez-Blanco J, Carazo JM, Sorzano CO (2016) Denoising of high-resolution single-particle electron-microscopy density maps by their approximation using three-dimensional gaussian functions. *J Struct Biol* 194:423-433.
  50. Pettersen EF, Goddard TD, Huang CC, Couch GS, Greenblatt DM, Meng EC, Ferrin TE (2004) Ucsf chimera--a visualization system for exploratory research and analysis. *J Comput Chem* 25:1605-1612.
  51. Sanchez Sorzano CO, Alvarez-Cabrera AL, Kazemi M, Carazo JM, Jonic S (2016) Structmap: Elastic distance analysis of electron microscopy maps for studying conformational changes. *Biophys J* 110:1753-1765.
  52. Tirion MM (1996) Large amplitude elastic motions in proteins from a single-parameter,

- atomic analysis. *Phys Rev Lett* 77:1905-1908.
53. Tama F, Gadea FX, Marques O, Sanejouand YH (2000) Building-block approach for determining low-frequency normal modes of macromolecules. *Proteins* 41:1-7.
54. Bruschweiler R (1995) Collective protein dynamics and nuclear spin relaxation. *J Chem Phys* 102:3396-3403.
55. Delarue M, Dumas P (2004) On the use of low-frequency normal modes to enforce collective movements in refining macromolecular structural models. *Proc Natl Acad Sci U S A* 101:6957-6962.
56. Suhre K, Navaza J, Sanejouand YH (2006) Norma: A tool for flexible fitting of high-resolution protein structures into low-resolution electron-microscopy-derived density maps. *Acta Crystallogr D Biol Crystallogr* 62:1098-1100.
57. Tama F, Miyashita O, Brooks CL, 3rd (2004) Flexible multi-scale fitting of atomic structures into low-resolution electron density maps with elastic network normal mode analysis. *J Mol Biol* 337:985-999.
58. Wang Y, Rader AJ, Bahar I, Jernigan RL (2004) Global ribosome motions revealed with elastic network model. *J Struct Biol* 147:302-314.
59. Humphrey W, Dalke A, Schulten K (1996) Vmd: Visual molecular dynamics. *J Mol Graph* 14:33-38, 27-38.
60. Powell MJD (2002) Uobyqa: Unconstrained optimization by quadratic approximation. *Mathematical Programming* 92:555-582.
61. Sorzano CO, Jonic S, El-Bez C, Carazo JM, De Carlo S, Thevenaz P, Unser M (2004) A multiresolution approach to orientation assignment in 3d electron microscopy of single particles. *J Struct Biol* 146:381-392.
62. Jonic S, Sorzano CO, Thevenaz P, El-Bez C, De Carlo S, Unser M (2005) Spline-based image-to-volume registration for three-dimensional electron microscopy.

Ultramicroscopy 103:303-317.

63. Van der Maaten L, Postma E, Van den Herik J. 2009. Dimensionality reduction: A comparative review. Tilburg centre for Creative Computing, Tilburg university. No. 2009-005.

## FIGURE LEGENDS

**Figure 1:** HEMNMA project tree for the case of analyzing images using an atomic structure as the reference model.

**Figure 2:** HEMNMA project tree for the case of analyzing images using a density volume as the reference model. The volume is converted into 3D Gaussian functions (the so-called pseudoatoms).

**Figure 3:** Volume-to-pseudoatoms conversion. (A) Dialog box. (B) Superposition of an input density volume (transparent grey) and its pseudoatomic representation (spheres where the amplitudes of Gaussian functions are color-coded from white (smallest) to red (largest)). In this example, the Gaussian-function standard deviation (pseudoatom radius) and the desired volume approximation error are 2 voxels and 5%, respectively. The desired volume approximation error parameter is hidden by default (“Expert Level” is set to “Normal”) and can be visualized by setting “Expert Level” to “Advanced” (panel A). The volume-to-pseudoatoms conversion is parallelized using threads (parallel processes sharing the same memory) and the number of threads to use can be specified in the “Threads” field (panel A). The masking options to remove background noise from the input volume are explained in the help message under the corresponding question mark (panel A).

**Figure 4:** Normal mode analysis and visualization. (A) Dialog box (note that “Cut-off mode” selected here is “relative” and that the interface allows selecting the cut-off mode “absolute”

as well). (B) Normal modes viewer allowing visualizing animations of normal modes (using “Display mode animation with VMD?” for the mode specified in “Mode number”) and checking their collectivities and scores (using “Display output Normal Modes?”).

**Figure 5:** Image analysis with normal modes. (A)-(B) Two sections of the dialog box, regarding input (panel A) and alignment (panel B). (C) Example of visualization of image analysis results in 3 dimensions (amplitudes along 3 normal modes). The image analysis task is parallelized using MPI protocol and the number of MPI cores (parallel processes that do not necessarily share the same memory) can be specified in the MPI field (panel A). In panel B, the default choice for “Rigid-body alignment method” is “wavelets & splines” (the alternative method is “projection matching”) and the elastic alignment is performed by Powell’s trust region method. In panel C, each point corresponds to a particle image with assigned orientations, translations, and normal-mode displacement amplitudes with respect to the reference model (close points correspond to similar conformations and vice versa).

**Figure 6:** Dimension reduction. (A) Dialog box for projecting normal-mode amplitudes computed by image analysis onto a space of lower dimension using principal component analysis (here, 2 in “Reduced dimension” means projecting onto a 2D space) or one of several other dimension reduction methods selected via “Dimensionality reduction method” (the default values of parameters of these methods are provided in the help message displayable by clicking on the corresponding question mark). (B) Dimension reduction viewer allowing visualizing normal-mode amplitudes in the low-dimensional space (here, a 2D space specified by axes 1 and 2) as well as opening the clustering (grouping) and trajectories tools. (C) Example of projecting normal-mode amplitudes onto a low-dimensional space (here, 2D space specified in panel B). In panel A, the methods available for the dimension reduction are Principal Component Analysis (PCA), Kernel PCA, Probabilistic PCA, Local Tangent Space Alignment (LTSA), Linear LTSA, Diffusion Map, Linearity Preserving Projection, Laplacian

Eigenmap, Hessian Locally Linear Embedding, Stochastic Proximity Embedding, and Neighborhood Preserving Embedding.

**Figure 7:** Animations after image analysis. (A) Trajectories Tool. (B) Example of trajectory (10 red points, see the text for more details). (C) Trajectory animation using VMD.

**Figure 8:** 3D reconstructions after image analysis. (A) Clustering Tool for grouping close points (images with similar conformations). (B) Example of selecting a group of points (circled in yellow) for 3D reconstruction. (C) Chimera superposition of three 3D reconstructions (yellow, cyan and gray isosurfaces) from the corresponding groups of images denoted in panel B (points surrounded by yellow, cyan and gray ellipses, panel B), with the most dominant motion shown by red arrows. The group of points saved using “Create Cluster” button (panel A) appears as a new box in the project tree (panel C), possibly after using “Refresh” button from the project window. The saved group of points can be inspected using “Analyze Results” button from the project window, which shows images in this group (not shown here) and slices of the volume reconstructed from these images (panel C). The volume isosurface can be visualized by clicking on Chimera icon in the slices display menu (panel C).

View

HEMNMA (Plugin)

HEMNMA (Plugin)

1. Reference model

★ a. Import PDB

★ b1. Import volume

P b2. Convert volume to PDB

2. Normal mode analysis

P Modes analysis & visualization

3. Stop here or continue

4. Images

★ Import particles

P Resize particles (optional)

5. Conformational distribution

P Image analysis with normal modes

6. Dimension reduction, clusters, and trajectories

P 3D reconstructions from image clusters, animated trajectories

Copy

Delete

Export

Export & upload

Labels

View: Tree

Refresh

PROJECT

scipion - import atomic structure finished

scipion - import particles finished

continuousflex - nma analysis finished

continuousflex - nma alignment finished

continuousflex - nma dimred finished

To extend the tree with the boxes containing 3D reconstruction results

To see the results of a selected box

Analyze Results

ViewHEMNMA (Plugin)

HEMNMA (Plugin)

1. Reference model

a. Import PDB

b1. Import volume

b2. Convert volume to PDB

2. Normal mode analysis

Modes analysis & visualization

3. Stop here or continue

4. Images

Import particles

Resize particles (optional)

5. Conformational distribution

Image analysis with normal modes

6. Dimension reduction, clusters, and trajectories

3D reconstructions from image clusters, animated trajectories

CopyDeleteExportExport & uploadLabels

PROJECT

scipion - import volumes finished

continuousflex - convert to pseudoatoms finished

continuousflex - nma analysis finished

continuousflex - nma alignment finished

continuousflex - nma dimred finished

scipion - import particles finished

View:TreeRefresh

To extend the tree with the boxes containing 3D reconstruction results

To see the results of a selected box

Analyze Results



A

Protocol Run: FlexProtConvertToPseudoAtoms

Protocol: continuousflex - convert to pseudoatoms finished Cite Help

**Run**

Run name continuousflex - convert to pse  Comment

Run mode ☒ Continue ☐ Restart ? Host localhost

Parallel Threads 4 ? Use queue? ☐ Yes ☒ No ?

Wait for ?

Expert Level ☐ Normal ☒ Advanced

**Input**

Input structure scipion - import volumes.outputVolume

Mask mode none ?

Pseudoatom radius (vox) 2.0 ?

Volume approximation error(%) 5.0 ?

Close Save Execute

B




A

## Protocol: continuousflex - nma analysis

finished  Cite  Help

## Run

Run name continuousflex - nma analysis Comment 

Run mode

☒ Continue☐ Restart 

Host

localhost 

Use queue?

☐ Yes☒ No Wait for 








Expert Level

☐ Normal☒ Advanced

Normal Mode Analysis

Animation

## Normal Mode Analysis

Input structure continuousflex - convert to pseudoatoms.outputPdb    Number of modes 20 Cut-off mode relative  Cut-off percentage 95.0 Threshold on collectivity 0.15 Number of residues per RTB block (for atomic structures) 10  Close Save Execute

B

## Protocol: continuousflex - viewer nma

 Cite  Help

Visualization

## Visualization

Display output Normal Modes? Plot max distance profile?  

Single mode

Mode number

7

Display mode animation with VMD? Plot mode distance profile?  

A

## Protocol: continuousflex - nma alignment

finished [Cite](#) [Help](#)

## Run

Run name continuousflex - nma alignmen

Comment

Run mode

☒ Continue☐ Restart

?

Host

localhost

Parallel

MPI 8

?

Use queue?

☐ Yes☒ No

?

Wait for

?

Expert Level

☐ Normal☒ Advanced

Input

Combined elastic and rigid-body alignment

Input

Normal modes continuousflex - nma analysis.outputModes

Q ?

Modes selection 7-9

?

Input particles scipion - import particles.outputParticles

Q ?

Precomputed results (for development)

?

Close

Save

Execute

B

Expert Level

☐ Normal☒ Advanced

Input

Combined elastic and rigid-body alignment

## Combined elastic and rigid-body alignment

Elastic-alignment trust region scale

1

?

Rigid-body alignment method

wavelets &amp; splines

?


Discrete angular sampling (deg)

10.0

?

C

Normal-mode amplitudes: Mode 7 Mode 8 Mode 9

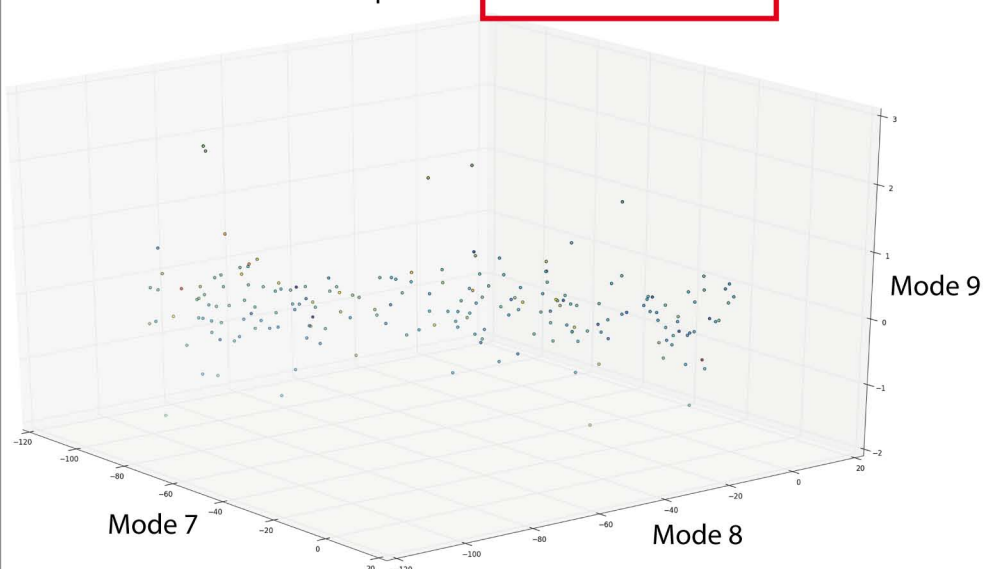

A

## Protocol Run: FlexProtDimredNMA

Protocol: continuousflex - nma dimred

finished



Run

Run name continuousflex - nma dimred

Comment

Run mode

☒ Continue☐ Restart

Host

localhost

Use queue?

☐ Yes☒ No

Wait for

Input

Input

Conformational distribution continuousflex - nma alignment

Dimensionality reduction method Principal Component Analysis (PCA)

Extra params

Reduced dimension

2

Close

Save

Execute

B

## Protocol Viewer: FlexDimredNMAViewer

Protocol: continuousflex - viewer nma dimred



Visualization

Visualization

Display normal-mode amplitudes in the low-dimensional space


12

Open clustering tool?


Open trajectories tool?

C


Normal-mode amplitudes in low-dimensional space: Axis 1 vs Axis 2

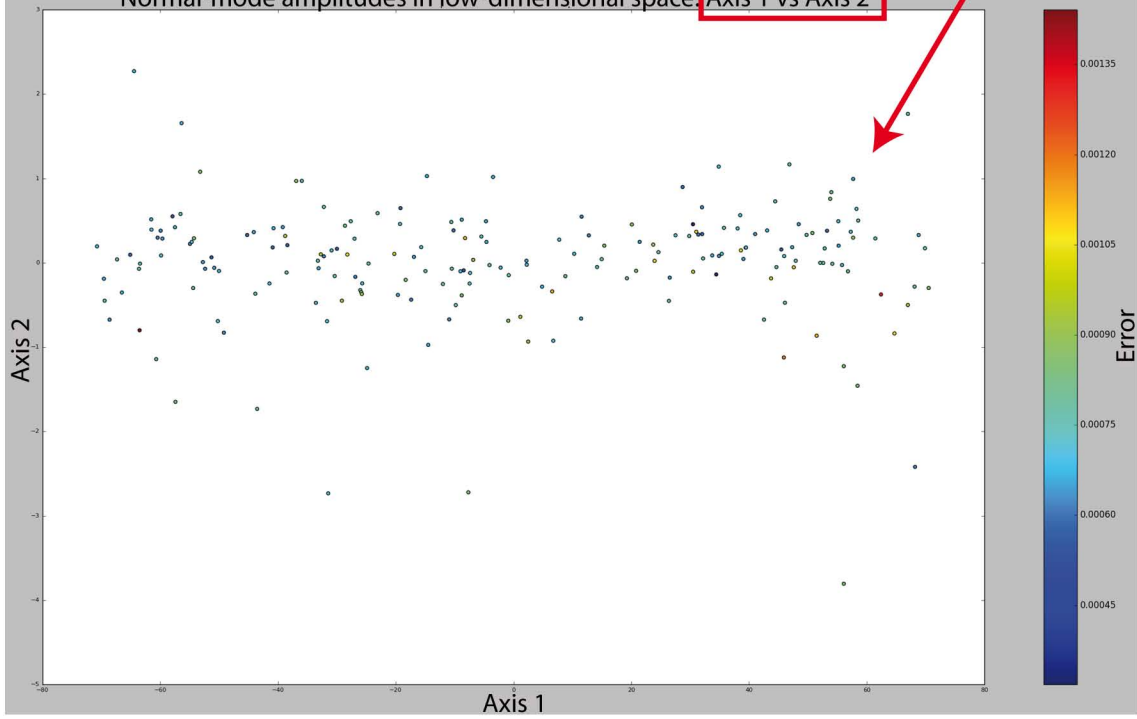



A

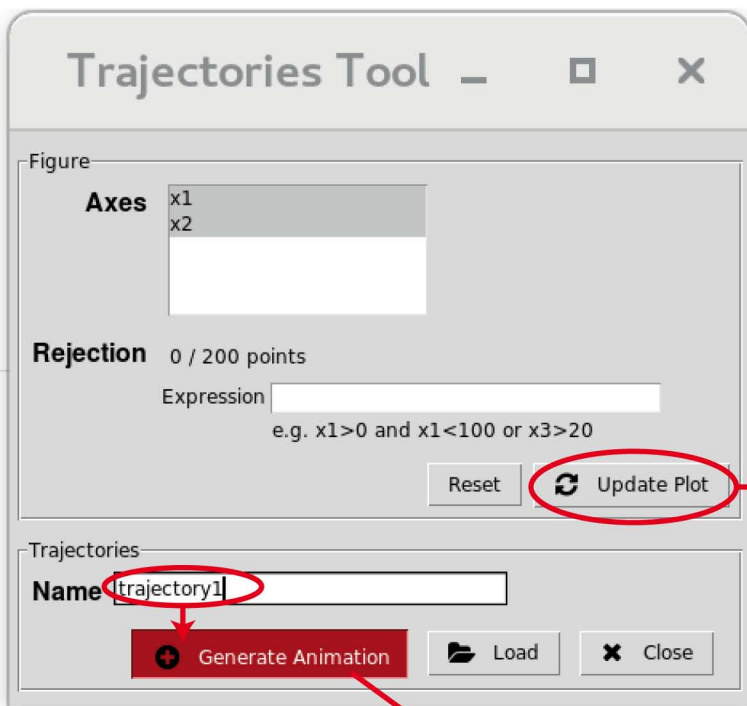


B

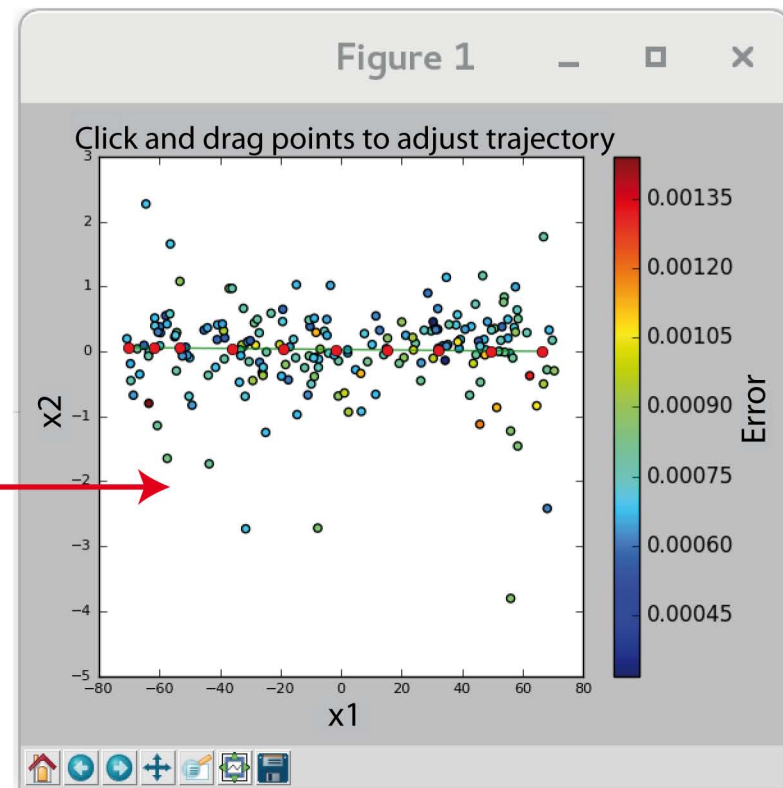


C

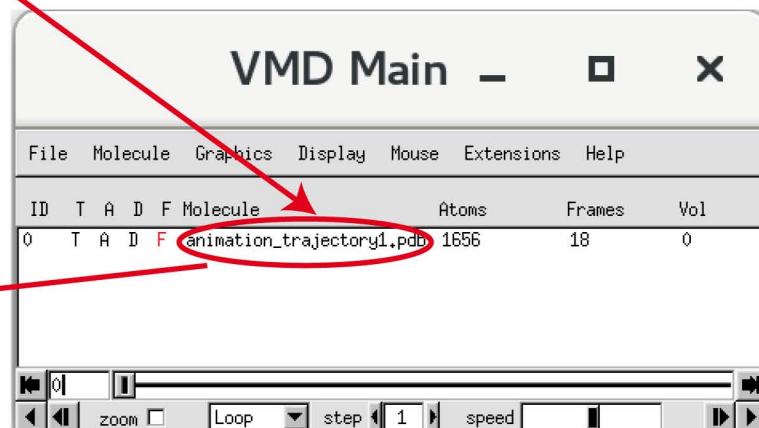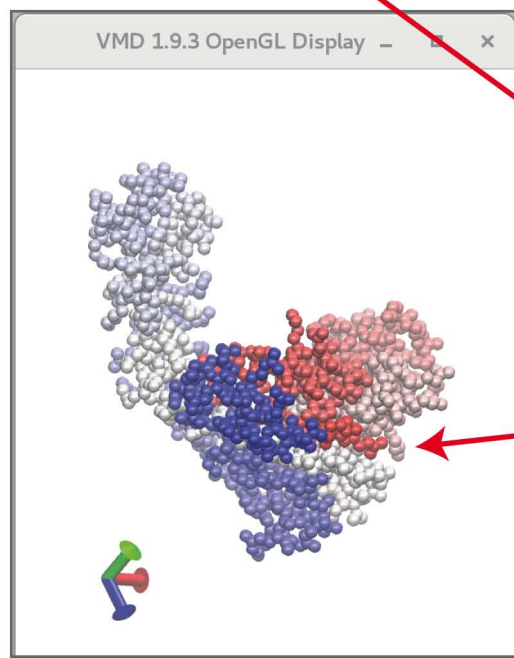

A

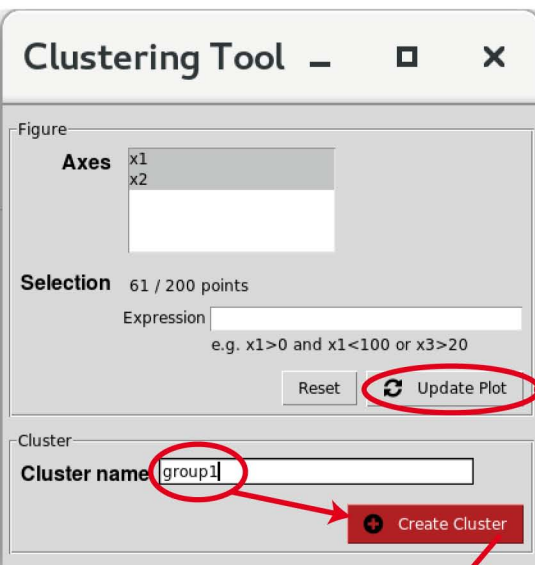

B

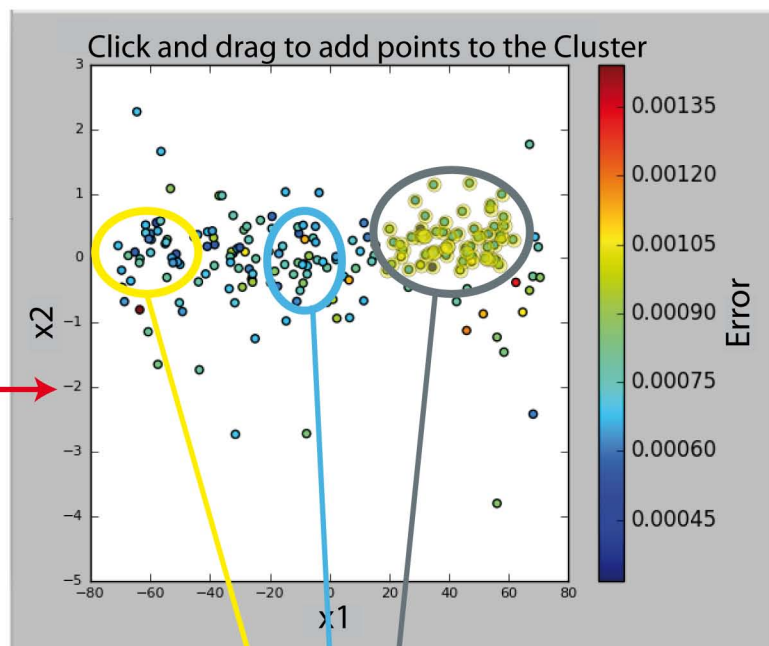

C

