



**HAL**  
open science

# Flash Chromatography System: A Practical Tool for Demonstrating the Influence of Column Characteristics on Chromatographic Resolution

Amaury Kasprowiak, Francine Cazier-Dennin, Pierre-Edouard Danjou

► **To cite this version:**

Amaury Kasprowiak, Francine Cazier-Dennin, Pierre-Edouard Danjou. Flash Chromatography System: A Practical Tool for Demonstrating the Influence of Column Characteristics on Chromatographic Resolution. *Journal of Chemical Education*, 2020, 97 (4), pp.1145-1150. 10.1021/acs.jchemed.9b00929 . hal-02865856

**HAL Id: hal-02865856**

**<https://hal.science/hal-02865856v1>**

Submitted on 21 Oct 2020

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

---

# Flash chromatography system: a practical tool for demonstrating the influence of column characteristics on chromatographic resolution.

5 Amaury Kasprowiak<sup>a</sup>, Francine Cazier-Dennin<sup>a,b</sup>, Pierre-Edouard Danjou<sup>a,b,\*</sup>

<sup>a</sup> Département de Chimie, Université du Littoral Côte d'Opale, 220 avenue de l'université, Dunkerque, France

<sup>b</sup> Unité de Chimie Environnementale et Interactions sur le Vivant, EA 4492, Université du Littoral Côte  
10 d'Opale, 145 Avenue Maurice Schumann, MREI 1, Dunkerque, France


\* danjou@univ-littoral.fr

## ABSTRACT

A flash purification system using prepacked cartridges containing different silica weight and particle sizes was employed to demonstrate the influence of column characteristics on  
15 chromatographic resolution. This chromatographic experiment was a helpful introduction to chromatographic courses where the influence of particle sizes was addressed only theoretically and it allows students to experience the concept in an active learning style experiment. After running a predefined set of experiments to understand concepts, students were asked to design their own experiment to achieve the best chromatographic resolution. Moreover, the difference of mass versus  
20 concentration detectors was highlighted using Ultraviolet (UV) and Evaporative Light Scattering Detectors (ELSD). Questionnaires were also used to define the students' initial knowledge and to assess on what they learned following this lab.

---

## GRAPHICAL ABSTRACT

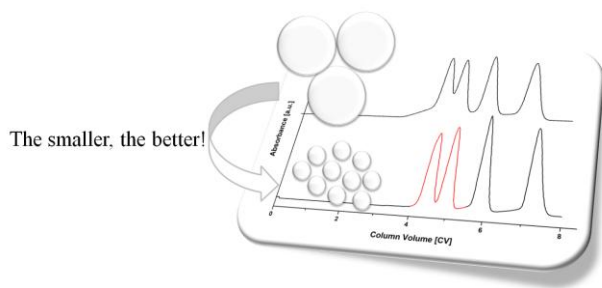

### 25 KEYWORDS


Chromatography, Analytical Chemistry, Laboratory Instruction, Hands-On Learning/Manipulatives, Upper-Division Undergraduate, Laboratory Equipment / Apparatus

The concept of flash chromatography was first introduced in 1978 with the promise of separating  
30 several grams of compound in a very short time<sup>1</sup> (15 minutes instead of 2-3 hours). While this  
technique is commonly used by organic chemists, it has not undergone any major evolution over  
several decades and remains the method of choice for the purification of organic compounds. Some  
minor modifications were proposed to facilitate access to this technique in teaching laboratories<sup>2,3,4</sup>.  
However, the emergence of automated flash chromatography systems like Reveleris (Buchi)<sup>5</sup>, Isolera  
35 (Biotage)<sup>6</sup> or Puriflash (Interchim)<sup>7</sup>, has deeply changed the ways compounds are purified and open  
new teaching opportunities in the organic synthesis as well as in analytical chemistry. It is obvious  
that such devices will find their usefulness in the organic chemistry laboratory in  
complement/replacement of traditional column chromatography<sup>8</sup> since it became increasingly popular  
in companies as well as in academic research labs. Moreover, the use of prepacked silica cartridges  
40 prevents the handling of hazardous silica gel<sup>9</sup> by the students, leading to a safer experiment.

Nevertheless, it could also be interesting to introduce flash chromatography systems in a  
chromatography module as it gives access to a large variety of prepacked cartridges of different  
sorbent weights, particles diameters, and phase nature (normal, reverse, grafted, chiral phases etc.).  
Sorbent weights could easily be related to HPLC column length and/or diameter since in both cases it  
45 presents an increase in active site number. It also gives the students notions regarding HPLC column  
choice. This variety of cartridges allows the study of other parameters than those commonly studied in  
typical liquid chromatography labs that are mainly focused on the influence of mobile phase  
composition on chromatographic separation. Moreover, some flash chromatography systems possess  
two different detectors, namely Ultraviolet (UV) and Evaporative Light Scattering (ELSD) Detectors,  
50 allowing students to compare the influence of the detection mode on chromatograms acquired during  
the same run.

The influence of some parameters like silica weight (i.e. column loading), particle sizes<sup>10</sup>, type of detectors<sup>11</sup> is often taught theoretically but students are rarely offered the possibility to demonstrate their influences during a lab session. This lack of practical experience is due to the fact that it is not common to own traditional HPLC columns in different lengths and silica particle sizes. Moreover, most of the time, HPLC only have one detector ranging from the most common UV-visible detector<sup>7,12,13</sup> to the most advanced (and expensive) tandem mass spectrometer.<sup>14,15</sup>

To address this lack of practical demonstration, we have designed a laboratory experiment (Figure 1) in which students are invited to experience different chromatographic parameters in an active learning style experiment thanks to an automated flash chromatography system. To the best of our knowledge, it is the first time that a flash chromatography device is used for teaching purposes.

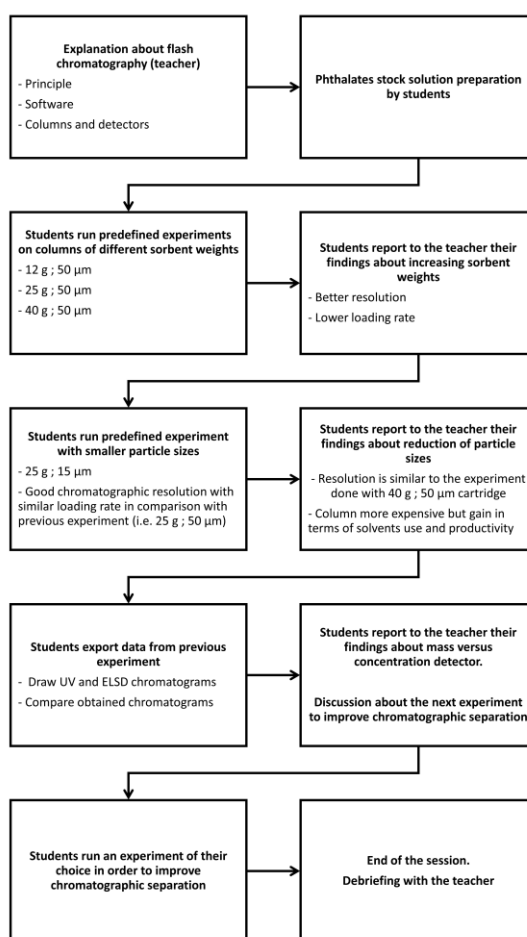

**Figure 1 : Flow diagram of the lab session**

---

## LEARNING OUTCOMES

At the end of the lab session, students are able to:

- 70
- (1) Explain the operation of a flash chromatography system
  - (2) Separate a mixture of compounds using flash chromatography system
  - (3) Explain the influence of sorbent weight as well as silica particle sizes on chromatographic resolution
  - (4) Conduct a scale change using Column Volume (CV)
- 75
- (5) Predict elution order on normal phase chromatography and compare with reverse phase HPLC
  - (6) Explain the difference between mass and concentration detector<sup>16</sup>

## EXPERIMENTAL PROCEDURE

### Apparatus description

The flash chromatography device used in this study is a classical commercial apparatus  
80 incorporating a quaternary gradient pump system, diode array detector (200 – 600 nm) allowing to monitor up to 4 wavelengths simultaneously and a full UV spectrum as well as an Evaporative Light Scattering Detector (ELSD) for universal detection.

### Column choice

Gradient optimization was performed on prepacked silica cartridges loaded with 25 g of 50  $\mu\text{m}$   
85 silica particles. Using the method of Column Volumes (CV), the gradient was transferred to cartridges with different dimensions and sorbent loadings (i.e., 12, 25 and 40 g). The same gradient was applied on the 15  $\mu\text{m}$  silica particles at a 25 g loading. Cartridges were reused several times (more than 5 times) during the same day without efficiency loss. As stated by the supplier and experienced by us, shelf-life of a column under hexane does not exceed one day.

### 90 Stock phthalate solution

Using a volumetric flask and ethyl acetate (analytical grade), students prepared 10 mL of stock solution containing 2g of each dialkyl phthalates: dimethylphthalate (DMP, >99%), diethylphthalate (DEP, 99.5%), dipropylphthalate (DPP, 98%), dibutylphthalate (DBP, 99%) This solution was used for each experiment.

---

95 **Elution conditions**

To avoid bias in comparison between column, students were asked to apply the following gradient (Table 1) expressed in Column Volume (CV) instead of time. During optimization, it was found that an adequate but not perfect resolution could be achieved in about 8 CV.

**Table 1. Elution conditions**

CV	n-hexane	Ethyl acetate
0	100	0
5	80	20
8	80	20

**HAZARDS**

100 Since the experiment requires manipulation of organic solvents and compounds as well as a potentially high-pressure device, students must wear appropriate protective equipment i.e. lab coat, goggles, and gloves. Ethyl acetate and n-hexane are highly flammable, irritant liquids, and should be manipulated with care in a well-ventilated fume hood. n-hexane is neurotoxic, suspected of damaging fertility, and toxic to aquatic life. Solvent should be disposed properly. Phthalates are potentially toxic  
105 for human reproduction<sup>17</sup> and should be manipulated with care.

**RESULTS**


The first set of experiments was performed using cartridges containing 50  $\mu\text{m}$  spherical silica in order to demonstrate the influence of column loading on the resolution (Table 2, entry 1-3). Gradient, injection volume and concentration of phthalates in solution remained constant for all experiments.  
110 The optimal flow rate for 40 g cartridge was chosen according to manufacturer data (entry 3) instead of the same 15 mL/min flow for all cartridges since no differences were evidenced between 26 and 15 mL/min (see supplementary data). All chromatograms were expressed in Column Volume (CV) as it allows to compare data recorded on cartridges of different sorbent weights. To be precise, retention times are dependent on flow rate and sorbent weight while CV are not affected by such parameters  
115 and stay relatively stable. Another advantage of the use of CV is that it makes it easy to conduct a change of scale, either up or down, in flash chromatography. From Figure 2, it can be seen that 12 g of 50  $\mu\text{m}$  silica led to a poor separation of all phthalates, particularly DBP and DPP which are clearly co-eluted. Increasing silica weight to 25 g allows us to discern the four different analytes with a resolution of only 0.76 for the first two peaks (DBP and BPP). This gain in resolution can be related to the

120 increase in number of theoretical plate which accompany increase in column length. As expected, a far better resolution (1.08) was obtained using a 40 g silica cartridge (entry 3) which results in a decrease in column loading. Resulting chromatograms were presented on Figure 2.

**Table 2. Experiment Overview**

Entry	Cartridge Conditions			CV <sub>0</sub> , mL <sup>a</sup>	Optimal Flow, mL/min	Column Loading, % <sup>b</sup>	Resolutions <sup>c</sup>		
	Sorbent Weight, g	Particle Size, μm	Internal Diameter × Length, mm				DBP/DPP	DPP/DEP	DEP/MP
1	12	50	21 × 78	20.5	15	6.6	ND <sup>d</sup>	ND <sup>d</sup>	1.16
2	25	50	21 × 127	32.5	15	3.2	0.76	1.41	1.92
3	40	50	27 × 127	52.0	26	2.0	1.08	1.85	2.35
4	25	15	21 × 127	32.5	15	3.2	1.13	2.68	2.95


<sup>a</sup>Column Volume according to manufacturer data. <sup>b</sup>Weight of analyte/weight of silica. <sup>c</sup>Resolutions were determined graphically. <sup>d</sup>Not determined due to DBP and DPP coelution.

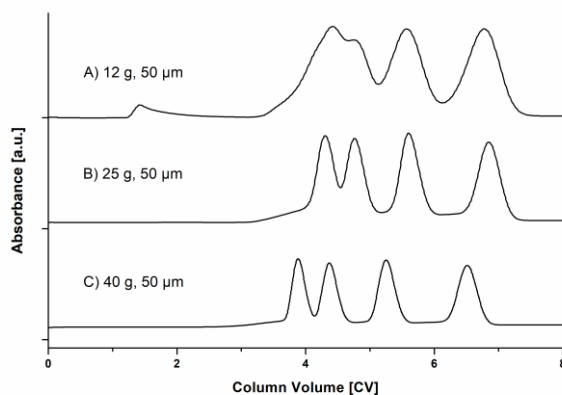

**Figure 2 : Separation of phthalates using different silica sorbent weights with 50 μm silica particles. A) 12 g, B) 25 g, and C) 40 g. Conditions: Cartridges and flow rates shown in Table 2, ambient temperature, 1 mL injection, 254 nm detection. Elution order: DBP, DPP, DET, DMP.**


125  
130 Afterwards, an experiment was run using 25g of 15 μm spherical silica in order to visualize the effect of particle sizes diminution on the chromatographic resolution. **Erreur ! Source du renvoi introuvable.**, entry 2 and 4 experiments were presented Figure 3. Students can then easily notice a better resolution (1.13) between DBP and DPP on this cartridge (25 g ,15μm silica), similar to the one obtained with 40 g cartridge (50 μm). At this point, the students were able to understand that if the CV  
135 remains fairly constant for all the compounds, the volume of solvent used is 40% lower on the

cartridge containing 25 g of silica (260 mL against 416 mL) leading to a more practical and economical separation with smaller silica particles.

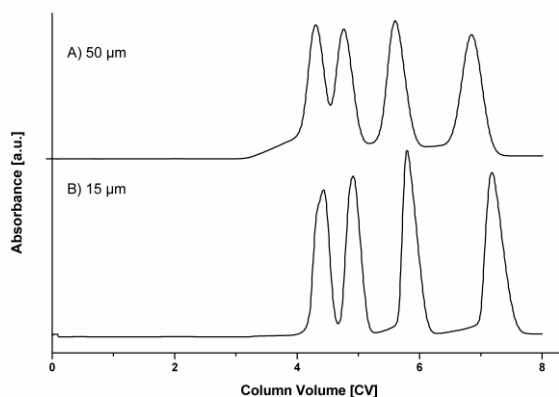

140 **Figure 3 : Separation of phthalates using different silica particles sizes at constant sorbent weight (25g). A) 50  $\mu\text{m}$ , B) 15  $\mu\text{m}$ . Conditions: 21 x 130 mm column, 15 mL/min flow, ambient temperature, 1 mL injection, 254 nm detection. Elution order: DBP, DPP, DET, DMP.**

145 A comparison of responses between the UV-detector (254 nm) and ELSD recorded during run 4 from table 2 was presented in Figure 4. It can be observed that the UV-detector provides a similar intensity of response for the four phthalates analytes while a great disparity was observed with ELSD leading to a poor signal for the smallest compound (DMP). This can be explained by the fact that in UV detection the signal is proportional to the analyte's concentration according to the Beer-Lambert's law while ELSD is sensitive to molecules sizes after solvents' evaporation, leading to, in general, intense signal for larger compounds.<sup>18</sup>

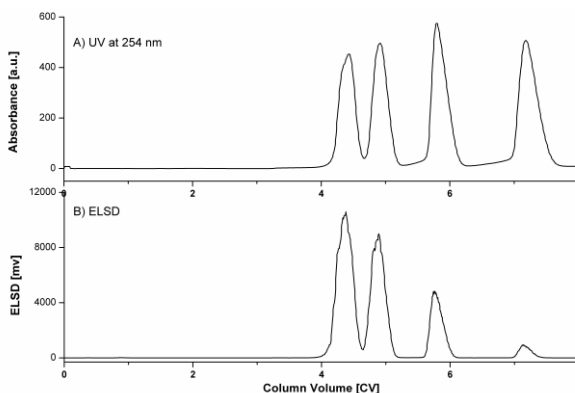

150 **Figure 4 : Response of detectors during the separation of phthalates. A) UV at 254 nm, B) ELSD. Conditions: 21 x 130 mm column, 25 g silica, 15  $\mu\text{m}$  particle sizes, 15 mL/min flow, ambient temperature, 1 mL injection, Elution order: DBP, DPP, DET, DMP.**



155

**DISCUSSION**

160

165

This 3 hours laboratory experiment was performed in autumn 2018 and autumn 2019 with master's degree students in analytical chemistry in addition to their traditional chromatography laboratory classes dealing with paraben identification by reverse-phase HPLC-UV and PAH quantification in polluted soil by GC-FID. The education and curricula of the students are diverse since some of them only know chromatography through TLC, others spend about 10 to 20 h in GC-FID, HPLC-UV and ion chromatography labs while a few of them spend a 3-month training period in the industry on a chromatograph. Those differences in student's knowledge explain why the chromatography program in our formation starts from scratch and why we choose to perform this module with graduate students and not undergraduate students. We believe that undergraduate courses dealing with chromatography theory could also benefit from this experiment.

170

The pedagogical staff decided to offer this lab session to students at the beginning of the chromatography courses in order to promote discovery learning of some chromatographic parameters that will be taught later in the semester. It enriched the students' curricula by incorporating notions in normal phase chromatography, flash chromatography and ELSD detector. Students were asked to answer the same questionnaire (see supplementary data) before and after the lab session to evaluate their initial knowledge and learning. It is important to mention that the first questionnaire was completed while students didn't yet receive instructions for the lab session in order to truly evaluate their initial knowledge. A compilation of their answers was presented in table 3.

**Table 3. Distribution of Students' Responses to the Pre-Lab Questionnaire**

Item	Questions for Response	Student Responses, <i>N</i> = 20			Remarks
		Yes, %	No, %	No Answer, %	
1	Had you ever used a liquid chromatograph?	60	40	0	
2	Had you ever heard of ELSD?	5	95	0	None could explain
3	Does the size of the silica particles in a column affect the chromatogram?	50	0	50	Only 10% could explain
4	Before the lab had you ever heard of flash chromatography?	40	60	0	Only 30% could explain
5	What type (s) of detector (s) did you use?	45 <sup>a</sup>	5 <sup>b</sup>	50	
6	What is the magnitude measured on the abscissa?	70 <sup>c</sup>	5 <sup>d</sup>	20	

<sup>a</sup>Percentage of students responding with "UV". <sup>b</sup>Percentage of students responding with "Conductimetry".

<sup>c</sup>Percentage of students responding with "Time". <sup>d</sup>Percentage of students responding with "Absorption".

---

175 Before the lab session, 60% of students never performed HPLC (Table 3 entry 1) and only 45% of  
students were able to mention UV as HPLC detector. The student who proposed the conductimetry  
was certainly confused with ion chromatography. The students probably guessed from the question  
that particle size affects the chromatogram (entry 3) but only 10% can explain how. No students were  
aware of the difference between mass and concentration detectors nor the existence or principle of the  
180 ELSD (entry 2). However, most of the students knew that the abscissa value on a chromatogram was  
the time (entry 6).

At the end of the lab session, all the students were able to mention in the questionnaire i) the type of  
detectors they employed (UV and ELSD) and if it was a mass or concentration detector ii) that particle  
sizes influence the chromatogram and how iii) that the abscissa value on a chromatogram could be  
185 expressed in time or Column Volume. Those results allow us to conclude that all the learning  
outcomes previously described were reached. A parallel was done during lecture between students'  
skill acquired during this lab session and the advantages of particle sizes diminution in UPLC as well  
as the difference between mass and concentration detectors. This parallel allows students to better  
understand those parameters that may seem abstract at first glance.

190 During the lab session students were asked to find a way to prove the elution order using common  
techniques. Among the responses, students proposed the following solutions:

- i) Inject standards compounds one by one
- ii) Add more weight of a particular phthalate in the existing solution (standard addition)
- iii) Collect pure fraction from the column and evaporate the solvent under reduced pressure to  
195 obtained pure compound that can be analyzed by infrared spectroscopy, NMR  
spectroscopy, refractive index or mass spectrometry.

Students were also asked how they can improve the chromatographic separation of the same mixture  
and to experiment with their proposal. The following solutions were proposed by students and tested  
during the lab session (see supplementary data):

- 200 i) Increase sorbent weight with 15 $\mu$ m silica particles (i.e. 40 g silica cartridge with 15  $\mu$ m  
particle sizes)

- 
- ii) Modifying the gradient (i.e. Ethyl acetate gradient in 8 CV instead of 5 CV)
  - iii) Inject less solution (i.e. 0.5 mL instead of 1 mL)

Among those solutions, the best resolutions of 1.50 and 1.73 between DBP and DPP were obtained  
205 by injecting 0.5 mL of the phthalates solution on cartridge containing 15  $\mu\text{m}$  silica particle sizes with a  
silica weight of 25 g and 40 g respectively.

## CONCLUSION

In conclusion, students were successfully introduced to automated flash chromatography, normal  
phase chromatography and the ELSD detector. They were given the opportunity to experience in an  
210 active learning style experiment the influence of column characteristics on chromatographic resolution  
as well as the different detectors responses on the same sample during the same chromatographic  
run. This lab experiment allows for a better understanding of i) column loading ii) particle sizes  
influence iii) mass vs concentration detector sensitivity.

## ASSOCIATED CONTENT

215 Supporting Information

Instructions for students; Additional information for instructors (DOCX)

## AUTHOR INFORMATION

Corresponding Author

220 \*E-mail: danjou@univ-littoral.fr

## ACKNOWLEDGMENTS

Authors would like to thanks the 2018-19 and 2019-20 promotions of the “Analyse Chimique,  
Contrôle Industriel, Environnement” master for testing the protocol, enriching the reproducibility data  
and their critical point of view on this work. Authors are also grateful to FCU Côte d'Opale and  
225 Université du Littoral Côte d'Opale (ULCO) for financial support.

- <sup>1</sup> Still, W. C.; Kahn, M.; Mitra, A. Rapid Chromatographic Technique for Preparative Separations with Moderate Resolution. *J. Org. Chem.* **1978**, 43 (14), 2923–2925.
- <sup>2</sup> Jacobson, B. M. An Inexpensive Way to Do Flash Chromatography. *J. Chem. Educ.* **1988**, 65 (5), 459.
- <sup>3</sup> Pontén, F.; Ellervik, U. Safe and Efficient Flash Chromatography Equipment for the Research/Teaching Lab. *J. Chem. Educ.* **2001**, 78 (3), 363.
- <sup>4</sup> Butler, J. D.; Choung, W.; Kurth, M. J. Flash Chromatography: A Novel Pressurization Apparatus. *J. Chem. Educ.* **2010**, 87 (11), 1265–1265.
- <sup>5</sup> Buchi English-language homepage. <https://www.buchi.com/> (accessed Jan 31, 2020).
- <sup>6</sup> Biotage English-language homepage. <https://www.biotage.com> (accessed Jan 31, 2020).
- <sup>7</sup> Interchim English-language homepage. <http://www.interchim.com> (accessed Jan 31, 2020).
- <sup>8</sup> Naumiec, G. R.; Del Padre, A. N.; Hooper, M. M.; St. Germaine, A.; DeBoef, B. A Modern Apparatus for Performing Flash Chromatography: An Experiment for the Organic Laboratory. *J. Chem. Educ.*, **2013** 90 (3), 376-378
- <sup>9</sup> Husgafvel- Pursiainen, K.; Kannio, A.; Oksa, P.; Suitiala, T.; Koskinen, H.; Partanen, R.; Hemminki, K.; Smith, S.; Rosenstock- Leib, R.; Brandt- Rauf, P. W. Mutations, tissue accumulations, and serum levels of p53 in patients with occupational cancers from asbestos and silica exposure. *Environ. Mol. Mutagen.*, **1997**, 30, 224-230
- <sup>10</sup> Chester, T. L. Recent Developments in High-Performance Liquid Chromatography Stationary Phases *Anal. Chem.*, **2013**, 85, 579-589
- <sup>11</sup> Swartz, M. HPLC detectors: a brief review. *J. Liq. Chromatogr. Relat. Technol.*, **2010**, 33:9-12, 1130-1150
- <sup>12</sup> Miranda, B.; Lawton, N. M.; Tachibana, S. R.; Swartz, N. A.; Hall, W. P. Titration and HPLC Characterization of Kombucha Fermentation: A Laboratory Experiment in Food Analysis. *J. Chem. Educ.*, **2016**, 93 (10), 1770-1775
- <sup>13</sup> Purcell, S. C.; Pande, P.; Lin, Y.; Rivera, E. J.; Paw U, L.; Smallwood, L. M.; Kerstiens, G. A.; Armstrong, L. B.; Robak, M.; Baranger, A. M.; Douskey, M. C. Extraction and Antibacterial Properties Of Thyme Leaf Extracts: Authentic Practice of Green Chemistry. *J. Chem. Educ.*, **2016**, 93 (8), 1422-1427
- <sup>14</sup> Parker, P. D.; Beers, B.; Vergne, M. J. What Is in Your Wallet? Quantitation of Drugs of Abuse on Paper Currency with a Rapid LC–MS/MS Method. *J. Chem. Educ.*, **2017**, 94 (10), 1522-1526
- <sup>15</sup> Betts, T. A.; Palkendo, J. A. Teaching Undergraduates LC–MS/MS Theory and Operation via Multiple Reaction Monitoring (MRM) Method Development. *J. Chem. Educ.*, **2018**, 95 (6), 1035-1039
- <sup>16</sup> Urban, P. L. Clarifying Misconceptions about Mass and Concentration Sensitivity. *J. Chem. Educ.*, **2016**, 93 (6), 984-987
- <sup>17</sup> Mariana, M.; Feiteiro, J.; Verde, I.; Cairrao, E. The effects of phthalates in the cardiovascular and reproductive systems: A review. *Environ. Int.*, **2016**, 94, 758-776
- <sup>18</sup> Waters Corporation. 2424 Evaporative Light Scattering Detector Operator's Guide. <https://www.waters.com/webassets/cms/support/docs/71500121802rb.pdf> (accessed Jan 31, 2020).