



Narrowband Modulation Two-Dimensional Mass Spectrometry and Label-Free Relative Quantification of Histone Peptides

Matthias Halper, Marc-André Delsuc, Kathrin Breuker, Maria A van Agthoven

► To cite this version:

Matthias Halper, Marc-André Delsuc, Kathrin Breuker, Maria A van Agthoven. Narrowband Modulation Two-Dimensional Mass Spectrometry and Label-Free Relative Quantification of Histone Peptides. *Analytical Chemistry*, 2020, 92 (20), pp.13945-13952. 10.1021/acs.analchem.0c02843 . hal-03090931

HAL Id: hal-03090931

<https://hal.science/hal-03090931>

Submitted on 30 Dec 2020

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

Narrowband Modulation Two-dimensional Mass Spectrometry and Label-free Relative Quantification of Histone Peptides

Matthias Halper¹, Marc-André Delsuc^{2,3}, Kathrin Breuker¹, Maria A. van Agthoven^{1#}

¹ Institute for Organic Chemistry, University of Innsbruck, Innrain 80/82, 6020 Innsbruck, Austria

² Institut de Génétique et de Biologie Moléculaire et Cellulaire, INSERM U596, UMR 7104, Université de Strasbourg, 1 rue Laurent Fries, 67404 Illkirch-Graffenstaden, France

³ CASC4DE, Pôle API, 300 Bd. Sébastien Grant, 67400 Illkirch-Graffenstaden, France

KEYWORDS Two-dimensional mass spectrometry, histones, post-translational modifications, label-free relative quantification, FT-ICR mass spectrometry

ABSTRACT: Two-dimensional mass spectrometry (2D MS) on an FT-ICR mass analyzer allows for tandem mass spectrometry without requiring ion isolation. In the ICR cell, precursor ion radii are modulated before fragmentation, which results in modulation of the abundance of their fragments. The resulting 2D mass spectrum enables correlation between precursor and fragment ions. In standard broadband 2D MS, the range of precursor ion cyclotron frequencies is determined by the lowest mass-to-charge (m/z) ratio to be fragmented in the 2D MS experiment, which leads to precursor ion m/z ranges that are much wider than necessary, thereby limiting the resolving power for precursor ions and the accuracy of the correlation between precursor and fragment ions. We present narrowband modulation 2D MS, which increases the precursor ion resolving power by reducing the precursor ion m/z range, with the aim of resolving the fragment ion patterns of overlapping isotopic distributions. In this proof-of-concept study, we compare broadband and narrowband modulation 2D mass spectra of an equimolar mixture of histone peptide isoforms. In narrowband modulation 2D MS, we were able to separate the fragment ion patterns of all ^{13}C isotopes of the different histone peptide forms. We further demonstrate the potential of narrowband 2D MS for label-free quantification of peptides.

INTRODUCTION

Histones are 11-21 kDa proteins that make up the chief components of chromatin, the complex around which DNA is wound in the nucleus.¹ Histones frequently carry numerous post-translational modifications (PTM), which have been found to play a central role in gene regulation, DNA repair, chromosome condensation, and spermatogenesis. A functional understanding of histone modifications and the proposed “histone code” requires far more experimental data than are available so far.²

Mass spectrometry (MS) is ideally suited to address this challenge as it can directly detect all mass-altering modifications and does not require laborious biochemical techniques for analysis. Specifically, MS can provide mass values of both the biomolecules under study and their fragments from dissociation in tandem mass spectrometry (MS/MS) experiments, for which pre-fractionation by chromatography is typically required.^{3, 4} However, it is often challenging - if not impossible - to fractionate very similar compounds, e.g., post-translationally modified proteins that differ only in their modification patterns. Likewise, standard MS/MS experiments, in which each ionized compound (precursor ion) in a sample is isolated before dissociation and detection of its fragment ions in the mass spectrometer, can be limited by overlapping isotopic distributions of different peptide or proteins forms. Although the label-free, direct localization and relative quantitation of histone⁵⁻⁷ and ribonucleic acid (RNA)⁸ modifications by MS/MS have been demonstrated, a more general approach that does not rely on precursor

ion isolation would significantly advance all fields of PTM research.

Two-dimensional mass spectrometry (2D MS) is a method for tandem mass spectrometry that enables the correlation between precursor and fragment ions without ion isolation.⁸⁻¹⁰ 2D MS is a useful technique for samples that are difficult to separate through chromatographic methods, like protein isoforms with relatively small modifications (e.g. methylations). Another advantage of 2D MS is the fact that all analytes are subjected to the same experimental conditions, which makes 2D MS a good candidate to perform quantification studies.

2D MS can be performed on FT-ICR or linear ion trap mass spectrometers, and has been successfully applied to the analysis of small molecules, polymers, and in bottom-up and top-down proteomics.¹¹⁻²¹ In most studies using 2D MS, the analytes were within an m/z range of several hundred a.m.u. In the case of modified histones, the m/z range is typically less than 10 a.m.u., depending on the charge state of the protein.

High-resolution mass spectrometers, i.e., the Fourier transform ion cyclotron resonance mass spectrometer (FT-ICR MS) and the Orbitrap are used to resolve complex fragmentation patterns of, for example, those of large peptides or proteins, with high mass accuracy.²²⁻²⁴ Efficient two-dimensional mass spectrometry in an FT-ICR MS requires a gas-free fragmentation method such as infrared multiphoton dissociation (IRMPD), which preferentially fragments the most labile bonds including those of PTMs of peptides and proteins, or electron-capture

dissociation (ECD), which preferentially fragments the C- α backbone bonds of peptides and proteins.²⁵⁻²⁷ These characteristics make ECD the method of choice for PTM analysis of peptides and proteins in 2D MS.²⁸

In standard broadband 2D MS, the cyclotron frequency range for precursor ions is determined by the lowest m/z ions to be fragmented in the 2D MS experiment, which leads to m/z ranges that are much wider than necessary, thereby limiting the resolving power of the precursor ions (corresponding to the isolation m/z window in standard tandem mass spectrometry). In one-dimensional FT-ICR mass spectrometry, heterodyne or “narrowband” techniques have been developed to enable ultrahigh resolution.²⁹⁻³¹ We present narrowband modulation 2D MS, which is conceptually different from one-dimensional narrowband MS and does not use heterodyne detection. Narrowband modulation 2D MS increases the precursor ion resolving power with the aim of differentiating the fragmentation patterns of overlapping isotopic distributions (e.g. peptides and proteins with and without disulphide bonds).

In this study, we compare the analysis of histone peptide isoforms with broadband and narrowband modulation 2D MS. We further discuss the use of 2D MS for the localization and assignment of PTMs, as well as their relative quantification. In order to simplify the manuscript, narrowband modulation 2D MS will be referred to as “narrowband 2D MS”, although this is not meant to indicate an equivalence to one-dimensional narrowband MS techniques.

EXPERIMENTAL METHODS

Sample Preparation

MiliQ® water (Merck Millipore, Darmstadt), methanol and acetic acid (VWR, Vienna, Austria) were used in all experiments. As model peptides, C-terminal GK-biotinylated histone H3 sequences (amino acid residues 21 to 44) without and with native modifications (mono-, di- and tri-methylation) at K27 were purchased (AnaSpec, Fremont, CA, USA) with a purity >95%, in which K27 of the full-length histone corresponds to K7 of the model peptides.

Peptides were desalted using MWCO 2000 Vivaspın centrifugal concentrators (Sartorius, Göttingen, Germany) at 7900 rcf, 6× ammonium acetate 100 mM (Sigma, Vienna, Austria), 6× H₂O. Peptide concentration was determined by UV-absorption at 280 nm using an Implen Nano PhotometerTM (Implen, München, Germany). For electrospray ionization (ESI), equimolar mixtures (0.1 μ M each) of unmodified (referred to as “wild type” in the text), mono-, di-, and trimethylated peptides (referred to as “K7 1m”, “K7 2m”, and “K7 3m”, respectively) in 50:50 H₂O/CH₃OH and 1% vol. CH₃COOH, pH ~3.0 were prepared from 100 μ M stock solutions of each peptide in H₂O.

Instrument Parameters


2D mass spectra were acquired on a 7 T Apex Ultra FT-ICR mass spectrometer (Bruker Daltonik, GmbH, Bremen, Germany) with an electrospray ion source operated in positive mode and direct injection at a flow rate of 90 μ L/h. Ions were accumulated for 0.1 s in the first hexapole and 0.1 s in the second hexapole.³² In the quadrupole, ions were isolated at m/z 490 with an isolation window of m/z 30 and transferred to the Infinity ICR cell through a series of focusing lenses.³³

The experimental script for the acquisition of 2D mass spectra is shown in Scheme 1. All excitation pulses (P_1 , P_2 , and P_3)

have a frequency range of 357,144.31 – 71,669.99 Hz (corresponding to m/z 301.08 – 1500).

The frequency list was composed of 458 frequencies with a 624.71 Hz decrement. For the encoding pulses P_1 and P_2 , the power level (corresponding to the amplitude of the pulses) was set at 106 V_{pp} with a pulse length of 1.0 μ s per frequency. For the excitation pulse P_3 , the power level was set at 80 V_{pp} with a

Scheme 1. Experimental script of two-dimensional mass spectrometry.

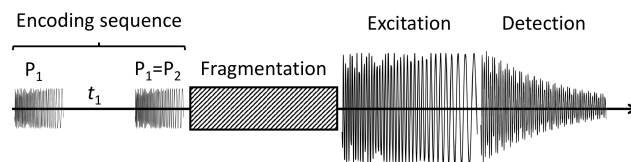

pulse length of 20 μ s per frequency (default settings in the Apex Control software, Bruker Daltonik, GmbH, Bremen, Germany).³⁴ The fragmentation mode used was electron capture dissociation, using a hollow cathode.³⁵ The heater was at 1.3 A, the lens at 20 V, the ECD bias at 2.0 V, and the irradiation pulse at 0.03 s. The transients were recorded with 1 Mwords (16 bits) for 1.4680 s, with a 714,288.62 Hz sampling frequency.

For the broadband 2D mass spectrum, the encoding period t_1 was incremented 4096 times with a 2.0 μ s increment, leading to a 250 kHz Nyquist frequency (mass range: m/z 334.439 – 1500 with the high m/z limit set by the frequency range of P_1). The longest encoding period was 8.192 ms. For the narrowband 2D mass spectrum, the encoding period t_1 was incremented 4096 times with a 8.0 μ s increment, leading to a 62.5 kHz Nyquist frequency (mass range: m/z 415.00 – 546.825). The longest encoding period was 32.768 ms. For neither 2D mass spectrum was the detection heterodyne.

Data Processing

The 2D mass spectra were processed and visualized using the SPIKE software developed independently by the University of Strasbourg and CASC4DE (Illkirch-Graffenstaden, France) in 64-bit Python 3.7 programming language on an open source platform distributed by the Python Software Foundation (Beaverton, OR, SA).³⁶ Processed datafiles were saved using the HDF5 file format. Both 2D mass spectra were digitally demodulated with a demodulation frequency of 71,669.99 Hz and denoised using the Support Selection for Noise Elimination (SANE) algorithm with a SANE rank of 50.^{9, 37-39}

For the broadband 2D mass spectrum, cyclotron frequencies in the vertical precursor ion dimension were calculated with the following equation:

$$f_{ICR} = f + f_0 \quad (\text{eq. 1})$$

in which f_{ICR} is the reduced cyclotron frequency of the precursor ion, f is the measured modulation frequency of the precursor ion, and f_0 is the lowest frequency in the excitation pulses (71,669.99 Hz). For the narrowband 2D mass spectrum, cyclotron frequencies in the vertical precursor ion dimension were calculated with the following equation:


$$f_{ICR} = f + f_0 + 2f_N \quad (\text{eq. 2})$$

in which f_N is the Nyquist frequency in the vertical precursor ion dimension (62,500.00 Hz).

Each 2D mass spectrum was internally calibrated in the horizontal fragment ion m/z dimension using a quadratic frequency-to-mass conversion.⁴⁰ For each precursor ion species, 4-

5 fragment ion spectra with the highest fragment ion intensities were added before peak-picking.

Scheme 2. Schematics of two-dimensional mass spectra obtained in broadband (left) and narrowband mode (right).

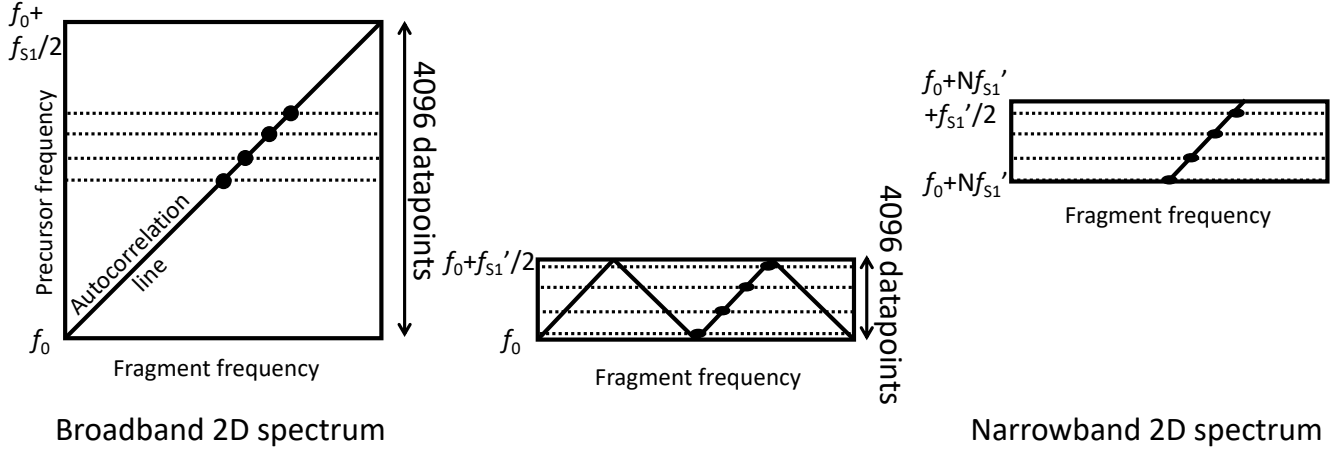

THEORY OF NARROWBAND TWO-DIMENSIONAL MASS SPECTROMETRY

Scheme 1 shows the experimental script for 2D FT-ICR MS. A detailed description of the mechanisms of the experiment can be found in previous articles.^{9, 41, 42} Scheme 2 (left) shows a schematic of a broadband 2D spectrum in the frequency domain. The modulation frequency in the vertical precursor frequency is $f_{ICR}-f_0$, in which f_{ICR} is the cyclotron frequency of the precursor ion and f_0 is the lowest frequency in the P_1 and P_2 pulses (see scheme 1). The highest modulation frequency f_{max1} in the vertical dimension is determined by the encoding period (t_1 in Scheme 1) increment Δt_1 :

$$f_{max1} = \frac{1}{2} f_{s1} = \frac{1}{2\Delta t_1} \quad (\text{eq. 3})$$

In which f_{s1} is the sampling frequency in the vertical precursor dimension. As a result, the range of precursor cyclotron frequencies in the 2D spectrum is $[f_0, f_{s1}/2]$.

According to the Shannon-Nyquist theorem, periodically sampling a signal at a sampling frequency of f_s produces an aliasing of the signal frequencies such that all frequencies f_N separated by f_s of a multiple of f_s produce the same measurement f :⁴³

$$f_N = \pm(f + Nf_s), N \in \mathbb{Z} \quad (\text{eq. 4})$$

The \pm sign is due to the cosine in the modulation equation for 2D MS.^{9, 41} Signals in 2D MS experiments that have frequencies higher than $f_{s1}/2$ can therefore be represented in aliased form in the 2D spectrum.

In Scheme 2 (left), the autocorrelation line has modulations signals up to $f_{s1}/2$ (precursor ion cyclotron frequencies are between f_0 and $f_0 + f_{s1}/2$). If the sampling frequency is reduced in the experiment by increasing the encoding period increment Δt_1 (Scheme 2 center), then the autocorrelation line is aliased and can be folded back into the spectrum multiple times and appears to zig-zag.

The sampling frequency f_{s1}' and the lowest excitation frequency f_0 can be chosen so that the modulation frequency range of interest is folded back into the spectrum an even number of times. The precursor ion cyclotron frequencies can be recovered by offsetting all frequencies in the precursor ion dimension by N

f_{s1}' , where N is the number of times that the signal of interest has been folded back into the spectrum (Scheme 2 right). Subsequent frequency-to-mass conversion yields the narrowband 2D mass spectrum.

In cases where only a small m/z range is necessary in the vertical precursor ion dimension, such as a single charge state of a protein and its isoforms, narrowband 2D MS significantly reduces the vertical precursor m/z range, which significantly increases the resolving power of each peak in the vertical precursor ion dimension and enables more accurate precursor-fragment correlation.

The total acquisition time of a 2D MS experiment is proportional to the number of t_1 increments. As a result, narrowband 2D MS can either be used to increase the vertical resolving power by keeping a constant number of data points, or decreasing experiment times and sample consumption and retaining the same resolving power.

RESULTS AND DISCUSSION

Figure 1 shows the full broadband (Fig. 1a) and narrowband (Fig. 1b) 2D mass spectra of a mixture of the 4 isoforms of the H3 histone peptide.

In both 2D mass spectra, fragment ion m/z ratios are plotted horizontally. Precursor ion m/z ratios, calculated after conversion from measured frequencies to reduced cyclotron frequencies with equation 1 (broadband) and equation 2 (narrowband), are plotted vertically. In both 2D mass spectra, the autocorrelation line has been extracted and shows the isotopic distribution of the $[M+6H]^{6+}$ ions of the 4 isoforms, evenly spaced by m/z 2.335942 (corresponding to 14.015652 Da for CH_2 divided by $z=6$). Because of the quadrupole isolation, these species are the only ones detected on the autocorrelation line. The resolving power on the autocorrelation line is approximately 147,000 at m/z 490 in both 2D mass spectra.

In both 2D mass spectra, the ECD fragmentation patterns of the peptide isoforms are represented in the horizontal fragment ion spectra. The sequence coverage of the unmodified (wild type) peptide is shown in Fig. 1a and 1b for the respective 2D mass spectra. In both cases, the sequence coverage is nearly complete. The peak assignments for all fragment ion spectra can be found in the Supporting Information for the broadband 2D

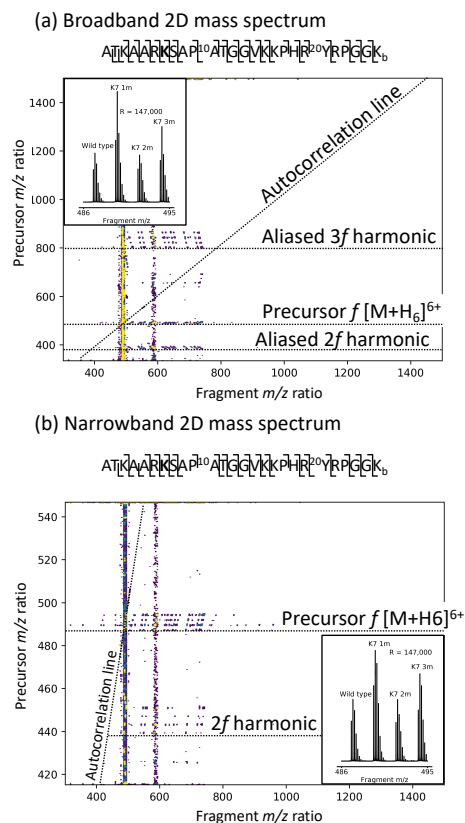


Figure 1. (a) Broadband 2D ECD mass spectrum of histone peptide isoforms (wild type, ^7K monomethylated or K7 1m, ^7K dimethylated or K7 2m, ^7K trimethylated or K7 3m). Inset: extracted autocorrelation line, sequence coverage of wild type histone peptide. (b) Narrowband 2D ECD mass spectrum of histone peptide isoforms (wild type, ^7K monomethylated, ^7K dimethylated, ^7K trimethylated). Inset: extracted autocorrelation line, sequence coverage of wild type histone peptide.

mass spectrum (Tables S1 to S4) and the narrowband 2D mass spectrum (Tables S5 to S8).

In both 2D mass spectra, the ECD fragmentation patterns of the peptide isoforms are represented in the horizontal fragment ion spectra. The sequence coverage of the unmodified (wild type) peptide is shown in Fig. 1a and 1b for the respective 2D mass spectra. In both cases, the sequence coverage is nearly complete. The peak assignments for all fragment ion spectra can be found in the Supporting Information for the broadband 2D mass spectrum (Tables S1 to S4) and the narrowband 2D mass spectrum (Tables S5 to S8). The sequence coverage for all four isoforms can be found in Figure S1 of the Supporting Information for the broadband 2D mass spectrum and Figure S2 for the narrowband 2D mass spectrum. For all isoforms, in both 2D mass spectra, the sequence coverage is nearly complete (over 20 out of 25 cleavages).

In the broadband 2D mass spectrum, the vertical range is m/z 334.439 – 1500. Only 0.77% of the vertical m/z range contains analytically useful information (m/z 486 – 495 for precursor ions), because it is determined by the Nyquist frequency, which has to be higher than the cyclotron frequency of the lowest m/z ratio of interest. In the narrowband 2D mass spectrum, the vertical mass range is m/z 415.00 – 546.825, which means that 6.8% contains useful analytical information. In the narrowband 2D mass spectrum, the vertical m/z range is still determined by

the Nyquist frequency, but the foldover effect is used to reduce it.

In both 2D mass spectra, harmonics in the vertical are clearly visible. Vertical harmonics, as well as residual scintillation noise, are clearly visible in the 2D mass spectra in Fig. 1 because of the low relative abundance of fragment ion peaks compared to precursor ion peaks in ECD.^{38,44} Vertical harmonics in 2D mass spectra arise from the non-sinusoidal modulation of signals during the 2D MS experiment.^{34,41,42,45} Their relative intensity can be minimized by adjusting the modulation zone in the ICR cell to the fragmentation zone during the encoding sequence (Scheme 1) but they cannot be eliminated.³⁴ A previous study has shown that a careful examination of the harmonics in precursor ion spectra can uncover additional peaks for those fragment ions whose fragmentation zones differ significantly from a Gaussian profile (e.g. in cases of secondary fragmentation or combined fragmentation methods, such as IR-ECD).⁴⁶

In the narrowband 2D mass spectrum (Fig. 1b), the $2f$ harmonic is visible (the $3f$ harmonic has too low a signal-to-noise ratio to be visible). The absence of precursor ion peaks on the autocorrelation line clearly marks the fragment ion spectra as harmonics. In the broadband 2D mass spectrum (Fig. 1a), the aliased $2f$ and $3f$ harmonics are visible. The foldover effect is marked by the fact that, in both cases, the slopes of the dissociation lines are negative. The signal from vertical harmonics can be used, in a similar way to the signal from harmonics in one-dimensional FT-ICR mass spectra, to increase the resolving power in the vertical dimension, and therefore to improve the accuracy of the precursor-fragment correlation.^{47–49} Although the smaller mass range in the precursor ion dimension of narrowband 2D mass spectra may lead to a decrease in harmonic signals, they may be folded back into the spectrum by aliasing.

Figure 2 shows vertical precursor ion spectra extracted from both the broadband and the narrowband 2D mass spectrum. Fig. 2a shows the vertical precursor ion spectrum of the monoisotopic c_5 fragment at m/z 460.2879. Fragment c_5 does not contain any modification, and therefore the precursor ion spectrum shows peaks for all 4 isoforms. In the vertical precursor ion spectrum extracted from the broadband 2D mass spectrum, the isotopic peaks of the four precursor ions of c_5 are not resolved (black circles show the theoretical isotopic peaks in the vertical dimension). In contrast, the isotopic distribution of the precursor ions is resolved in the precursor ion spectrum extracted from the narrowband 2D mass spectrum, albeit not to the baseline.

The vertical resolving power in both spectra can be calculated from the extracted precursor ion spectra shown in Fig. 2b, which correspond to the monoisotopic z_{24}^{4+} fragment of the wild type peptide at m/z 682.8849 (to help compare between the two peaks, the spectra have been normalized, which makes the noise level of the spectrum extracted from the narrowband 2D mass spectrum appear higher). The z_{24}^{4+} fragment contains the K7 residue on which all modifications (methylation, dimethylation, trimethylation) are located, so the precursor ion spectrum only shows a peak for the wild type peptide. In addition, the z_{24}^{4+} fragment contains 24 out of 26 residues. Any ^{13}C isotope in the precursor ion has a probability close to 100% to remain in the z_{24}^{4+} fragment instead of its complement c_2 .^{9,50} Therefore, the monoisotopic z_{24}^{4+} fragment is the product of the monoisotopic peak of the wild type precursor with negligible contribution from any ^{13}C isotopes of the precursor ion, and the extracted precursor ion spectra in Fig. 2b can be used to accurately determine the vertical resolving power and the accuracy of the

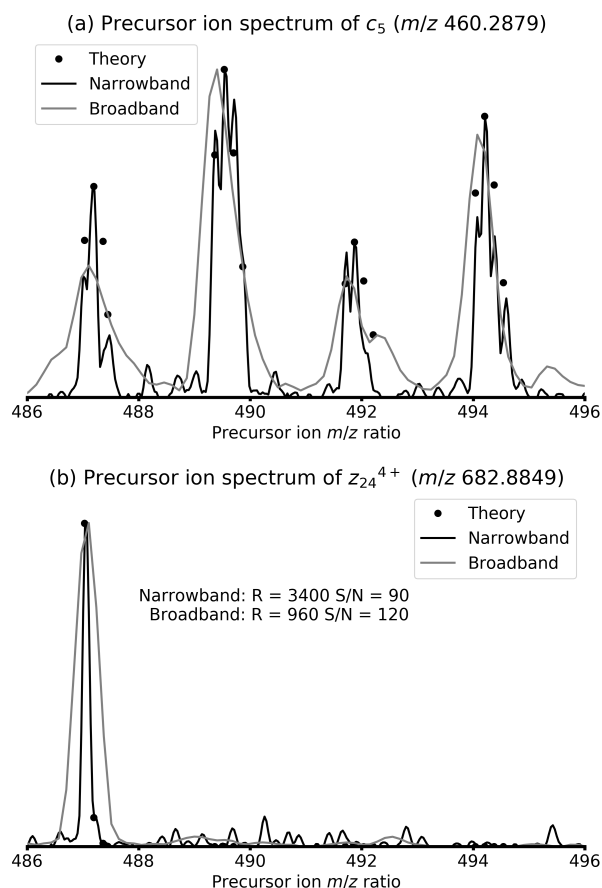


Figure 2. (a) Normalized precursor ion spectrum of c_5 extracted from the broadband and the narrowband 2D mass spectra (black dots: theoretical isotopic distribution). (b) Normalized precursor ion spectrum of z_{24}^{4+} extracted from the broadband and the narrowband 2D mass spectra at the first harmonic (black dots: theoretical isotopic distribution).

precursor-fragment correlation. The isotopic distributions of fragment ions in 2D mass spectra are further illustrated in Fig. S4 in the Supporting Information.

As expected, the vertical resolving power in the narrowband 2D mass spectrum (3,400 at m/z 487) is almost 4 times the vertical resolving power in the broadband 2D mass spectrum (960 at m/z 487). This result is consistent with the fact that the frequency range in the broadband 2D mass spectrum is 4 times the frequency range in the narrowband 2D mass spectrum with an equal number of datapoints in the vertical dimension.

Figure 3 shows the vertical precursor ion spectra of all four isotopes of the z_{24}^{4+} fragment of the $[M+6H]^{6+}$ charge state of the wild type histone peptide. Figure 3a shows the first harmonic peak (i.e. corresponding to the modulation frequency of the precursor ion radius) and figure 3b shows its $2f$ harmonic peak (or second harmonic) with an adapted frequency-to-mass conversion as extracted from the broadband 2D mass spectrum. Figure 3c shows the first harmonic peak and figure 3d shows its second harmonic peak as extracted from the narrowband 2D mass spectrum.

In Fig. 3a (first harmonic of the broadband 2D mass spectrum), the isotopic distribution is not resolved ($R = 960$ for the ^{12}C isotopic peak). In Fig. 3b, the use of the second harmonic enables an increase in resolving power ($R = 2,200$), but the isotopic distribution is still not resolved. The signal-to-noise ratio

is slightly lower for the second harmonic than in the first harmonic. The third harmonic is shown in the Supporting Information Fig. S5.

In Fig. 3c, while the isotopic distribution of the $[M+6H]^{6+}$ precursor is resolved ($R = 3,400$ for the ^{12}C isotopic peak), the peaks are not baseline-resolved. In Fig. 3d, the isotopic distribution is resolved ($R = 7,700$ for the ^{12}C isotopic peak). In the first harmonic, the signal-to-noise ratio of the $1\times^{13}\text{C}$ isotopic peak was calculated to be 110, but in the second harmonic, it was calculated to be 50. The horizontal fragment ion spectra for the second harmonic of all four isoforms (visible between vertical m/z 435-450 in Fig. 1b) were peak-picked. The peak assignments are shown in the Supporting Information (Tables S9-S12).

The sequence coverage of all 4 isoforms in the second harmonic (Supporting Information Figure S3) was approximately 50%, compared to approximately 90% in the first harmonic. Unlike in one-dimensional FT-ICR mass spectrometry, where the relative intensity of $2f$ - $6f$ harmonics can be maximized by adapting the design of ICR cells,⁴⁷⁻⁴⁹ increasing the relative intensity of higher harmonics in 2D FT-ICR MS comes at the expense of overall signal-to-noise ratio and thus analytical information.³⁴ However, in 2D MS, higher harmonics in the vertical dimension can be useful to establish precursor-fragment correlations when first harmonic fragment ion peaks have insufficient resolving power.

Fig. 3 also highlights the necessity to add up fragment ion spectra in the narrowband 2D mass spectrum to obtain complete isotopic distributions. Further illustration can be found in Fig. S6 in the Supporting Information.

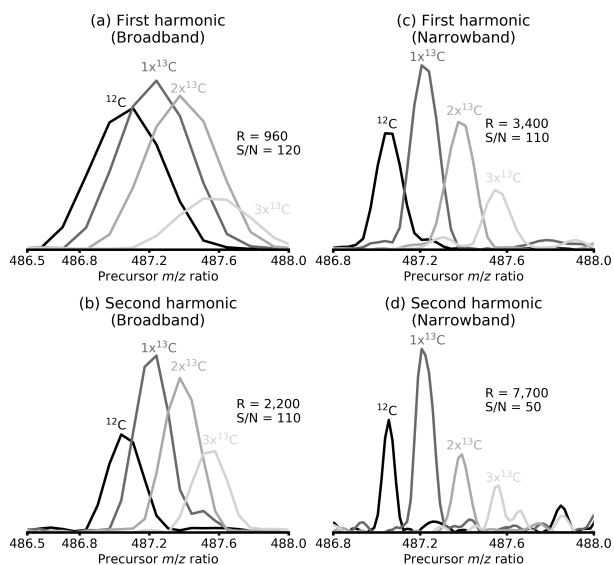


Figure 3. Precursor ion spectra of the isotopes of the z_{24}^{4+} fragment of the wild type histone peptide, with frequency-to-mass conversion extracted from the broadband 2D mass spectrum for (a) the first harmonic peak, and (b) the second harmonic peak and extracted from the narrowband 2D mass spectrum for (c) the first harmonic peak, and (d) the second harmonic peak. The resolving power is measured on the ^{12}C isotopic peak, and the signal-to-noise ratio is measured for the $1\times^{13}\text{C}$ isotopic.

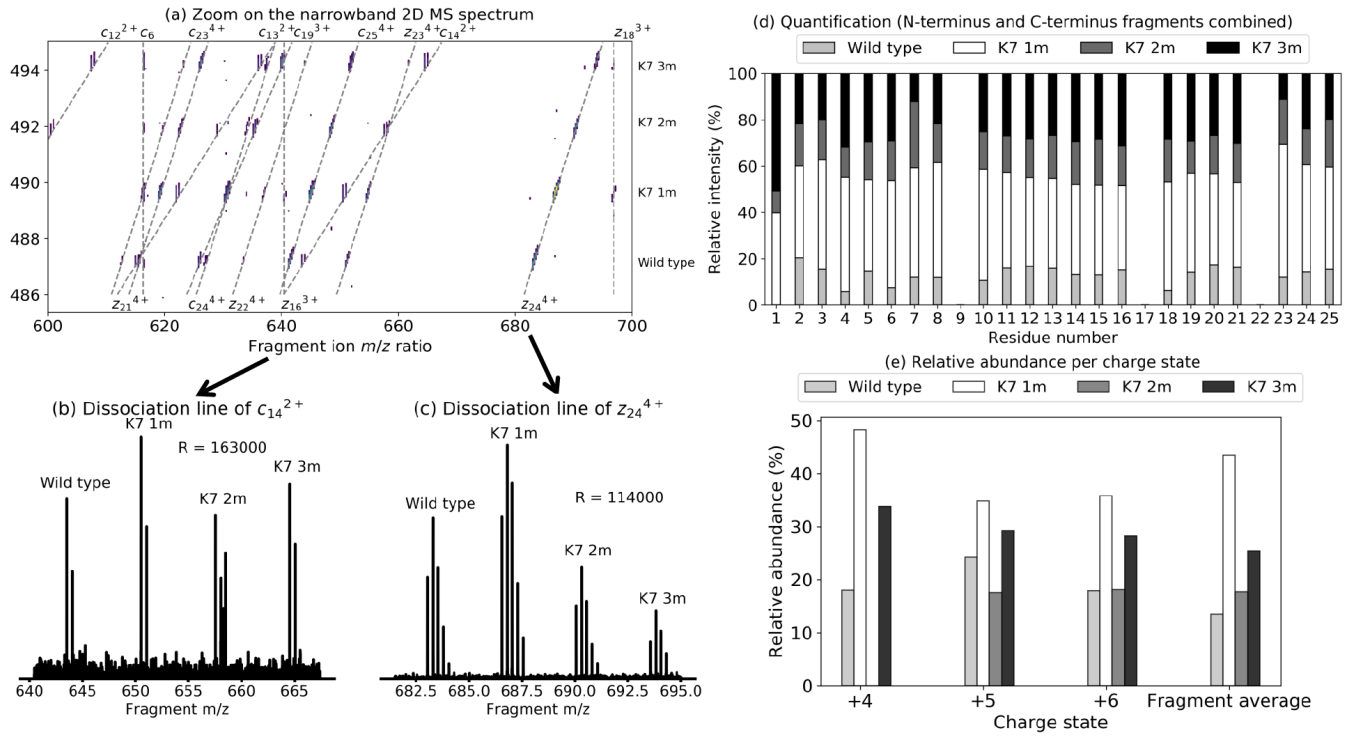


Figure 4. (a) Zoom on the narrowband 2D mass spectrum of the histone peptide mixture with various dissociation lines highlighted. (b) Extracted dissociation line of the c_{14}^{2+} fragment. (c) Extracted dissociation line of the z_{24}^{4+} fragment. (d) Label-free relative quantification of the four histone peptide isoforms using fragment abundances for each residue (N-terminus and C-terminus combined). (e) Label-free relative quantification of the four histone peptide isoforms using the relative abundances of each charge state (6+ precursor ion, 5+ and 4+ charge-reduced states) and average of the quantification using fragment ion abundances.

Figure 4 shows how the narrowband 2D mass spectrum enables the sequencing of the four peptide isoforms, the identification and location of the modifications, and their relative quantification. Figure 4a shows a zoom-in of the narrowband 2D mass spectrum between m/z 600–700. For each peptide isoform, horizontal fragment ion spectra were extracted for each isotope and added up before peak-picking and peak assignment (Tables S5–S8 and Fig. S2 in the Supporting Information).

Fig. 4a illustrates how to locate peptide modifications. In the 2D mass spectrum, precursor ions of different m/z with fragments of identical m/z (i.e. fragments not containing the modified lysine) can all be found on the same vertical precursor ion spectrum. In Fig. 4a, c_6 (m/z 616.3888), z_{16}^{3+} (m/z 640.3481), and z_{18}^{3+} (m/z 696.3790) all form vertical precursor ion spectra. Since the methylations are situated on residue number 7 counting from the N-terminus or residue number 19 counting from the C-terminus, these vertical precursor ion spectra enable the accurate localization of the modifications.

The complementary ions for fragments with identical m/z ratios can be found on diagonal dissociation lines following equation 9:^{9, 15, 46, 51}

$$(m/z)_p = \frac{p}{n}(m/z)_f + \frac{\Delta m}{n} \quad (\text{eq. 9})$$

in which $(m/z)_p$ is the m/z ratio of the precursor ion, $(m/z)_f$ the m/z ratio of the fragment, n the charge state of the precursor, p the charge state of the fragment, and Δm the mass difference between the precursor and the fragment. In Fig. 4a, c_{12}^{2+} , c_{13}^{2+} , and c_{14}^{2+} are observed on dissociation lines with a slope of $2/6=0.33$, c_{19}^{3+} on a dissociation line with a slope of $3/6=0.5$, and finally z_{21}^{4+} , z_{22}^{4+} , z_{23}^{4+} , z_{24}^{4+} , c_{23}^{4+} , c_{24}^{4+} , and c_{25}^{4+} on dissociation lines with a slope of $4/6=0.67$ (linear regressions

available in Table S13 in the Supporting Information). The use of vertical precursor ion spectra and diagonal dissociation lines can therefore be used to locate the sites of peptide modifications.

Fig. 4b and 4c show the extracted dissociation lines of c_{14}^{2+} (Fig. 4b) and z_{24}^{4+} (Fig. 4c). The resolving power on both dissociation lines is sufficient to calculate exact mass differences between the isoforms for correct assignments of the peptide modifications. In Fig. 4b, the modifications are measured to have masses of 14.0155 Da, 28.0313 Da, and 42.0465 Da, corresponding to mono-, di-, and trimethylation, respectively. The unusual pattern for the K7 2m isotopic distribution results from overlap of the c_{14}^{2+} and z_{23}^{4+} fragment ion distributions (see Fig. 4a). Figs. 4a–c therefore show that a visual analysis of narrowband 2D mass spectra can be used for peptide sequencing, and to identify and locate the modifications.

Figure 4d and 4e show different methods to quantify the four peptide isoforms. In label-free relative quantification studies with one-dimensional MS/MS, relative intensities of ions with given mass differences are compared to locate modifications and quantify isoforms on one graph.^{7, 52, 53} With 2D MS, locating modifications can be accomplished visually (cf. Fig. 4a) and relative intensities can be plotted without distinguishing m/z ratios. In Fig. 4d, the relative intensities are given for both N-terminus and C-terminus fragments.

Fig. 4e compares the quantification obtained by using peaks on the autocorrelation line (6+), the one-electron capture line (5+), the two-electron capture line (4+), and the average quantification used with the fragment ion peaks (fragment average). The proportions of the four peptide isoforms are consistent between all four methods. The only discrepancy was in the

proportions of the two-electron capture line, where K7 2m was not detected, due to the low abundance of the $[M+6H]^{4+}$ charge state.

CONCLUSION

This study shows that decreasing the sampling frequency along the t_1 axis (i.e. precursor ion modulation) in 2D MS can substantially increase the resolving power in the vertical precursor ion m/z dimension for an unchanged acquisition time. Narrowband 2D MS can resolve isotopic distributions of the $[M+6H]^{6+}$ ions in a mixture of four histone peptide isoforms (26 residues). The precursor-fragment correlation was high enough to enable the detailed examination of fragment isotopic distributions in the 2D mass spectrum. Finally, the study shows the capability of 2D MS to identify and locate histone modifications, and to quantify their relative abundances.

As this study demonstrates, narrowband 2D MS has potential for top-down mass spectrometry of biomolecules. The high precursor-fragment correlation can be used to distinguish the fragmentation patterns of overlapping isotopic distributions (e.g. for proteins with and without disulphide bonds) and easily identify and locate modifications by extracting dissociation lines. The capacity of 2D MS to analyse all analytes in a sample with identical experimental parameters in a single experiment enable label-free relative quantification. These characteristics of 2D MS enable the analysis of heavily modified biomolecule, such as histones, for the elucidation of the histone code.⁵⁴

ASSOCIATED CONTENT

Supporting Information

The Supporting Information is available free of charge on the ACS Publications website.

Peak assignment tables (broadband and narrowband 2D mass spectra), sequence coverage schemes (broadband and narrowband 2D mass spectra), slope and offset of dissociation lines, isotopic distributions in the horizontal dimension (Word file, .docx)

AUTHOR INFORMATION

Corresponding Author

Maria A. van Agthoven, Institute for Organic Chemistry, University of Innsbruck, Innrain 80/82, 6020 Innsbruck, Austria

Email: maria.van-agthoven@uibk.ac.at

Author Contributions

The manuscript was written through contributions of all authors.

ACKNOWLEDGMENT

M.v.A. dedicates this article to the memory of Christina M. Olsthoorn-de Munck. The authors thank Ms. Sarah Heel, Mr. Giovanni Calderisi and Mr. Michael Palasser for helpful conversations. M.v.A. and K.B. thank Der Wissenschaftsfonds (Austrian Science Fund, FWF) for the Lise Meitner Fellowship project M 2757-B.

REFERENCES

- Peterson, C. L.; Laniel, M.-A., Histones and histone modifications. *Curr. Biol.* **2004**, *14* (14), R546-R551.
- Moradian, A.; Kalli, A.; Sweredoski, M. J.; Hess, S., The top-down, middle-down, and bottom-up mass spectrometry approaches for characterization of histone variants and their post-translational modifications. *Proteomics* **2014**, *14* (4-5), 489-497.

- Taucher, M.; Breuker, K., Characterization of Modified RNA by Top-Down Mass Spectrometry. *Angew. Chem., Int. Ed.* **2012**, *51* (45), 11289-11292.
- Tran, J. C.; Zamborg, L.; Ahlf, D. R.; Lee, J. E.; Catherman, A. D.; Durbin, K. R.; Tipton, J. D.; Vellaichamy, A.; Kellie, J. F.; Li, M.; Wu, C.; Sweet, S. M. M.; Early, B. P.; Siuti, N.; LeDuc, R. D.; Compton, P. D.; Thomas, P. M.; Kelleher, N. L., Mapping intact protein isoforms in discovery mode using top-down proteomics. *Nature* **2011**, *480* (7376), 254-258.
- Zheng, Y.; Fornelli, L.; Compton, P. D.; Sharma, S.; Canterbury, J.; Mullen, C.; Zabrouskov, V.; Fellers, R. T.; Thomas, P. M.; Licht, J. D.; Senko, M. W.; Kelleher, N. L., Unabridged Analysis of Human Histone H3 by Differential Top-Down Mass Spectrometry Reveals Hypermethylated Proteoforms from MMSET/NSD2 Overexpression. *Mol. Cell. Proteomics* **2016**, *15* (3), 776-790.
- Zheng, Y.; Huang, X.; Kelleher, N. L., Epiproteomics: quantitative analysis of histone marks and codes by mass spectrometry. *Curr. Opin. Chem. Biol.* **2016**, *33*, 142-150.
- Dang, X.; Singh, A.; Spetman, B. D.; Nolan, K. D.; Isaacs, J. S.; Dennis, J. H.; Dalton, S.; Marshall, A. G.; Young, N. L., Label-Free Relative Quantitation of Isobaric and Isomeric Human Histone H2A and H2B Variants by Fourier Transform Ion Cyclotron Resonance Top-Down MS/MS. *J. Proteome Res.* **2016**, *15* (9), 3196-3203.
- van Agthoven, M. A.; Delsuc, M.-A.; Bodenhausen, G.; Rolando, C., Towards analytically useful two-dimensional Fourier transform ion cyclotron resonance mass spectrometry. *Anal. Bioanal. Chem.* **2013**, *405* (1), 51-61.
- van Agthoven, M. A.; Lam, Y. P. Y.; O'Connor, P. B.; Rolando, C.; Delsuc, M.-A., Two-dimensional mass spectrometry: new perspectives for tandem mass spectrometry. *Eur. Biophys. J.* **2019**, *48* (3), 213-229.
- van Agthoven, M. A.; Delsuc, M.-A.; Rolando, C., Two-dimensional FT-ICR/MS with IRMPD as fragmentation mode. *Int. J. Mass Spectrom.* **2011**, *306* (2-3), 196-203.
- van Agthoven, M. A.; O'Connor, P. B., Two-dimensional mass spectrometry in a linear ion trap, an in silico model. *Rapid Comm. Mass Spectrom.* **2017**, *31* (8), 674-684.
- van Agthoven, M. A.; Barrow, M. P.; Chiron, L.; Coutouly, M.-A.; Kilgour, D.; Wootton, C. A.; Wei, J.; Soulby, A.; Delsuc, M.-A.; Rolando, C.; O'Connor, P. B., Differentiating Fragmentation Pathways of Cholesterol by Two-Dimensional Fourier Transform Ion Cyclotron Resonance Mass Spectrometry. *J. Am. Soc. Mass Spectrom.* **2015**, *26* (12), 2105-2114.
- Simon, H. J.; van Agthoven, M. A.; Lam, P. Y.; Floris, F.; Chiron, L.; Delsuc, M. A.; Rolando, C.; Barrow, M. P.; O'Connor, P. B., Uncoiling collagen: a multidimensional mass spectrometry study. *Analyst* **2016**, *141* (1), 157-165.
- Floris, F.; van Agthoven, M.; Chiron, L.; Soulby, A. J.; Wootton, C. A.; Lam, Y. P. Y.; Barrow, M. P.; Delsuc, M.-A.; O'Connor, P. B., 2D FT-ICR MS of Calmodulin: A Top-Down and Bottom-Up Approach. *J. Am. Soc. Mass Spectrom.* **2016**, *27* (9), 1531-1538.
- van Agthoven, M. A.; Wootton, C. A.; Chiron, L.; Coutouly, M.-A.; Soulby, A.; Wei, J.; Barrow, M. P.; Delsuc, M.-A.; Rolando, C.; O'Connor, P. B., Two-Dimensional Mass Spectrometry for Proteomics, a Comparative Study with Cytochrome c. *Anal. Chem.* **2016**, *88* (8), 4409-4417.
- Floris, F.; Vallotto, C.; Chiron, L.; Lynch, A. M.; Barrow, M. P.; Delsuc, M.-A.; O'Connor, P. B., Polymer Analysis in the Second Dimension: Preliminary Studies for the Characterization of Polymers with 2D MS. *Anal. Chem.* **2017**, *89* (18), 9892-9899.
- Floris, F.; van Agthoven, M. A.; Chiron, L.; Wootton, C. A.; Lam, P. Y. Y.; Barrow, M. P.; Delsuc, M.-A.; O'Connor, P. B., Bottom-Up Two-Dimensional Electron-Capture Dissociation Mass Spectrometry of Calmodulin. *J. Am. Soc. Mass Spectrom.* **2018**, 207-210.
- Snyder, D. T.; Szalwinski, L. J.; Wells, J. M.; Cooks, R. G., Logical MS/MS scans: a new set of operations for tandem mass spectrometry. *Analyst* **2018**, *143* (22), 5438-5452.

19. Paris, J.; Morgan, T. E.; Wootton, C. A.; Barrow, M. P.; O'Hara, J.; O'Connor, P. B., Facile Determination of Phosphorylation Sites in Peptides Using Two-Dimensional Mass Spectrometry. *Anal. Chem.* **2020**, *92* (10), 6817-6821.
20. Snyder, D. T.; Szalwinski, L. J.; St. John, Z.; Cooks, R. G., Two-Dimensional Tandem Mass Spectrometry in a Single Scan on a Linear Quadrupole Ion Trap. *Analytical Chemistry* **2019**, *91* (21), 13752-13762.
21. Szalwinski, L. J.; Holden, D. T.; Morato, N. M.; Cooks, R. G., 2D MS/MS Spectra Recorded in the Time Domain Using Repetitive Frequency Sweeps in Linear Quadrupole Ion Traps. *Analytical Chemistry* **2020**, *92* (14), 10016-10023.
22. Marshall, A. G.; Hendrickson, C. L.; Jackson, G. S., Fourier transform ion cyclotron resonance mass spectrometry: a primer. *Mass. Spectrom. Rev.* **1998**, *17* (1), 1-35.
23. Amster, I. J., Fourier transform mass spectrometry. *J. Mass Spectrom.* **1996**, *31* (12), 1325-1337.
24. Denisov, E.; Damoc, E.; Lange, O.; Makarov, A., Orbitrap mass spectrometry with resolving powers above 1,000,000. *Int. J. Mass Spectrom.* **2012**, *325-327*, 80-85.
25. Little, D. P.; Speir, J. P.; Senko, M. W.; O'Connor, P. B.; McLafferty, F. W., Infrared multiphoton dissociation of large multiply charged ions for biomolecule sequencing. *Anal. Chem.* **1994**, *66* (18), 2809-15.
26. Zubarev, R. A.; Kelleher, N. L.; McLafferty, F. W., Electron Capture Dissociation of Multiply Charged Protein Cations. A Nonergodic Process. *J. Am. Chem. Soc.* **1998**, *120* (13), 3265-3266.
27. Hakansson, K.; Cooper, H. J.; Emmett, M. R.; Costello, C. E.; Marshall, A. G.; Nilsson, C. L., Electron capture dissociation and infrared multiphoton dissociation MS/MS of an N-glycosylated tryptic peptide to yield complementary sequence information. *Anal. Chem.* **2001**, *73* (18), 4530-4536.
28. Domon, B.; Aebersold, R., Mass Spectrometry and Protein Analysis. *Science* **2006**, *312* (5771), 212-217.
29. Drader, J. J.; Shi, S. D. H.; Blakney, G. T.; Hendrickson, C. L.; Laude, D. A.; Marshall, A. G., Digital Quadrature Heterodyne Detection for High-Resolution Fourier Transform Ion Cyclotron Resonance Mass Spectrometry. *Anal. Chem.* **1999**, *71* (20), 4758-4763.
30. Fujiwara, M.; Katakura, H.; Inoue, M., Quadrature heterodyne method for fourier-transform ion cyclotron resonance mass spectrometry. *Rapid Comm. Mass Spectrom.* **1990**, *4* (7), 237-238.
31. Wei, J.; Antzutkin, O. N.; Filippov, A. V.; Iuga, D.; Lam, P. Y.; Barrow, M. P.; Dupree, R.; Brown, S. P.; O'Connor, P. B., Amyloid Hydrogen Bonding Polymorphism Evaluated by ¹⁵N{(17)O}REAPDOR Solid-State NMR and Ultra-High Resolution Fourier Transform Ion Cyclotron Resonance Mass Spectrometry. *Biochemistry* **2016**, *55* (14), 2065-2068.
32. Taucher, M.; Rieder, U.; Breuker, K., Minimizing base loss and internal fragmentation in collisionally activated dissociation of multiply deprotonated RNA. *J. Am. Soc. Mass Spectrom.* **2010**, *21* (2), 278-285.
33. Caravatti, P.; Allemann, M., The infinity cell: a new trapped-ion cell with radiofrequency covered trapping electrodes for Fourier transform ion cyclotron resonance mass spectrometry. *Org. Mass Spectrom.* **1991**, *26* (5), 514-18.
34. van Agthoven, M. A.; Chiron, L.; Coutouly, M.-A.; Sehgal, A. A.; Pelulessy, P.; Delsuc, M.-A.; Rolando, C., Optimization of the discrete pulse sequence for two-dimensional FT-ICR mass spectrometry using infrared multiphoton dissociation. *Int. J. Mass Spectrom.* **2014**, *370* (0), 114-124.
35. Tsybin, Y. O.; Witt, M.; Baykut, G.; Kjeldsen, F.; Hakansson, P., Combined infrared multiphoton dissociation and electron capture dissociation with a hollow electron beam in Fourier transform ion cyclotron resonance mass spectrometry. *Rapid Commun. Mass Spectrom.* **2003**, *17* (15), 1759-1768.
36. Chiron, L.; Coutouly, M.-A.; Starck, J.-P.; Rolando, C.; Delsuc, M.-A., SPIKE a processing software dedicated to Fourier spectroscopies. *arXiv.org, e-Print Arch., Phys.* **2016**, 1-13.
37. Chiron, L.; van Agthoven, M. A.; Kieffer, B.; Rolando, C.; Delsuc, M.-A., Efficient denoising algorithms for large experimental datasets and their applications in Fourier transform ion cyclotron resonance mass spectrometry. *Proc. Natl. Acad. Sci. U. S. A.* **2014**, *111* (4), 1385-1390.
38. van Agthoven, M. A.; Coutouly, M.-A.; Rolando, C.; Delsuc, M.-A., Two-dimensional Fourier transform ion cyclotron resonance mass spectrometry: reduction of scintillation noise using Cadzow data processing. *Rapid. Commun. Mass Spectrom.* **2011**, *25* (11), 1609-1616.
39. Bray, F.; Bouclon, J.; Chiron, L.; Witt, M.; Delsuc, M.-A.; Rolando, C., Nonuniform Sampling Acquisition of Two-Dimensional Fourier Transform Ion Cyclotron Resonance Mass Spectrometry for Increased Mass Resolution of Tandem Mass Spectrometry Precursor Ions. *Anal. Chem.* **2017**, *89* (17), 8589-8593.
40. Ledford, E. B., Jr.; Rempel, D. L.; Gross, M. L., Space charge effects in Fourier transform mass spectrometry. II. Mass calibration. *Anal. Chem.* **1984**, *56* (14), 2744-8.
41. Guan, S.; Jones, P. R., A theory for two-dimensional Fourier-transform ion cyclotron resonance mass spectrometry. *J. Chem. Phys.* **1989**, *91* (9), 5291-5.
42. Pfändler, P.; Bodenhausen, G.; Rapin, J.; Walser, M. E.; Gäumann, T., Broad-band two-dimensional Fourier transform ion cyclotron resonance. *J. Am. Chem. Soc.* **1988**, *110* (17), 5625-5628.
43. Jerri, A. J., The Shannon sampling theorem; Its various extensions and applications: A tutorial review. *Proceedings of the IEEE* **1977**, *65* (11), 1565-1596.
44. Zubarev, R. A.; Horn, D. M.; Fridriksson, E. K.; Kelleher, N. L.; Kruger, N. A.; Lewis, M. A.; Carpenter, B. K.; McLafferty, F. W., Electron Capture Dissociation for Structural Characterization of Multiply Charged Protein Cations. *Anal. Chem.* **2000**, *72* (3), 563-573.
45. van Agthoven, M. A.; Kilgour, D. P. A.; Lynch, A. M.; Barrow, M. P.; Morgan, T. E.; Wootton, C. A.; Chiron, L.; Delsuc, M.-A.; O'Connor, P. B., Phase relationships in two-dimensional mass spectrometry. *J. Am. Soc. Mass Spectrom.* **2019**, *30* (12), 2594-2607.
46. van Agthoven, M. A.; Lynch, A. M.; Morgan, T. E.; Wootton, C. A.; Lam, Y. P. Y.; Chiron, L.; Barrow, M. P.; Delsuc, M.-A.; O'Connor, P. B., Can Two-Dimensional IR-ECD Mass Spectrometry Improve Peptide de Novo Sequencing? *Anal. Chem.* **2018**, *90* (5), 3496-3504.
47. Grosshans, P. B.; Marshall, A. G., Can Fourier transform mass spectral resolution be improved by detection at harmonic multiples of the fundamental ion cyclotron orbital frequency? *Int. J. Mass Spectrom. Ion Process.* **1991**, *107* (1), 49-81.
48. Vorobyev, A.; Gorshkov, M. V.; Tsybin, Y. O., Towards data acquisition throughput increase in Fourier transform mass spectrometry of proteins using double frequency measurements. *Int. J. Mass Spectrom.* **2011**, *306* (2), 227-231.
49. Park, S.-G.; Anderson, G. A.; Bruce, J. E., Parallel Detection of Fundamental and Sixth Harmonic Signals Using an ICR Cell with Dipole and Sixth Harmonic Detectors. *J. Am. Soc. Mass Spectrom.* **2020**, *31* (3), 719-726.
50. O'Connor, P. B.; Little, D. P.; McLafferty, F. W., Isotopic Assignment in Large-Molecule Mass Spectra by Fragmentation of a Selected Isotopic Peak. *Anal. Chem.* **1996**, *68* (3), 542-5.
51. van Agthoven, M. A.; Chiron, L.; Coutouly, M.-A.; Delsuc, M.-A.; Rolando, C., Two-Dimensional ECD FT-ICR Mass Spectrometry of Peptides and Glycopeptides. *Anal. Chem.* **2012**, *84* (13), 5589-5595.
52. Glasner, H.; Riml, C.; Micura, R.; Breuker, K., Label-free, direct localization and relative quantitation of the RNA nucleobase methylations m6A, m5C, m3U, and m5U by top-down mass spectrometry. *Nucleic Acids Res.* **2017**, *45* (13), 8014-8025.
53. Yuan, Z.-F.; Lin, S.; Molden, R. C.; Cao, X.-J.; Bhanu, N. V.; Wang, X.; Sidoli, S.; Liu, S.; Garcia, B. A., EpiProfile Quantifies Histone Peptides With Modifications by Extracting Retention Time and Intensity in High-resolution Mass Spectra. *Mol. Cell. Proteomics* **2015**, *14* (6), 1696-1707.
54. Jenuwein, T.; Allis, C. D., Translating the histone code. *Science* **2001**, *293* (5532), 1074-1080.

Table of Content:

