



**HAL**  
open science

## **Petroleomics by Direct Analysis in Real Time-Mass Spectrometry**

Wanderson Romão, Lilian V Tose, Boniek G Vaz, Sara G Sama, Ryszard Lobinski, Pierre Giusti, Hervé Carrier, Brice Bouyssière

► **To cite this version:**

Wanderson Romão, Lilian V Tose, Boniek G Vaz, Sara G Sama, Ryszard Lobinski, et al.. Petroleomics by Direct Analysis in Real Time-Mass Spectrometry. *Journal of The American Society for Mass Spectrometry*, 2016, 27 (1), pp.182-185. 10.1007/s13361-015-1266-z . hal-03133641

**HAL Id: hal-03133641**

**<https://hal.science/hal-03133641>**

Submitted on 10 Feb 2021

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

# Petroleomics by Direct Analysis in Real Time-Mass Spectrometry

Wanderson Romão,<sup>1,2 †</sup> Lílian V. Tose,<sup>1</sup> Boniek G. Vaz,<sup>3</sup> Sara G. Sama,<sup>4,5</sup> R. Lobinski,<sup>4</sup> P. Giusti,<sup>5</sup> Hervé Carrier,<sup>6</sup> and Brice Bouyssiere<sup>4 ‡</sup>

<sup>1</sup> Laboratório de Petroleômica e Forense, Departamento de Química, Universidade Federal do Espírito Santo, 29075-910, Vitória, ES, Brazil.

<sup>2</sup> Instituto Federal de Educação, Ciência e Tecnologia do Espírito Santo, 29106-010, Vila Velha, ES, Brazil.

<sup>3</sup> Instituto de Química, Universidade Federal de Goiás, 74001-970, Goiânia, GO, Brazil.

<sup>4</sup> LCABIE-IPREM, Université de Pau et des Pays de l'Adour, Hélioparc, 2 Av. Pr. Angot, 64053 Pau CEDEX, France.

<sup>5</sup> TOTAL Raffinage Chimie, TRTG, BP 27, 76700 Harfleur, France .

<sup>6</sup> LFR-R, Université de Pau et des Pays de l'Adour, Av. de l'Université, BP 576, 64012 Pau CEDEX, France.

Corresponding author:

<sup>†</sup>wandersonromao@gmail.com / Phone: + 55-27-3149-0833

<sup>‡</sup>[brice.bouyssiere@univ-pau.fr](mailto:brice.bouyssiere@univ-pau.fr) / Phone: + 33 559 407 752

23

24 **Abstract**

25         The analysis of crude oil and its fractions applying ambient ionization techniques is  
26 yet under explored in mass spectrometry (MS). Direct Analysis in Real Time (DART) in  
27 positive ion mode detection was coupled to linear ion trap quadrupole (LTQ) and Orbitrap  
28 mass spectrometer and optimized to analyze crude oil and paraffin samples. The ionization  
29 and acquisition parameters of the DART-MS such as the template substrates (paper, TLC  
30 plate and glass), temperature (from 100 up to 400 °C), carrier gas (helium and nitrogen),  
31 concentration of analyte (from 0.33 to 6 mg mL<sup>-1</sup>) and acquisition time (from 1 to 10 μscans)  
32 were optimized for crude oil analyzes. DART-MS rendered the optimum conditions of  
33 operation using paper as a substrate, T = 400 °C, helium as a carrier gas, sample concentration  
34 ≥ 6 mg mL<sup>-1</sup>, and acquisition time < 2 μscans. For crude oils analyzes, DART(+)-Orbitrap  
35 mass spectra detected nitrogen-containing protonated species, whereas for paraffin samples,  
36 hydroxylated HCs species (O<sub>x</sub> classes, where x = 1-4) with DBEs of 1-4 were detected, being  
37 their structures and connectivity confirmed by CID experiments (MS<sup>2</sup>). The DART(+)-MS  
38 and CID experiments (MS<sup>2</sup> and MS<sup>3</sup>) were also able to identify porphyrin standard  
39 compounds as [M + H]<sup>+</sup> ions of *m/z* 615.2502 and 680.1763, where M = C<sub>44</sub>H<sub>30</sub>N<sub>4</sub> and  
40 C<sub>44</sub>H<sub>28</sub>N<sub>4</sub>OV, respectively.

41 **Key-words:** ambient mass spectrometry, DART-MS, crude oil, paraffin, porphyrin.

42

43 **1. Introduction**

44

45 A paragraph on interest of what you have done

46 A paragraph on other analytical techniques that were applied to what you have done  
47 and their limitations.

48 The objective of this work was to investigate the application of a direct analysis in  
49 real time (DART) for this purpose. DART developed by Cody and Laramée,<sup>1</sup> is a type of  
50 ambient pressure ionization technique. Its coupling with mass spectrometry has allowed  
51 many applications in different fields such as forensic,<sup>2-6</sup> pharmaceutical,<sup>7-9</sup> food,<sup>10-14</sup> analysis,  
52 in biology,<sup>15-17</sup> chemistry,<sup>18-20</sup> In crude oil analysis <sup>21</sup> there are no application of DART  
53 except the work of Rummel et al.,<sup>22</sup> where it was coupled to Fourier Transform Ion Cyclotron  
54 Resonance (FT-ICR) mass spectrometer . Here, the DART source was coupled to a hybrid  
55 mass analyzer (LTQ Orbitrap Velos Pro™) and the ionization and acquisition parameters  
56 (substrates, temperature and type of gas heater, concentration of analyte and acquisition time)  
57 were optimized for crude oil analyzes. The ability of DART source in paraffin and porphyrin  
58 compounds ionization was also evaluated.

59

## 60 ***2. Procedure***

### 61 ***2.1 Reagents and samples***

62 Dichloromethane, heptane, and tetrahydrofuran, THF, (analytical grades with purity  
63 higher than 99.5%) were supplied by Sigma–Aldrich Chemicals USA and were used to  
64 prepare solutions of crude oil, paraffin and porphyrin standard compounds.

65 Seven crude oil samples, named samples A to G, were supplied by PETROBRAS and  
66 characterized to determine the API degree (ASTM D1298-99) and saturates, aromatics,  
67 resins, and asphaltenes (SARA) content. For all the crude oils evaluated, API degree values

68 ranged from 26 to 30, classifying them as medium crude oils (API degree = 22-30). Saturates  
69 content ranged from 51 to 65 wt %. To evaluate the detection sensibility of DART technique,  
70 samples were diluted in dichloromethane in five different concentrations: 0.33; 1.0; 1.6; 3.3  
71 and 6 mg mL<sup>-1</sup>. A volume of 10µL was spotted on the paper surface and analysed by DART-  
72 MS. Two other surfaces were tested as substrates: glass and a thin layer chromatography  
73 (TLC) plate. The stationary phase of TLC is composed of silica gel.

74 Three saturated hydrocarbon samples were also studied in this work, named paraffin  
75 A, B and C. The first two samples, paraffins A and B, were purchased from Sigma–Aldrich  
76 Chemicals USA and Vetec Química Fina Ltda, respectively, whereas the paraffin C was  
77 obtained from a food grade process. The organic solutions were prepared in heptane at ≈ 10  
78 mg mL<sup>-1</sup> and a volume of 10µL was spotted on the TLC plate.

79 The porphyrin standard compounds such as 5,10,15,20-tetraphenyl-21H,23H-  
80 porphine and 5,10,15,20-tetraphenyl-21H,23H-porphine vanadium(IV) oxide, with  
81 molecular formula (M) equal to C<sub>44</sub>H<sub>30</sub>N<sub>4</sub> and C<sub>44</sub>H<sub>28</sub>N<sub>4</sub>OV and molecular weight of  
82 614.2471 and 679.1703, respectively, were also studied. Both porphyrin standard compounds  
83 were purchased from Sigma–Aldrich Chemicals USA. Solutions were prepared at 100 mg L<sup>-1</sup>  
84 in THF and after, a volume of 10µL was spotted on paper substrate and analysed by  
85 DART(+)-MS.

86

## 87 **2.2 DART(+)-MS**

88 For the experiments, DART–MS system consisting of a DART ion source (IonSense,  
89 Saugus, MA, USA) was coupled with a hybrid mass spectrometer: LTQ Orbitrap Velos Pro™  
90 (Thermo Fisher Scientific, Bremen, Germany). The operating conditions of the DART ion

91 source were as follows: positive ion mode; helium flow: 4.0 L min<sup>-1</sup>; discharge needle  
92 voltage: 3.0 kV; perforated and grid electrode potentials: +150 and +350 V, respectively. The  
93 distance between the DART gun exit and mass spectrometer inlet was about 5-10 mm. For  
94 glasses and paper surfaces, the sample introductions were carried out manually, whereas for  
95 TLC plate substrate it was automatically performed using Dip-it holder samplers. To assess  
96 the influence of the gas beam temperature on the signal intensity, crude oil spots were  
97 analyzed at different temperatures ranging from 100 to 400 °C using helium and nitrogen as  
98 the gas beam.

99 DART(+) mass spectra were acquired using both a LTQ and an Orbitrap mass  
100 analyzers. For high resolution experiments using Orbitrap mass analyzer, a mass resolving  
101 power of 100,000 (FWHM, at  $m/z$  400; 0.5 s scan cycle time) was reached. As consequence,  
102 a mass error about 3 to 5 ppm was measured, where the mass accuracy is determined from  
103 
$$\text{Error (ppm)} = ((m/z_{\text{measured}} - m/z_{\text{theoretical}})/m/z_{\text{theoretical}}) \times 10^6.$$

104 The maximum ion injection time was about 10 and 300 ms for LTQ and Orbitrap,  
105 respectively, with the automatic gain control (AGC, corresponding to the number of changes  
106 transferred from the front-stage ion trap to the orbitrap analyzer) target set at  $1 \times 10^5$ . The  
107 full scan mass spectra were acquired over the range of  $m/z$  150–1000. Tandem mass  
108 spectrometry (DART-MS/MS) was also performed by collision-induced  
109 dissociation with a collision energy of 15–30% (manufacturer's unit) using LTQ as  
110 mass analyzer.

111 For crude oil samples, the DART(+)-Orbitrap mass spectra were acquired and processed  
112 using Composer software (Sierra Analytics, Pasadena, CA, USA). The MS data were  
113 processed and the elemental composition of the sample was determined by measuring the

114  $m/z$  values. Class and DBE distributions, and carbon number (CN) *versus* DBE graphs were  
115 plotted to better analyze the results. DBE is defined as the number of rings added and the  
116 number of double bonds in each molecular structure. The unsaturation level of each  
117 compound can be deduced directly from its DBE value according to equation 1:<sup>23,24,25</sup>

$$118 \quad \text{DBE} = c - h/2 + n/2 + 1 \text{ (equation 1)}$$

119 Where  $c$ ,  $h$ , and  $n$  are the numbers of carbon, hydrogen, and nitrogen atoms,  
120 respectively, in the molecular formula.

121

122

123


### 124 **3. Results**

#### 125 **3.1 DART(+)-MS Optimization**

126 The sensibility of DART(+)-MS technique for crude oil analyzes as function of some  
127 parameters were evaluated such as: substrate (TLC plate, glass and paper, **Figure 2**);  
128 temperature (from 100 to 400 °C, **Figure 3**); type of gas heater (He and N<sub>2</sub>, **Figure 4**);  
129 concentration of analyte (from 0.33 to 6 mg mL<sup>-1</sup>, **Figure 5**); and acquisition time for  
130 Orbitrap experiments (number of microscans, **Figure 8**).

131 **Figures 2a-c** show DART(+)-LTQ mass spectra (using Helium as heater at 400 °C) of  
132 a typical Brazilian crude oil (sample A at 10 mg mL<sup>-1</sup>) as a function of different substrates:  
133 **(a)** TLC; **(b)** glass; and **(c)** paper. A higher amplitude and distribution of signals with profiles  
134 ranging from  $m/z$  200-800 was observed when paper substrate was applied. The sensibility  
135 of DART(+) source increased in the following order: TLC < glass < paper. As consequence,

136 the average molar mass distribution ( $M_w$ ) was shifted for higher values of  $m/z$  (from  $m/z$  280  
137 to 360). A lower chemical interaction between polar organic compounds and the cellulose  
138  $((C_6H_{10}O_5)_n)$  was the key contributor to have a better ionization efficiency using paper as  
139 substrate.

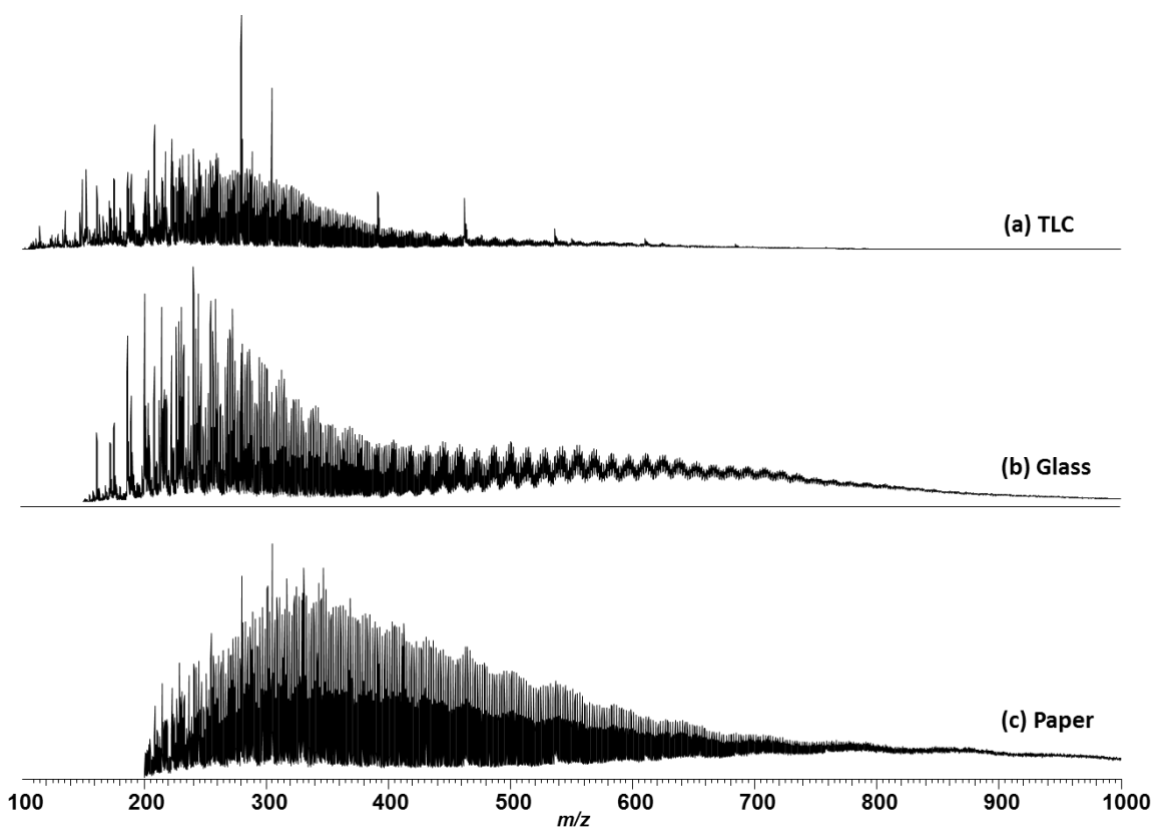


140  
141 **Figure 2.** DART(+)-LTQ MS for a typical Brazilian crude oil sample on the (a) TLC plate,  
142 (b) glass and (c) paper substrates.

143 The ionization efficiency of DART source was also tested as a function of  
144 temperature, **Figure 3**, and the type of the gas heater (Helium or  $N_2$ ), **Figure 4**. **Figure 3**  
145 shows the DART(+)-LTQ mass spectra as a function of temperature using Helium at (3a)  
146 100°C, (3b) 200°C, (3c) 300 °C, and (3d) 400 °C. It is possible to note that a higher sensibility  
147 and a higher amplitude of signals with  $m/z$  from 200-1000 and  $M_w = 400$  Da was observed




148 with 400 °C as the optimum temperature, **Figure 3d**. Changing the He gas heater to N<sub>2</sub>  
149 (**Figures 4a-b**), the DART(+)-LTQ mass spectrum at 400 °C, **Figure 4b**, showed a similar  
150 performance to that one with He at 300 °C, **Figure 4c**, thus proving the better efficiency of  
151 molecules ionization in presence of He due to its higher internal energy of ionization (He =  
152 19.8 eV *versus* N<sub>2</sub> = 15.6 eV).

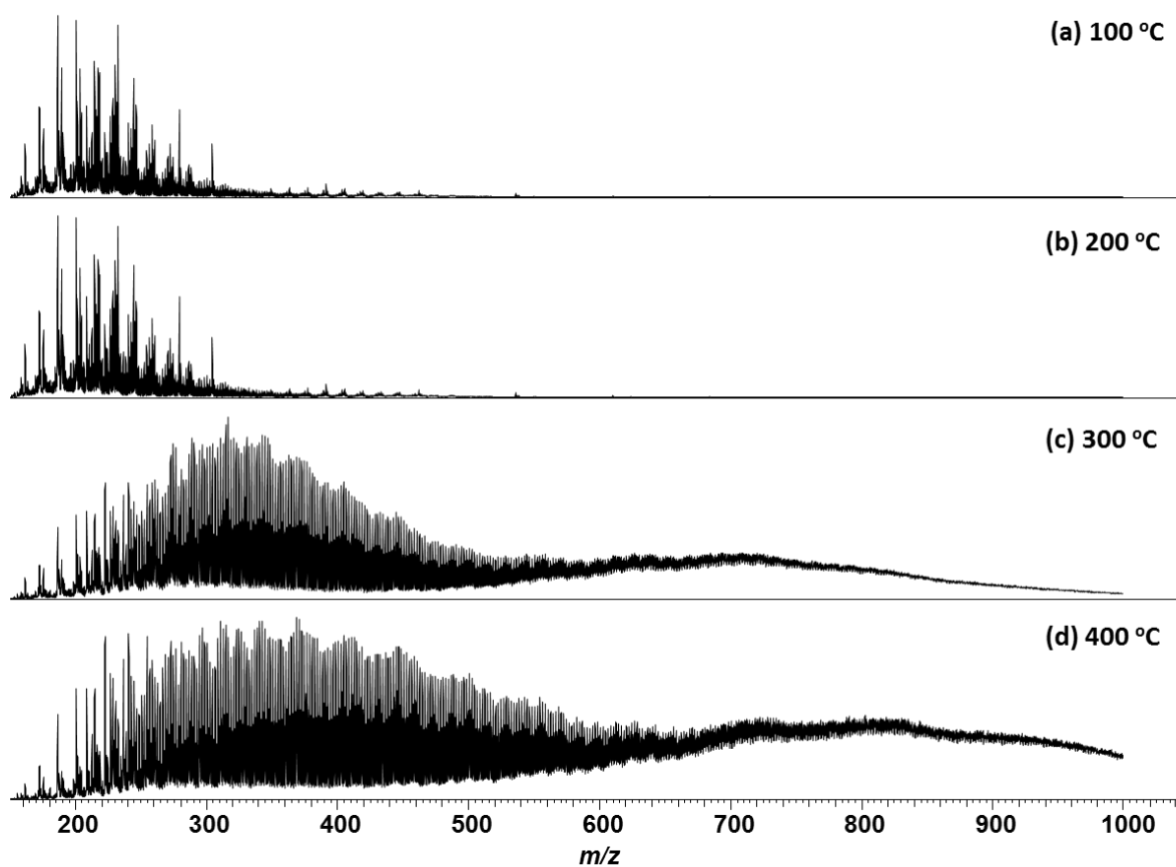

153  
154 **Figure 3.** DART(+)-LTQ MS on a Brazilian crude oil solution spot (concentration of 10 mg  
155 mL<sup>-1</sup>) as a function of temperature using He as gas heater at (a) 100, (b) 200, (c) 300 and (d)  
156 400 °C.

157

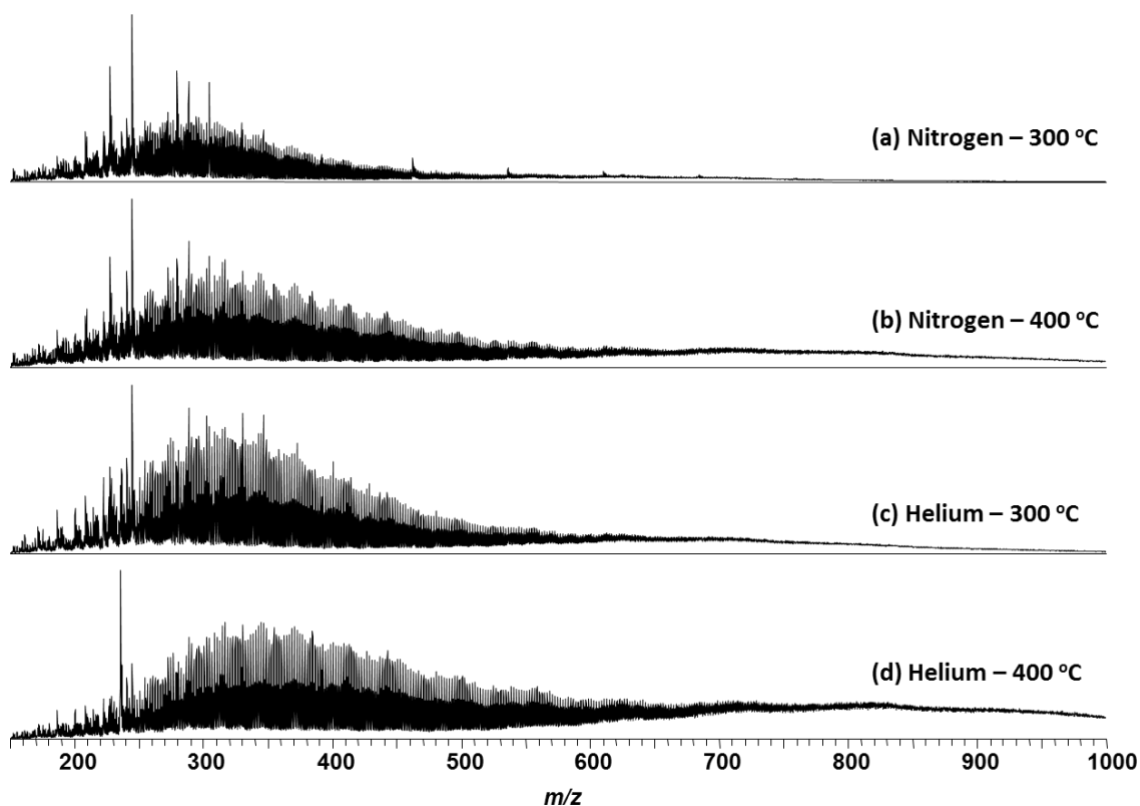


158

159 **Figure 4.** DART(+)-LTQ MS for a crude oil solution spot as a function of gas heater:  
 160 nitrogen and helium at 300 (a and c) and 400 °C (b and d), respectively.

161

162 The sensibility of DART(+)-MS technique was also evaluated as a function of  
 163 concentration of crude oil (from 0.33 to 6 mg mL<sup>-1</sup>). As consequence of increasing crude oil  
 164 concentration, a higher signal-to-noise ratio and amplitude of signals were easily observed  
 165 (see the decreasing of relative intensity of ions of  $m/z$  279, 329, 346 and 411 formed from  
 166 the paper substrate). Therefore, it is suggested that concentrations higher than 6 mg mL<sup>-1</sup>  
 167 must be used for crude oil analysis.

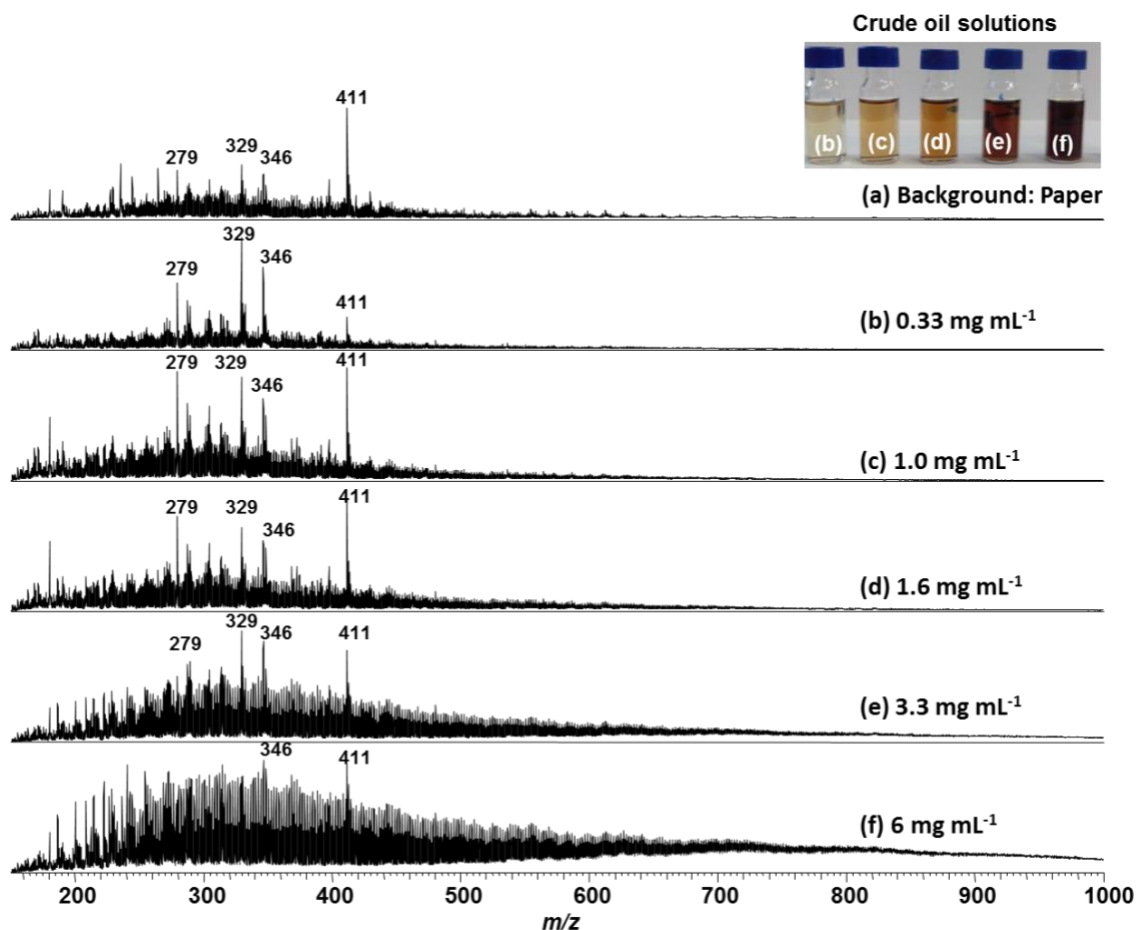

168

169 **Figure 5.** DART(+)-LTQ MS as a function of crude oil concentration: (a) 0, (b) 0.33 (c) 1.0,  
 170 (d) 1.6, (e) 3.3 and (f) 6 mg mL<sup>-1</sup>. He at 400 °C was used as gas heater.

171


### 172 **3.2 Crude oil analyzes**

173 After optimization of DART source ionization conditions as following: substrate  
 174 (paper), temperature (T = 400 °C), gas heater (Helium) and concentration ( $\geq 6$  mg mL<sup>-1</sup>),  
 175 DART(+) mass spectra at high resolution (FWHM  $\approx 100,000$  at  $m/z$  400) were acquired for  
 176 seven crude oil samples (samples A-G), using Orbitrap mass analyzer, **Figure 6.**

177           **Figures 6a-g** show the DART(+)-Orbitrap mass spectra of seven typical Brazilian  
178 crude oil samples, showing peaks profile from  $m/z$  200-700 with  $M_w$  centered from  $m/z$  374  
179 to 453. Nitrogen-containing species were detected as protonated molecules, that is,  $[M-H]^+$   
180 cations, according to the proposed mechanism of **equation 3**. The magnified region near  $m/z$   
181 310-312 indicates that DART(+) detected pyridine analogous compounds (N class):  
182  $[C_{23}H_{21}N + H]^+$ ,  $[C_{22}H_{31}N + H]^+$ ,  $[C_{23}H_{21}N + H]^+$  and  $[C_{22}H_{33}N + H]^+$  ions with  $m/z$   
183 310.1577, 310.2517, 312.1734 and 312.2674, and DBEs of 15, 8, 13 and 7, respectively. The  
184 theoretical  $m/z$  values for these ions are 310.1590, 310.2529, 312.1734 and 312.2674, thus  
185 providing a medium mass error of about -4.02 ppm. In Petroleomic, ultra-high resolution and  
186 accuracy mass spectrometry with FWHM > 400,000 and exact mass < 1 ppm is required for  
187 the identification of complex organic mixtures, thus ensuring an excellent recalibration data  
188 from Composer *software*. Accurate mass measurements define the unique elemental  
189 composition ( $C_cH_hN_nO_oS_s$ ) and DBE from singly charged ions.<sup>26,27</sup> To correct the mass  
190 deviation (Error > 1 ppm), DART(+)-Orbitrap mass spectra were further processed with the  
191 Composer *software*, especially designed for formula attribution via automatic recalibration  
192 for known homologous series from the measured  $m/z$  values of polar crude oil markers.<sup>28,29</sup>  
193 A medium of approximately 350 molecular formula (where  $n = 7$ ) were assigned from  
194 monoisotopic components for the Orbitrap mass spectra, corresponding to a medium  
195 percentage of 70 % of all assignments.

196

197

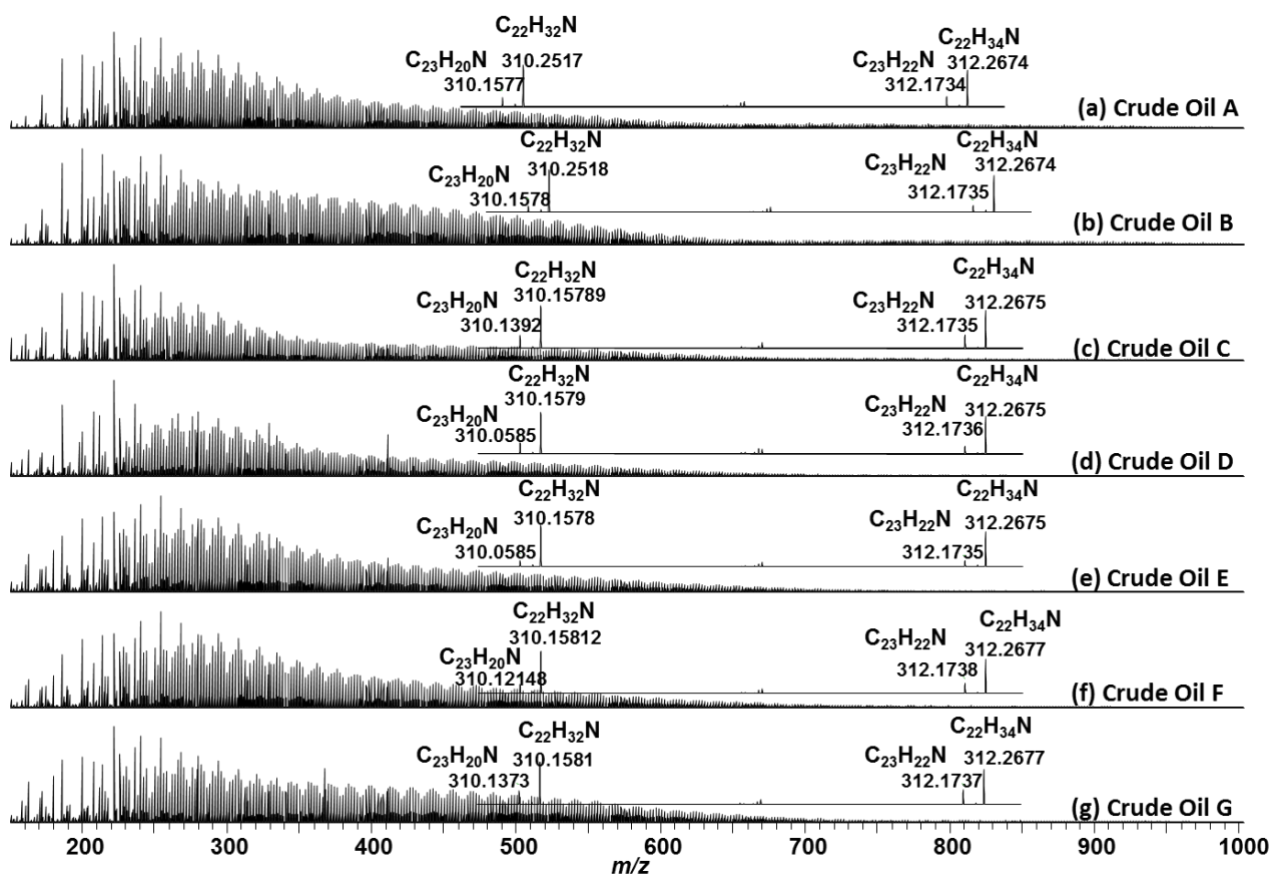

198


199 **Figure 6.** DART(+)-Orbitrap MS for seven crude oil samples (A-G). He was used as gas  
 200 heater at 400 °C. Acquisition time of 1 microscan.

201

202 One way to display the similarities or differences between the signal patterns of crude  
 203 oil samples is the construction of certain types of plots, such as the plots of the relative  
 204 abundances of different classes of compounds, DBE vs intensity and DBE vs CN,<sup>30</sup> **Figure**  
 205 **7a-c.** **Figure 7a** displays the distribution of polar compound classes (NH), NO[H], and O[H])  
 206 obtained from DART(+)-Orbitrap MS. In all cases, DART(+) seems to efficiently promoted  
 207 the ionization of polar compounds, as protonated cations ( $[M + H]^+$ ), with their magnitude  
 208 following the order: N[H] > O[H] > NO[H]. The DBE relative abundance distributions of

209 N[H] class for samples A-G were also evaluated, **Figure 7b**, in which a distribution ranging  
210 from 4 to 20 was observed. **Figure 7c** presents the DBE *versus* CN for the majority class,  
211 N[H] class, for the sample F. Carbon numbers ranged from C<sub>12</sub> to C<sub>45</sub> for pyridine compound  
212 species (DBEs = 4-20), with maximum abundance around C<sub>16</sub> and DBE = 8 were observed.


213 An attempt to improve the exact mass and consequently the signal-to-noise ratio in  
214 an Orbitrap analyzer was made by increasing its acquisition time. **Figure 8** shows the  
215 DART(+)-Orbitrap mass spectra as a function of the acquisition time (number of microscans)  
216 for sample F. The M<sub>w</sub> decreased as the number of microscan increased (1 → 10) as well as  
217 the population of nitrogen compounds assignments, depicted by the DBE vs CN plot (see the  
218 insert of **Figure 8**). The number of assigned molecular formulas decreases from 288 (for 1  
219 μscan) to 199 (for 10 μscans). Probably, the ions transmission from the LTQ to Orbitrap was  
220 affected by fact that desorption and ionization mechanism of DART acts only in a specific  
221 point. Hence, the ionic population was reduced as function of time.

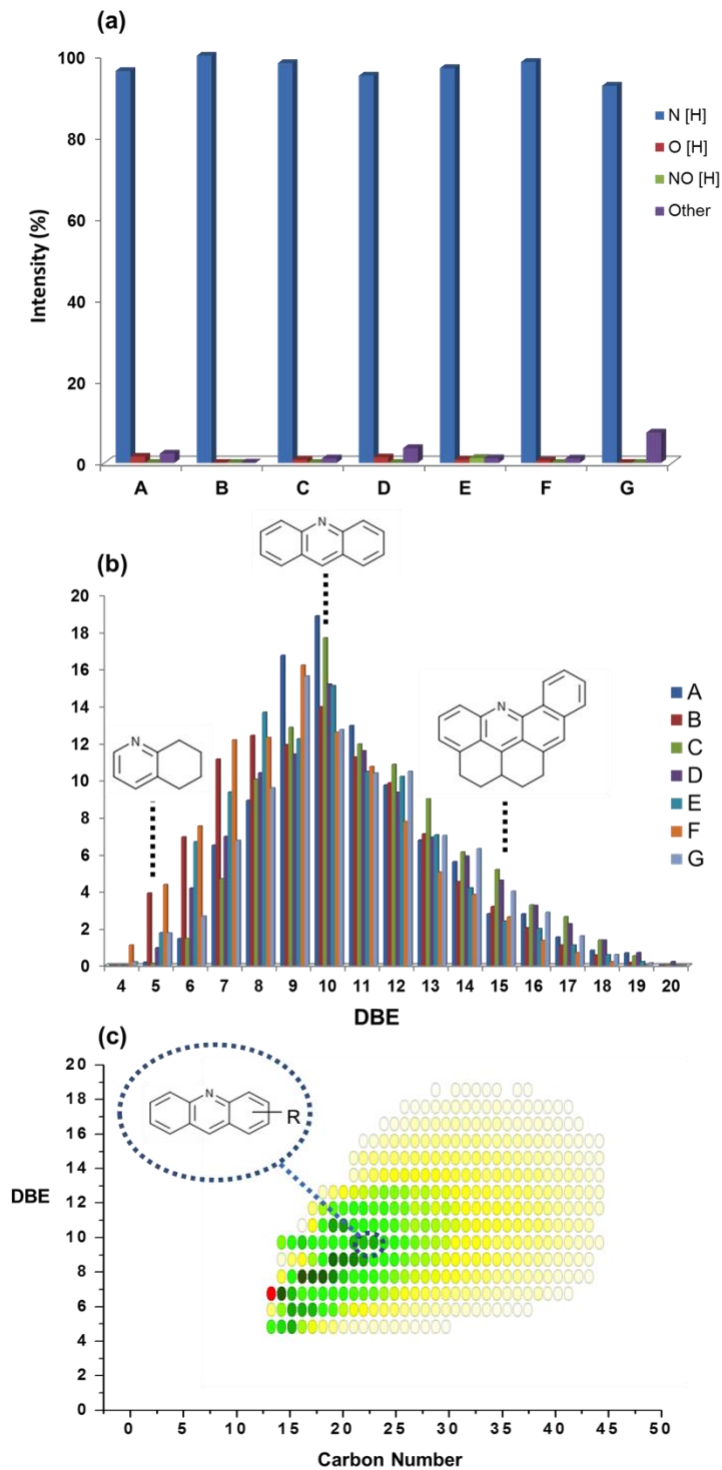

222

223 **Figure 7.** Class distribution, DBE *versus* intensity of seven crude oil samples (A-G) and

224 DBE *versus* CN plots of sample F generated from DART(+)Orbitrap data.

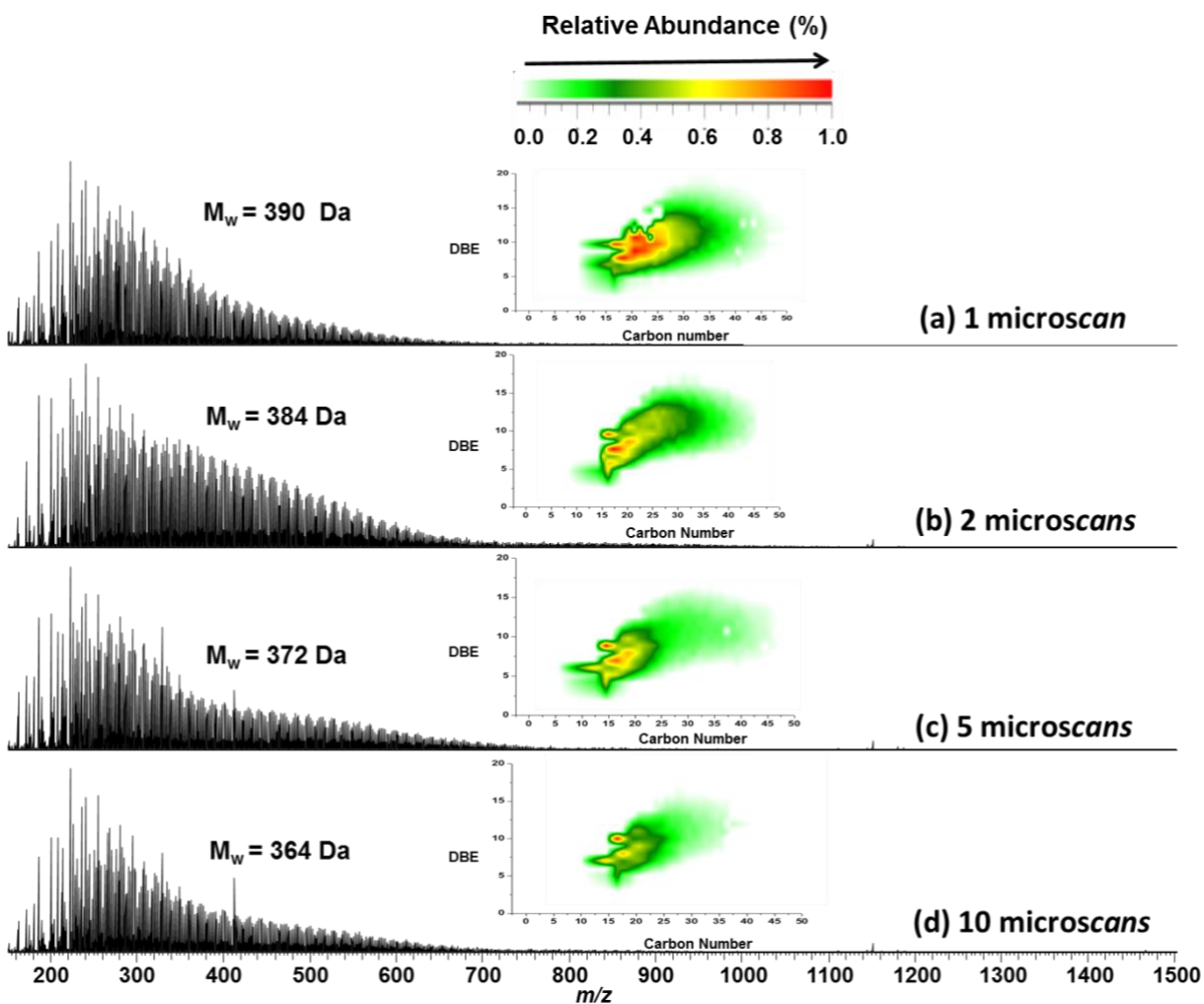

225

226 **Figure 8.** DART(+)-Orbitrap mass spectra as a function of acquisition time (number of  
 227 microscans) for crude oil sample A. Note that  $M_w$  decreases as the number of microscan  
 228 increases as well as the nitrogen compounds assignments decreases as showed by DBE vs  
 229 CN plot (see the insert of **Figure 8**).

230

### 231 3.3 Paraffin Analyzes

232 The analysis of hydrocarbons (HCs) using atmospheric or ambient ionization  
 233 techniques still remain a challenge in mass spectrometry.<sup>31</sup> **Figure 9** shows DART(+)-  
 234 Orbitrap mass spectra for the three paraffin samples evaluated. The right side inserts show




235 paraffin detection as oxygenated HCs species ( $O_x$  classes, where  $x = 1-4$ ) with DBEs of 1-4  
236 and a mass error of - 3-4 ppm. In all cases, a similar Gaussian profile from  $m/z$  250-600 was  
237 observed. The oxidized HCs species were generated at high temperature and atmospheric  
238 pressure from short-life time oxygen-based species such as OH and OOH radicals and also  
239  $H_3O^+$  ion in contact with HCs compound species, producing oxidized HCs ( $O_x$  classes).<sup>31,32</sup>

240 In 2009, Cooks et al.<sup>33</sup> has also analyzed saturated hydrocarbons ( $C_{15}H_{32}$  to  $C_{30}H_{62}$ )  
241 using discharge-induced oxidation in desorption electrospray ionization. Multiple oxidations  
242 and dehydrogenations occurred during the DESI discharge, but no hydrocarbon  
243 fragmentation was observed. DESI-Orbitrap mass spectrum of a petroleum distillate  
244 containing vacuum gas oil saturates (boiling point > 316 °C) showed HCs species containing  
245 two oxygen additions in alkanes structures from  $C_{21}H_{44}$  to  $C_{36}H_{74}$ , similar to observed for  
246 DART(+)-Orbitrap data.

247 To confirm the structures and the connectivity of some oxygenated HCs compound  
248 classes ( $O_x$  classes), which were identified using DART(+)-Orbitrap MS, the DART(+)-  
249 MS/MS spectra were acquired for the ions with  $m/z$  337, 351, 365, 379 and 393. This  
250 approach identifies the characteristic neutral loss and confirms the existence of hydroxylated  
251 HCs compounds such as alcohols from successive eliminations of 18 Da ( $H_2O$ ), and 28 Da  
252 ( $CH_2=CH_2$ ) along the molecule, **Figure 10**.


253

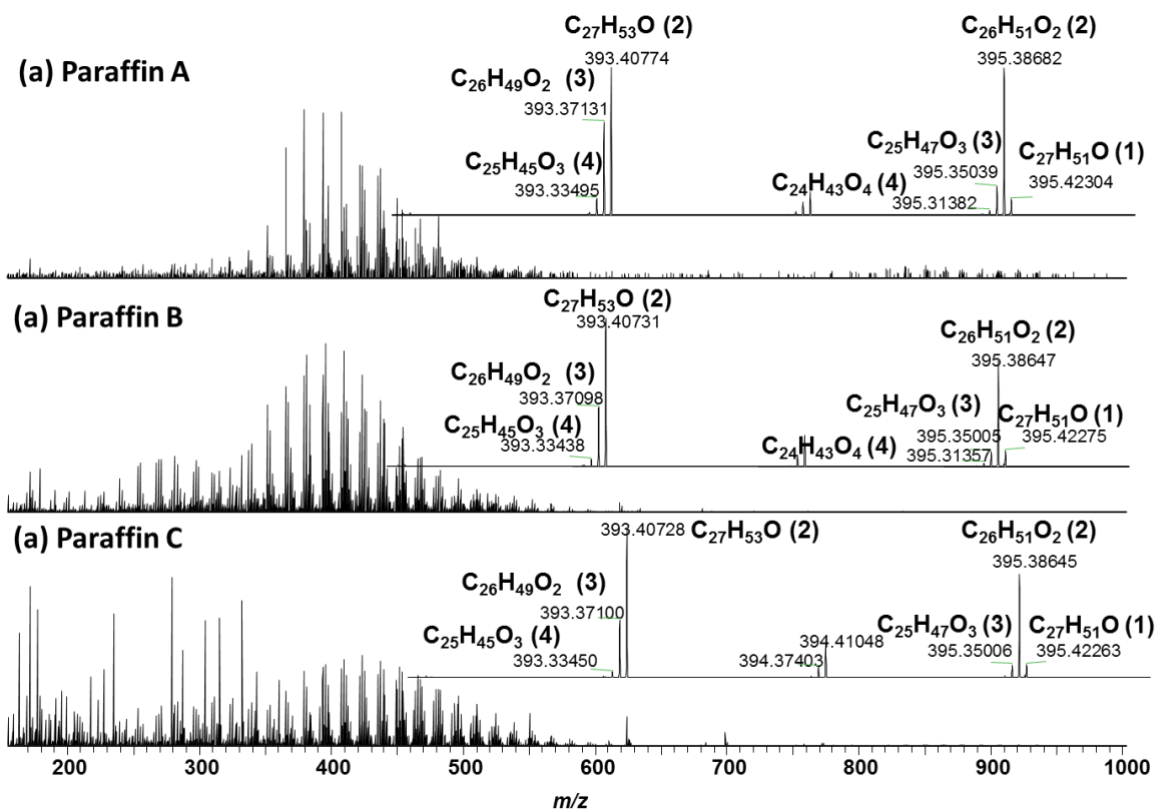

254

255 **Figure 9.** DART(+)-Orbitrap mass spectra for paraffin samples. The right side inserts show  
 256 paraffin detection in oxygenated form ( $O_x$  classes). The number between parentheses  
 257 corresponds to DBE value.

258

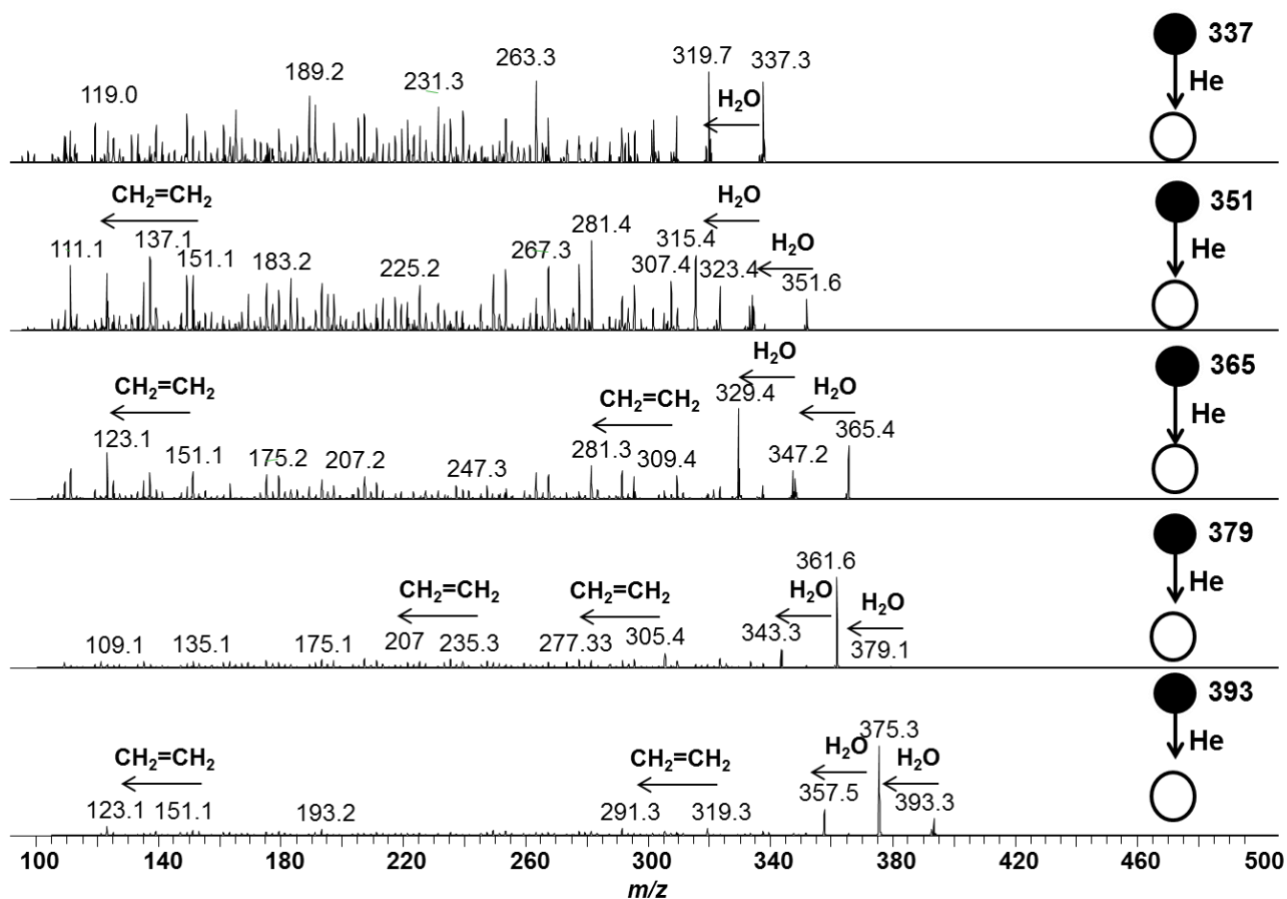

259


260 **Figure 10.** DART(+)-MS/MS for ions with  $m/z$  337, 351, 365, 379 and 393.

261

### 262 3.4 Porphrin compounds analyzes

263 The detection ability of DART(+)-MS technique regarding standard porphyrin and  
 264 metal porphyrin compounds were evaluated, **Figure 11** and **12**, respectively. The **Figure 11a**  
 265 shows the DART(+)-Orbitrap MS for the 5,10,15,20-tetraphenyl-21H,23H-porphine  
 266 compound, detected as  $[\text{M} + \text{H}]^+$  ion with  $m/z$  615.2502, where  $\text{M} = \text{C}_{44}\text{H}_{30}\text{N}_4$ . Its chemical  
 267 structure was confirmed from CID experiments, **Figure 11b**, in which the two neutral losses  
 268 of 77 Da and one of 16 Da identifies the presence of two phenyl rings and one amine group  
 269 ( $\text{C}_6\text{H}_5$  and  $\text{NH}_2$ ). The standard metal porphyrin compound was also detected by DART(+)-


270 Orbitrap MS, with lower ionization efficiency, as  $[M + H]^+$  ion with  $m/z$  680.1763, where M  
 271 =  $C_{44}H_{28}N_4O$ , **Figure 12a**. Its connectivity and its structure were confirmed from  
 272 DART(+)-MS<sup>2</sup> and DART(+)-MS<sup>3</sup> experiments, in which two neutral losses of 77 Da  
 273 happened simultaneously ( $m/z$  680  $\rightarrow$  603 and  $m/z$  603  $\rightarrow$  526 transitions **Figures 12b-c**).

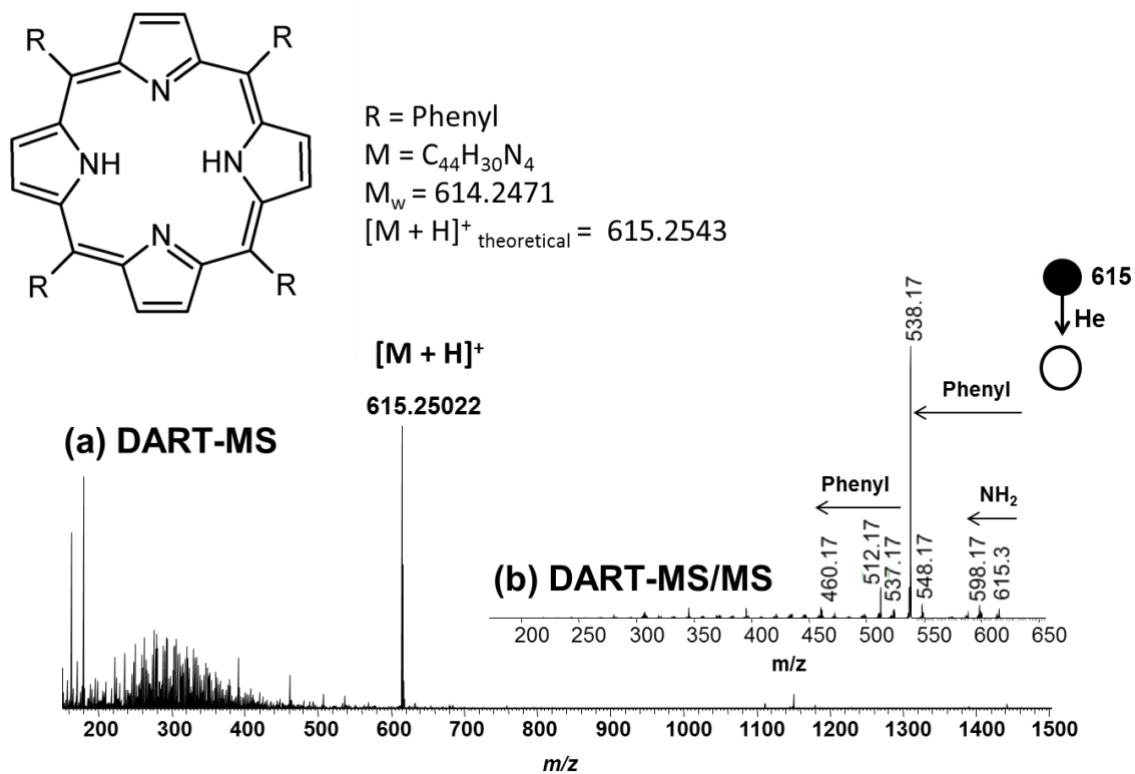

274

275 **Figure 11.** (a) DART(+) mass spectrum of porphyrin standard using He at 400 °C and (b)  
 276 DART(+)MS/MS for the ion with  $m/z$  615.

277

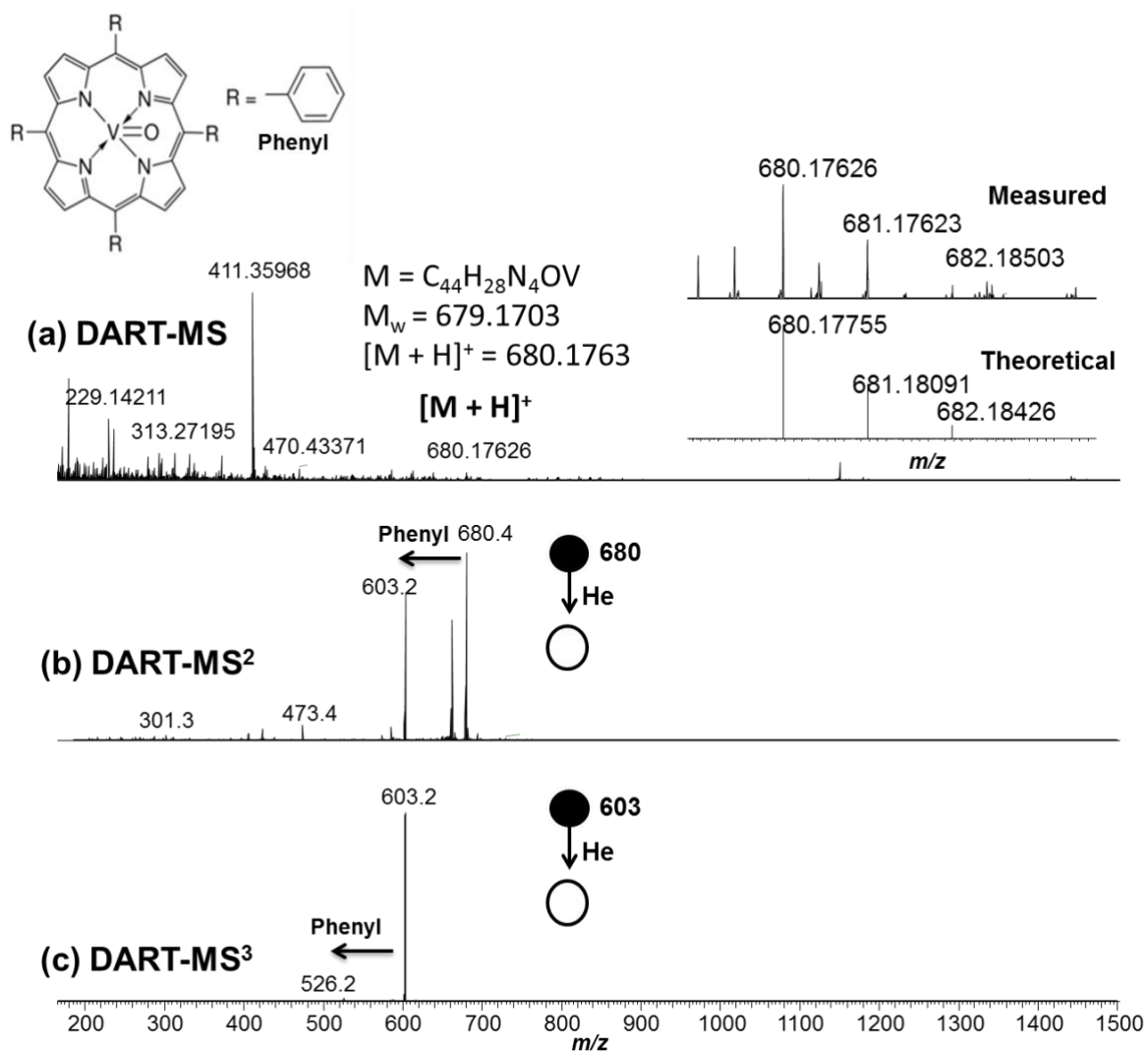

278

279 **Figure 12.** (a) DART(+) mass spectrum of vanadium porphyrin standard and (b)

280 DART(+)MS<sup>2</sup> and (c) DART(+)MS<sup>3</sup> for ions with  $m/z$  680 and 603, respectively.

281

282

283

284

285

286

#### 287 **4. Conclusion and Perspectives**

288           The DART(+)-hybrid ion-trap-Orbitrap MS is a powerful, simple, and easy  
289 analytical tool that can be applied to petroleum analysis to assess assessing chemical  
290 composition at molecular level. Cellulose-based substrates together with high temperature  
291 (400 °C), He as a gas heater and crude oil concentrations higher than 6 mg mL<sup>-1</sup> increased  
292 the sensibility of DART source. Nitrogen-containing species were detected as protonated  
293 molecules in crude oil samples, following by NO[H] and O[H] class species.

294           DART(+) also rendered hydroxylated HCs species (O<sub>x</sub> classes, where x ranged from  
295 1 to 4) with DBEs of 1 to 4 for paraffin samples. Oxidation reactions occur at high  
296 temperature and atmospheric pressure between HCs species and generated short-lifetime  
297 oxygen-based species such as OH, OOH radicals and H<sub>3</sub>O<sup>+</sup>. DART(+)-MS and CID  
298 experiments (MS<sup>2</sup> and MS<sup>3</sup>) were also able to identify porphyrin standard compounds as [M  
299 + H]<sup>+</sup> ions with *m/z* 615.2502 and 680.1763, where M = C<sub>44</sub>H<sub>30</sub>N<sub>4</sub> and C<sub>44</sub>H<sub>28</sub>N<sub>4</sub>OV,  
300 respectively.

#### 301 **5. Acknowledgments**

302           The authors thank FAPES, FAPEG, PETROBRAS, CNPq, and CAPES for their financial  
303 support. The financial support of the Conseil Régional d'Aquitaine (20071303002PFM) and  
304 FEDER (31486/08011464) is acknowledged.

305

306

---

## *References*

- <sup>1</sup> R. B. Cody, J. A. Laramée, H. D. Durst, Versatile New Ion Source for the Analysis of Materials in Open Air under Ambient Conditions, *Anal. Chem.* 77 (2005) 2297-03.
- <sup>2</sup> J. A. Laramée, Cody RB, Nilles JM, Durst HD (2007) In: *Forensic Analysis on the Cutting Edge*, Blackledge, RD, Ed. Forensic Applications of DART (Direct Analysis in Real Time) Mass Spectrometry. John Wiley and Sons, NJ.
- <sup>3</sup> R. W. Jones, R. B. Cody, J. F. McClelland, Differentiating writing inks using direct analysis in real time mass spectrometry, *J. Foren. Sci.* 51(4) (2006) 915-18.
- <sup>4</sup> R. W. Jones, J. F. McClelland, Analysis of writing inks on paper using direct analysis in real time mass spectrometry, *Foren. Sci. Int.* 231 (2013) 73–81.
- <sup>5</sup> R. A. Musah, R. B. Cody, M. A. Domin, A. D. Lesiak, A. J. Dane, J. R.E. Shepard, DART–MS in-source collision induced dissociation and high mass accuracy for new psychoactive substance determinations. *Foren. Sci. Int.* 244 (2014) 42–9.
- <sup>6</sup> A. D. Lesiak, R. B. Cody, A. J. Dane, R. A. Musah, Rapid detection by direct analysis in real time-mass spectrometry (DART-MS) of psychoactive plant drugs of abuse: The case of *Mitragyna speciosa* aka “Kratom”, *Foren. Sci. Int.* 242 (2014) 210–8.
- <sup>7</sup> F. M. Fernandez, R. B. Cody, M. D. Green, C. Y. Hampton, R. McGready, S. Sengaloundeth, N. J. White, P. N. Newton, Characterization of solid counterfeit drug samples by desorption electrospray ionization and direct-analysis-in-real-time coupled to time-of-flight mass spectrometry, *Chem. Med Chem.* 1(7) (2006) 702-05.

- 
- <sup>8</sup> R-M. Räsänen, P. Dwivedi, F. M. Fernández, T. J. Kauppila, Desorption atmospheric pressure photoionization and direct analysis in real time coupled with travelling wave ion mobility mass spectrometry, *Rapid Commun. Mass Spectrom.* 28 (2014) 2325–36.
- <sup>9</sup> Y. Zhao, M. Lam, D. Wu, R. Mak, Quantification of small molecules in plasma with direct analysis in real time tandem mass spectrometry, without sample preparation and liquid chromatographic separation, *Rapid Commun. Mass Spectrom.* 22 (2008) 3217–24.
- <sup>10</sup> O. P. Haefliger, N. Jeckelman, Direct mass spectrometric analysis of flavors and fragrances in real applications using DART, *Rapid Commun. Mass Spectrom.* 21 (2007) 1361-66.
- <sup>11</sup> T. Vail, P. R. Jones, O. D. Sparkman, Rapid and unambiguous identification of melamine in contaminated pet food based on mass spectrometry with four degrees of confirmation, *J Anal Toxicol* 31(6) (2007) 304-12.
- <sup>12</sup> L. Vaclavik, T. Cajka, V. Hrbek, J. Hajslova, Ambient mass spectrometry employing direct analysis in real time (DART) ion source for olive oil quality and authenticity assessment. *Anal. Chimica Acta* 645 (2009) 56–63.
- <sup>13</sup> V. Hrbek, L. Vaclavik, O. Elich, J. Hajslova, Authentication of milk and milk-based foods by direct analysis in real time ionization high resolution mass spectrometry (DARTeHRMS) technique: A critical assessment, *Food Control* 36 (2014) 138-45.
- <sup>14</sup> Lukas Vaclavik, Beverly Belkova, Zuzana Reblova, Katerina Riddellova, Jana Hajslova. Rapid monitoring of heat-accelerated reactions in vegetable oils using direct analysis in real time ionization coupled with high resolution mass spectrometry. *Food Chem.* 138 (2013) 2312–2320.
- <sup>15</sup> C. Y. Pierce, J. R. Barr, R. B. Cody, R. F. Massung, A. R. Woolfitt, H. Moura, H. A. Thompson, F. M. Fernandez, Ambient generation of fatty acid methyl ester ions from



---

bacterial whole cells by direct analysis in real time (DART) mass spectrometry, *Chem. Commun.* 28(8) (2007) 807-09.

<sup>16</sup> J. Y. Yew, R. B. Cody, E. A. Kravitz, Cuticular hydrocarbon analysis of an awake behaving fly using direct analysis in real-time time-of-flight mass spectrometry, *Proc. Nat. Acad. Sci. USA* 105(20) (2008) 7135-40.

<sup>17</sup> S. Yu, E. Crawford, J. Tice, B. Musselman, J-T. Wu, Bioanalysis without sample cleanup or chromatography: the evaluation and initial Implementation of Direct Analysis in Real Time Ionization Mass Spectrometry for the Quantification of Drugs in Biological Matrixes, *Anal. Chem.* 81 (2009) 193-02.

<sup>18</sup> G. Morlock, Y. Ueda, New coupling of planar chromatography with direct analysis in real time mass spectrometry, *J. Chromatogr. A* 1143 (2007) 243-51.

<sup>19</sup> R. B. Cody, Observation of Molecular Ions and Analysis of Nonpolar Compounds with the Direct Analysis in Real time Ion Source, *Anal. Chem.* 81 (2009) 1101-07.

<sup>20</sup> L. Song, A. B. Dykstra, H. Yao, J. E. Bartmess, Ionization Mechanism of Negative Ion-Direct Analysis in Real Time: A Comparative Study with Negative Ion-Atmospheric Pressure Photoionization, *J. Am. Soc. Mass Spectrom.* 20 (2009) 42-50.

<sup>21</sup> R. Haddad, T. Regiani, C. F. Klitzke, G. B. Sanvido, Y. E. Corilo, D. V. Augusti, V. M. D. Pasa, R. C. C. Pereira, W. Romão, B. G. Vaz, R. Augusti, M. N. Eberlin, Gasoline, Kerosene, and Diesel Fingerprinting via Polar Markers, *Energy Fuels* 26 (2012) 3542-7.

<sup>22</sup> J. L. Rummel, A. M. McKenna, A. G. Marshall, J. R. Eyler, D. H. Powell, The coupling of direct analysis in real time ionization to Fourier transform ion cyclotron resonance mass spectrometry for ultrahigh-resolution mass analysis, *Rapid Commun. Mass Spectrom.* 24 (2010) 784–90.

- 
- <sup>23</sup> G. P. Dalmaschio, M. M. Malacarne, V. M. D. L. Almeida, T. M. C. Pereira, A. O. Gomes, E. V. R. Castro, S. J. Greco, B. G. Vaz, W. Romão, Characterization of polar compounds in a true boiling point distillation system using electrospray ionization FT-ICR mass spectrometry, *Fuel*, 115 (2014) 190-202.
- <sup>24</sup> N. S. Tessarolo, R. C. Silva, G. Vanini, A. Pinho, W. Romão, E. V. R. Castro, D. A. Azevedo, Assessing the chemical composition of bio-oils using FT-ICR mass spectrometry and comprehensive two-dimensional gas chromatography with time-of-flight mass spectrometry. *Microchem. J.* 117 (2014) 68-76.
- <sup>25</sup> T. M. C. Pereira, G. Vanini, L. V. Tose, F. M. Cardoso, F. P. Fleming, P. T. V. Rosa, C. J. Thompson, E. V. R. Castro, B. G. Vaz, W. Romão, FT-ICR MS analysis of asphaltenes: Asphaltenes go in, fullerenes come out. *Fuel (Guildford)*, 131 (2014) 49-58.
- <sup>26</sup> L. A. Terra, P. R. Filgueiras, L. V. Tose, W. Romão, D. D. Souza, E. V. R. Castro, M. S. L. Oliveira, J. C. M. Dias, R. J. Poppi, Petroleomics by Electrospray Ionization FT-ICR Mass Spectrometry Coupled to Partial Least Squares with Variable Selection Methods: Prediction of the Total Acid Number of Crude Oils, *The Analyst* 139 (2014) 4908-16.
- <sup>27</sup> H. B. Costa, L. M. Souza, L. C. Soprani, B. G. Oliveira, E. M. Ogawa, A. M. N. Korres, J. A. Ventura, W. Romão, Monitoring the Physicochemical Degradation of Coconut Water Using ESI-FT-ICR MS, *Food Chem.* 174 (2014) 139-46.
- <sup>28</sup> C. F. Klitzke, Y. E. Corilo, K. Siek, J. Binkley, J. Patrick, M. N. Eberlin, Petroleomics by Ultrahigh-Resolution Time-of-Flight Mass Spectrometry, *Energy Fuels* 26 (2012) 5787-94.
- <sup>29</sup> M. Benassi, A. Berisha, W. Romão, E. Babayev, A. Rompp, B. Spengler, Petroleum crude oil analysis using low-temperature plasma mass spectrometry. *RCM. Rapid Commun. Mass Spectrom.* 27 (2013) 825-834.

- 
- <sup>30</sup> K. A. P. Colati, G. P. Dalmaschio, E. V. R. de Castro, A. O. Gomes, B. G. Vaz, W. Romão, Monitoring the liquid/liquid extraction of naphthenic acids in brazilian crude oil using electrospray ionization FT-ICR mass spectrometry (ESI FT-ICR MS), *Fuel* 108 (2013) 647-55.
- <sup>31</sup> L. V. Tose, F. M. R. Cardoso, F. P. Fleming, M. A. Vicente, S. R. C. Silva, E. V. R. Castro, G. M. F. V. Aquije, B. G. Vaz, W. Romão, Analyzes of Hydrocarbons by Atmosphere Pressure Chemical Ionization FT-ICR Mass Spectrometry using Isooctane as ionizing Reagent, *Fuel*, 2015, *in press*.
- <sup>32</sup> M. Schiorlin, E. Marotta, M. D. Molin, C. Paradisi, *Environ. Sci. Technol.* 47 (2013) 542–548.
- <sup>33</sup> C Wu, K Qian, M. Nefliu, R. G. Cooks, Ambient analysis of saturated hydrocarbons using discharge-induced oxidation in desorption electrospray ionization. *J. Am Soc Mass Spectrom.* 21 (2010) 261-267.