



HAL
open science

Detecting the elemental and molecular signatures of life: Laser-based mass spectrometry technologies

Niels Ligterink, Andreas Riedo, Marek Tulej, Rustam Lukmanov, Valentine Grimaudo, Coenraad De Koning, Peter Wurz, Christelle Briois, Nathalie Carrasco, France Ricardo Arevalo, et al.

► **To cite this version:**

Niels Ligterink, Andreas Riedo, Marek Tulej, Rustam Lukmanov, Valentine Grimaudo, et al.. Detecting the elemental and molecular signatures of life: Laser-based mass spectrometry technologies. *Bulletin of the American Astronomical Society*, 2021, 53 (4), <10.3847/25c2cfef.782ec4d0>. <hal-03991342>

HAL Id: hal-03991342

<https://hal.science/hal-03991342v1>

Submitted on 6 Jun 2023

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons CC BY 4.0 - Attribution - International License

Detecting the elemental and molecular signatures of life: Laser-based mass spectrometry technologies

A whitepaper for

The Planetary Science and Astrobiology Decadal Survey 2023–2032
National Academy of Sciences

Niels F.W. Ligterink^{1,*}, Andreas Riedo², Marek Tulej², Rustam Lukmanov², Valentine Grimaudo², Coenraad de Koning², Peter Wurz², Christelle Briois³, Nathalie Carrasco⁴, Ricardo Arevalo Jr.⁵, and William B. Brinckerhoff⁶

¹Center for Space and Habitability, University of Bern, Switzerland
(*niels.ligterink@csh.unibe.ch, [+41 31 631 4415](tel:+41316314415))

²Space Research and Planetary Sciences, Physics Institute, University of Bern, Switzerland

³Laboratoire de Physique et Chimie de l'Environnement et de l'Espace (LPC2E), Orléans, France

⁴Laboratoire Atmosphères, Milieux et Observations Spatiales (LATMOS), Guyancourt, France

⁵Department of Geology, University of Maryland, USA

⁶NASA Goddard Space Flight Center, USA

List of co-signatories: [LINK](#)

Recommendation: This whitepaper highlights the advantages that laser-based mass spectrometry offers for the *in situ* detection of chemical biosignatures on planets and moons. Future planetary lander and rover missions should be encouraged to make use of laser-based mass spectrometry instruments in their scientific payload.

Introduction

The detection and identification of signatures of life, extinct or extant, on Solar System bodies other than Earth, is nowadays one of the most challenging endeavors in space research. These signatures can be abundances of specific biogenic elements, fractionated isotopes due to metabolic processes, or biomolecules, such as amino acids, nucleobases, and more complex macromolecular networks¹. Past *in situ* space exploration missions have attempted to identify biosignatures, albeit with limited success. Some problems have been pointed out with traditional planetary exploration instruments, such as insufficient sensitivity, the possibility that faint biosignatures in micrometer-sized fossils get lost in bulk analysis techniques, or the chemical alteration of biosignatures during, e.g., sample processing. To counter these issues, a new type of analysis technique is required, which combines high detection sensitivity and spatial resolution measurement capabilities at micrometer level and below, and does not chemically alter the sample material of interest. An analysis technique that combines all these requirements is laser-based mass spectrometry, including Laser ablation Ionization and Laser Desorption Mass Spectrometry (LIMS, LD-MS). LIMS and LD-MS instruments have the advantage that they do not require carrier gasses, energy-consuming ovens, or extensive sample preparation steps. This sets them apart from to date traditional planetary exploration instruments, such as pyrolysis Gas Chromatography – Mass Spectrometry (pyr GC – MS) or other chromatographic methods involving extraction and preparation of the sample by wet or gaseous chemistry, biasing the analysis to certain chemical groups, before mass spectrometric analysis. It also means that LIMS and LD-MS instruments can be significantly more compact, lower in weight and power consumption, and are less complex. In this whitepaper, we give an overview of prototype LIMS and LD-MS systems for space exploration, their versatility, and their beneficial characteristics for future deployment for the detection of life. A Mars and a Europa scenario are used to show how the instruments that are developed for use on planetary science missions, can operate on rocky planets or ice-covered ocean worlds, respectively, and search for signatures of life. Finally, we will discuss other uses of LIMS and LD-MS instruments on future space exploration missions.

Laser-based mass spectrometry measurement techniques, LIMS and LD-MS

In its essence, all LIMS and LD-MS instruments operate following the same principle. A single laser pulse is fired at a solid material, where the laser-matter interaction results in the release of atoms or molecules to the gas-phase. During the same laser-matter interaction a fraction of the species is ionized, which makes it possible to detect and identify them using a mass analyzer. The main difference between LIMS and LD-MS is their applied laser irradiance ($W\text{ cm}^{-2}$), which induces material ablation or desorption and ionization, respectively. By tuning the laser irradiance, which in turn depends on

parameters such as laser pulse length, laser wavelength, focal spot size, and pulse energy, different laser-matter interactions can be achieved, which can be grouped in both laser ablation and laser desorption conditions outlined in the following in more detail.

Laser Desorption Mass Spectrometry (LD-MS). At moderate laser power densities (MW to GW/cm²) this interaction results in laser desorption; the breaking of physisorbed bonds of molecules to a surface, followed by the intact release of these molecules to the gas phase. Ionization is either realized directly at the desorption of molecules from the surface or via a molecule-laser interaction in the gas-phase. Typically, nanosecond laser systems are preferred for LD-MS, as, in comparison to short-pulsed femtosecond laser systems, the photon density is significantly lower, which induces only modest fragmentation of the molecules desorbed from the surface; a severe fragmentation would challenge their identification. Laser desorption has a number of advantages: in contrast to pyrolysis, it does not chemically alter molecules or induce reactions between organic molecules and salts or minerals, although compound molecules may form. At mild desorption conditions near the laser desorption threshold, molecular fragmentation is limited, which prevents isobaric interference and contrasts with, e.g., electron ionization, which results in extensive fragmentation.

Laser Ablation Ionization Mass Spectrometry (LIMS). At high laser irradiances (GW/cm² to TW/cm² level), the laser-matter interaction results in laser ablation conditions; a volatile matter-plasma conversion, accompanied by the near-complete atomization and, to a large extent, ionization of the material. Laser ablation is particularly well suited to investigate the element and isotope composition of solids of interest. By applying multiple laser pulses to the same surface spot, layers of material can be removed shot by shot, and subsurface material can be analyzed (chemical depth profiling). In particular femtosecond laser pulses at UV wavelengths are well suited for this because of efficient ionization, close to stoichiometric ion production, and minimized element fractionation effects. With well-focused laser beams, focal spot sizes of several micrometers can be achieved, which makes it possible to perform chemical analysis of micrometer-sized features (microfossils, finely laminated biopatterns, etc). The application of chemical imaging and depth profiles in close vicinity to each other allows, e.g., spatially resolved chemical analysis of the solids in three dimensions, and, e.g., to identify regions of interest. The laser ablation process can also help by the identification of intact neutral molecules when neutral molecules form clusters with metal ions in the plasma.

Subsequent analysis of the ions formed in the laser desorption or ablation process is performed with a mass analyzer. For *in situ* space applications, several miniature mass analyzers have been developed, such as Quadrupole Mass Spectrometers (QMS)^{2,3}, high mass resolution Orbitrap systems⁴⁻⁶, and time-of-flight Mass Spectrometers^{7,8} (TOF MS), some of which are optimized for laser ablation and laser desorption studies. Although in principle not necessary for LIMS and LD-MS systems, the mass analyzer can be outfitted with a post ionization stage, such as electron ionization or a secondary laser ionization

technique^{9,10} to ionize neutral species and significantly increase the systems' sensitivity, however, at the expense of higher complexity of the overall system (e.g., more optical elements).

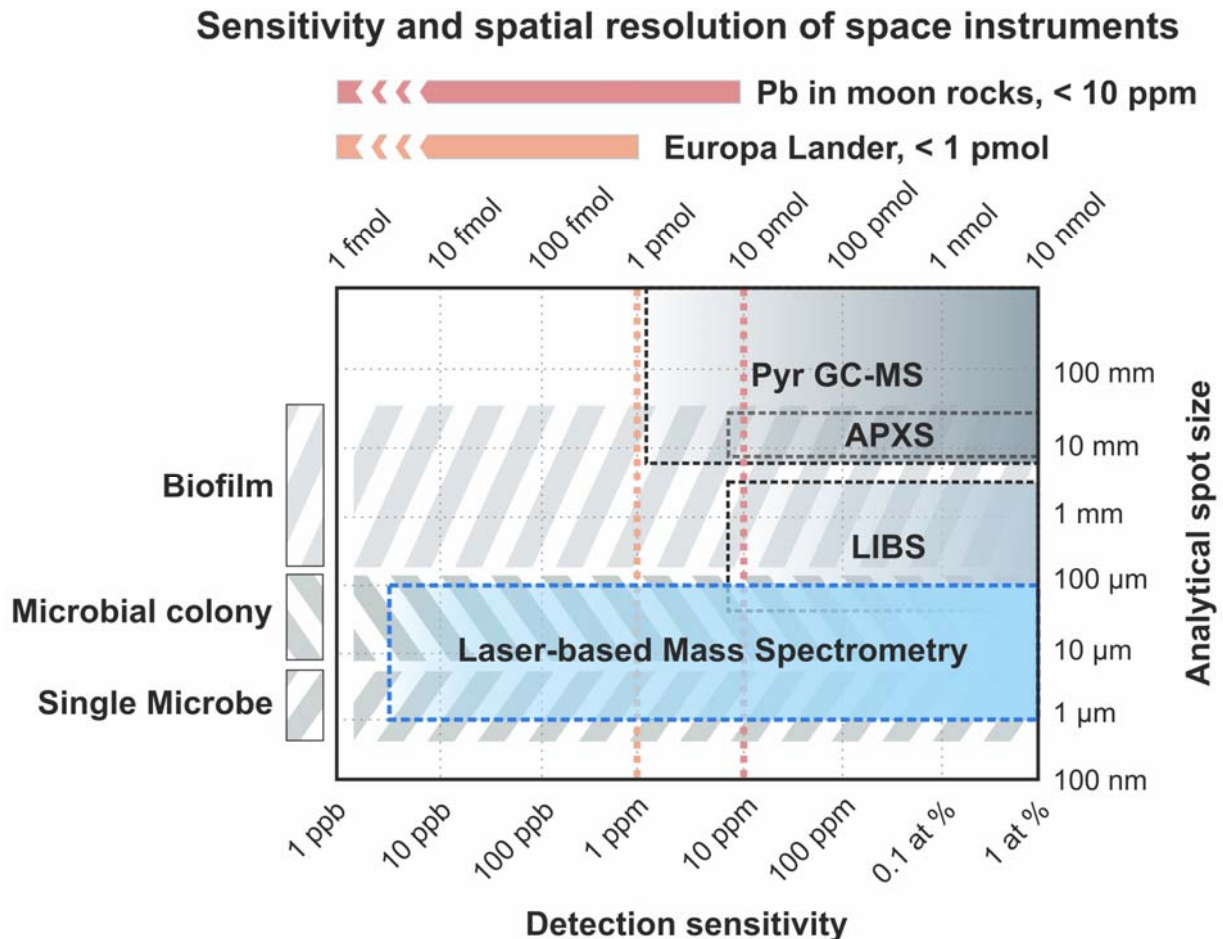


Figure 1. Sensitivity and spatial resolution of elemental and chemical analysis instruments developed for in situ space exploration missions. Laser-based Mass Spectrometry techniques (LIMS, LD-MS) are compared with Laser-Induced Breakdown Spectroscopy (LIBS), Alpha Particle X-ray Spectrometry (APXS), and pyrolysis Gas Chromatography – Mass Spectrometry (pyr GC – MS). Typical sizes of biological units are indicated, as well as quantities and concentrations required by space exploration missions or for element analysis. Note that the top and bottom sensitivity axis are related by the total sample mass. In both sensitivity and spatial resolution, laser-based instrumentation performs better than current space instruments.

The Mars scenario

While evidence points in the direction of a watery past for Mars, nowadays it is a mostly barren place, which has lost most of its liquid water. Therefore, the expectation is to find extinct life on Mars, most probably in the subsurface that protects the fragile signatures from the harsh ionizing environment. Recently, the most promising biosignatures were

reviewed in detail and allocated to six groups, ranging from biomolecules to microstructures, to solids where specific elements are isotope fractionated (e.g., sulfur)¹.

The high detection sensitivity of laser-based mass spectrometers (see Fig. 1 for a comparison of experimental techniques) coupled to spatially resolved chemical analysis, allows for the identification and chemical characterization of small scale features at the micrometer level. These advantages have now been recognized by the space exploration community and resulted in the inclusion of LD-MS in the Mars Organic Molecule Analyzer (MOMA) instrument on the ESA/Roscosmos ExoMars rover^{2,3}.

Recent measurement campaigns demonstrate that LIMS systems not only have the measurement capabilities to analyze the elemental composition of minerals¹¹, lunar regolith¹², and meteorites¹³ but can also identify putative fossil structures of micrometer dimensions embedded in a host material by monitoring elements relevant to life, such as carbon^{14,15}. Measurements conducted on Martian mudstone analog material allowed the identification of the elemental signatures of single microbes that were inoculated artificially at low number density¹⁶. Moreover, LIMS systems, such as LAb-CosmOrbitrap⁴, have been able to successfully characterize the chemical composition of organics and tholins mixed with minerals^{17,18}, making use of metal-organic clusters for identification and the high mass resolving power of Orbitrap to disentangle isobaric interferences of possible biogenic and non-biogenic molecular fragments. In combination with a microscope camera system, features of micrometer dimensions can be localized and targeted with the laser beam with high accuracy, which allows for the chemical investigation of the feature with minimal contribution of the host material that may obscure the faint life signature. Chemical imaging further refines the analysis of these microfossils by scanning the entire region of interest.

Sample preparation is limited for LIMS measurements on Mars. In a best-case scenario, the LIMS system can be placed on a robotic arm to make chemical measurements on the local environment. Alternatively, the LIMS system can be placed on the payload platform on a lander or rover and have samples introduced to it via a sample delivery system. To conduct a thorough analysis, it needs to be possible to manipulate samples in the x,y,z direction and potentially rotate them. Note that the two laser-based mass spectrometry systems that to date have been prepared for space flight, LAZMA¹² for the Roscosmos Phobos-Grunt (Russian: Phobos-Soil) mission, and MOMA³ for the ESA/Roscosmos ExoMars rover, both have a sample carousel that allows for successive analyses of an aliquot of the sample across the arc of the carousel motion.

The Europa scenario

While Europa and other ocean worlds have inhospitable ice shells, they do contain subsurface oceans of liquid water and in some cases conditions that are considered right for life to emerge and persist¹⁹. Through cracks in the ice crust or with plumes, signatures of this life can end up on its surface, form deposits, and become accessible for *in situ* analysis. However, the transfer of molecular biosignatures to the surface is likely inefficient and at the surface, these molecules can be processed by radiation. Instruments that

search for these signatures thus need to be able to detect trace amounts of biomolecules and identify associations between detected organic compounds and geological phases that may shield them from harsh environmental conditions. LIMS and LD-MS instruments commonly reach picomol mm⁻² detection sensitivity^{3,9,18,20-24}, and in some cases even achieve detection limits down to several femtomol mm⁻², such as with the LD-MS system called ORIGIN²⁵. Space LIMS and LD-MS systems have demonstrated that they can detect several classes of biomolecules, such as amino acids and nucleobases, even in complex mixtures, or when mixed with minerals and salts^{4,9,18,20,21,25,26}. No signs for the chemical alteration of biomolecules by laser-induced reactions between minerals and organic molecules are seen during LIMS or LD-MS analysis of molecules mixed with minerals or salts.

LIMS and LD-MS systems can perform measurements directly on ice samples, but this may not be the most efficient way to detect biomolecules. It is better to make use of the physical characteristics of Europa and its ice constituents, by simply warming up an ice sample until the water ice sublimates and a high concentration residue of biomolecules, salts, and minerals remain, which can be considered as a simple enrichment process. As shown before, these residues do not have to be filtered since the LIMS/LD-MS measurement does not chemically alter its components. Sample preparation is therefore limited and simple and only consists of warming up an ice sample after its introduction into an instrument, which is for free by using the thermal dissipation of the instrument payload²⁵.

Other LIMS/LD-MS applications on space exploration missions

Due to a large number of technical benefits and its measurement capabilities, LIMS and LD-MS instruments are an excellent addition to the scientific payloads of planetary exploration missions when searching for signatures of life. However, the unique properties of LIMS and LD-MS systems also make other uses possible. Since these systems are compact and lightweight they are suited for miniature landers or handheld applications as well. The Hayabusa-2 mission to asteroid Ryugu deployed multiple small landers, such as MASCOT, to investigate the asteroid *in situ*²⁷. MASCOT weighs around 10 kg and has a volume of 15,500 cm³. This type of lander can therefore easily fit a space LIMS/LD-MS instrument and be deployed to perform *in situ* chemical analysis of asteroids or comets. Due to their compact and lightweight nature, future space exploration missions could carry multiple miniature landers with LIMS/LD-MS systems and investigate multiple science targets in a single mission. Similarly, a handheld LIMS system can be used during future Moon landings to perform easy *in situ* analysis of Moon regolith. On sample return missions, LIMS and LD-MS can be employed to perform preliminary chemical analysis of sample sites, so the scientifically most interesting samples are collected and brought to Earth.

Concluding remarks

Laser-based mass spectrometer systems for use on space exploration missions have been developed and tested for over two decades, demonstrating the maturity of this

technology. Due to the versatility of LIMS and LD-MS systems, they can be applied in a wide variety of applications where *in situ* chemical analysis is desired and are in particular well suited to find chemical indicators of life. Because chemical processing of samples is not required, LIMS/LD-MS instruments are compact, require no consumables, and above all are not biased to specific chemical groups. Chemical imaging and depth profiling on micrometer scales makes it possible to detect the chemical biosignatures of microfossils. A number of concept planetary exploration missions, such as the Europa Lander, Lunar landers, future Mars missions, and various sample return missions, can benefit from including a Laser Desorption/Ablation Mass Spectrometer to meet their scientific goals.

- 1 Hays, L. E. *et al.* Biosignature Preservation and Detection in Mars Analog Environments. *Astrobiology* **17**, 363-400 (2017).
- 2 Li, X. *et al.* Mars Organic Molecule Analyzer (MOMA) laser desorption/ionization source design and performance characterization. *International Journal of Mass Spectrometry* **422**, 177-187 (2017).
- 3 Goesmann, F. *et al.* The Mars Organic Molecule Analyzer (MOMA) instrument: characterization of organic material in martian sediments. *Astrobiology* **17**, 655-685 (2017).
- 4 Briois, C. *et al.* Orbitrap mass analyser for in situ characterisation of planetary environments: Performance evaluation of a laboratory prototype. *Planetary and Space Science* **131**, 33-45 (2016).
- 5 Arevalo, R. *et al.* in *2019 Astrobiology Science Conference*. (AGU).
- 6 Willhite, L. *et al.* Characterization of Regolith and Trace Economic Resources (CRATER) Instrument: Integration into the Artemis Program. *LPICo* **2241**, 5104 (2020).
- 7 Rohner, U. *et al.* A miniature laser ablation time-of-flight mass spectrometer for in situ planetary exploration. *Measurement Science and Technology* **14**, 2159 (2003).
- 8 Rohner, U. *et al.* Highly miniaturized laser ablation time-of-flight mass spectrometer for a planetary rover. *Review of Scientific Instruments* **75**, 1314-1322 (2004).
- 9 Getty, S. A. *et al.* Compact two-step laser time-of-flight mass spectrometer for in situ analyses of aromatic organics on planetary missions. *Rapid Communications in Mass Spectrometry* **26**, 2786-2790 (2012).
- 10 Anderson, F. S. *et al.* in *2012 IEEE Aerospace Conference*. 1-18 (IEEE).
- 11 Neubeck, A. *et al.* Mineralogical determination in situ of a highly heterogeneous material using a miniaturized laser ablation mass spectrometer with high spatial resolution. *International Journal of Astrobiology* **15**, 133-146 (2015).
- 12 Managadze, G. *et al.* Study of the main geochemical characteristics of Phobos' regolith using laser time-of-flight mass spectrometry. *Solar System Research* **44**, 376-384 (2010).
- 13 Frey, S. *et al.* Chemical analysis of a lunar meteorite by laser ablation mass spectrometry. *Planetary and space science* **182**, 104816 (2020).
- 14 Tulej, M. *et al.* Chemical Composition of Micrometer-Sized Filaments in an Aragonite Host by a Miniature Laser Ablation/Ionization Mass Spectrometer. *Astrobiology* **15**, 669-682 (2015).
- 15 Wiesendanger, R. *et al.* Chemical and optical identification of micrometre sized 1.9 billion-year-old fossils by combining a miniature LIMS system with an optical microscope. (2018).
- 16 Stevens, A. H. *et al.* Detectability of biosignatures in a low-biomass simulation of martian sediments. *Scientific Reports* **9**, 9706 (2019).
- 17 Selliez, L. *et al.* High-resolution mass spectrometry for future space missions: Comparative analysis of complex organic matter with LAb-CosmOrbitrap and laser desorption/ionization

- Fourier transform ion cyclotron resonance. *Rapid Communications in Mass Spectrometry* **34**, e8645 (2020).
- 18 Arevalo Jr, R. *et al.* An Orbitrap-based laser desorption/ablation mass spectrometer designed for
spaceflight. *Rapid Communications in Mass Spectrometry* **32**, 1875-1886 (2018).
- 19 Lunine, J. I. Ocean worlds exploration. *Acta Astronautica* **131**, 123-130 (2017).
- 20 Uckert, K. *et al.* IR resonance-enhanced organic detection with two-step laser desorption time-
of-flight mass spectrometry. *Icarus* **299**, 15-21 (2018).
- 21 Selliez, L. *et al.* Identification of organic molecules with a laboratory prototype based on the
Laser Ablation-CosmOrbitrap. *Planetary and Space Science* **170**, 42-51 (2019).
- 22 Riedo, A. *et al.* Highly accurate isotope composition measurements by a miniature laser ablation
mass spectrometer designed for in situ investigations on planetary surfaces. *Planetary and Space
Science* **87**, 1-13 (2013).
- 23 Riedo, A. *et al.* Performance evaluation of a miniature laser ablation time-of-flight mass
spectrometer designed for in situ investigations in planetary space research. *Journal of Mass
Spectrometry* **48**, 1-15 (2013).
- 24 Riedo, A. *et al.* Coupling of LMS with a fs-laser ablation ion source: elemental and isotope
composition measurements. *Journal of Analytical Atomic Spectrometry* **28**, 1256-1269 (2013).
- 25 Ligterink, N.F.W. *et al.* ORIGIN: a novel and compact Laser Desorption–Mass Spectrometry
system for sensitive in situ detection of amino acids on extraterrestrial surfaces. *Scientific reports*
10, 1-10 (2020).
- 26 Moreno-García, P. *et al.* Towards matrix-free femtosecond-laser desorption mass spectrometry
for in situ space research. *Rapid Communications in Mass Spectrometry* **30**, 1031-1036 (2016).
- 27 Ho, T.-M. *et al.* MASCOT—the mobile asteroid surface scout onboard the HAYABUSA2 mission.
Space Science Reviews **208**, 339-374 (2017).