

## Application of a Widely Tunable Near-Infrared Laser Instrument for Stable Isotope Ratio Measurements

Julien Cousin, Stéphane Plus, Weidong Chen, Éric Fertein, Daniel Boucher

## ▶ To cite this version:

Julien Cousin, Stéphane Plus, Weidong Chen, Éric Fertein, Daniel Boucher. Application of a Widely Tunable Near-Infrared Laser Instrument for Stable Isotope Ratio Measurements. European Geosciences Union General Assembly 2007 (EGU 2007), Apr 2007, Vienne (Austria), Austria. pp.10773. hal-04026187

HAL Id: hal-04026187

https://hal.science/hal-04026187


Submitted on 13 Mar 2023

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

Geophysical Research Abstracts, Vol. 9, 10773, 2007

SRef-ID: 1607-7962/gra/EGU2007-A-10773

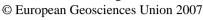



## 1 Application of a Widely Tunable Near-Infrared Laser Instrument for Stable Isotope Ratio Measurements

## J. COUSIN, S. PLUS, W. CHEN, E. FERTEIN, D. BOUCHER

Laboratoire de PhysicoChimie de l'Atmosphère, CNRS UMR 8101, Université du Littoral Côte d'Opale, 189 A, Av. Maurice Schumann, 59140 Dunkerque, France (chen@univ-littoral.fr)

In this paper, we report on the development and the application of a widely tunable near-infrared laser spectrometer for trace gas and isotope analysis. The spectrometer was based on fiber-coupled continuous-wave (cw) Telecom external cavity laser (ECDL, Tunics Plus) that is continuously tunable from 1500 to 1640 nm (C and L band) with an output power up to 3 mW and a tuning resolution of 0.001 nm ( $\sim 4\times 10^{-3}~\rm cm^{-1}$ ). The effective laser linewidth is less than 1 MHz. A multipass cell (New Focus – model 5612) in Herriott configuration with an optical path of 100 m was used to enhance the detection sensitivity.

The developed instrument has been used for isotopic composition analysis in woodbased combustion emission and in human breath. In a wood-based combustion, the measured  $^{13}\text{C}/^{12}\text{C}$  isotope ratio in CO<sub>2</sub> emission is found to be (1.1011±0.0040) % for the full burn operation regime. The corresponding  $\delta$ -value relative to PDB standard is (-20.17±3.53) %, that is in good agreement with the typical value of (-25±2) %, for wood.