



**HAL**  
open science

## Reflection Zone Plate Spectrometer For Li K Emission Detection In Lithium Compounds

Khalil Hassebi, Nicolas Rividi, Anne Verlaguet, Alexei Erko, Juergen Probst,  
Heike Loechel, Thomas Krist, Karine Le Guen, Philippe Jonnard

► **To cite this version:**

Khalil Hassebi, Nicolas Rividi, Anne Verlaguet, Alexei Erko, Juergen Probst, et al.. Reflection Zone Plate Spectrometer For Li K Emission Detection In Lithium Compounds. International Workshop on the Characterisation and quantification of lithium, Jun 2023, Paris, France. . hal-04056076

**HAL Id: hal-04056076**

**<https://hal.science/hal-04056076>**

Submitted on 20 Jun 2023

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution 4.0 International License

# REFLECTION ZONE PLATE SPECTROMETER FOR LI K EMISSION DETECTION IN LITHIUM COMPOUNDS

Khalil Hassebi<sup>1</sup>, Nicolas Rividi<sup>2</sup>, Anne Verlaguet<sup>3</sup>, Alexei Erko<sup>†4</sup>, Juergen Probst<sup>5</sup>, Heike Loechel<sup>5</sup>, Thomas Krist<sup>5</sup>, Karine Le Guen<sup>1</sup>, Philippe Jonnard<sup>1</sup>



1 Laboratoire de Chimie Physique—Matière et Rayonnement, Faculté des Sciences et Ingénierie, Sorbonne Université, UMR CNRS, 4 Place Jussieu, 75252 Paris Cedex 05, France  
 2 Service Camparis, UMS 7154—CNRS, OSU Ecce Terra, Sorbonne Université, 4 place Jussieu, 75005 Paris, France  
 3 Sorbonne Université, CNRS-INSU, Institut des Sciences de la Terre de Paris, ISTeP UMR 7193, F 75005 Paris, France  
 4 Institut für Angewandte Photonik e. V., Rudower Chaussee 29-31, D-12489 Berlin, Germany  
 5 Nano Optics Berlin GmbH, Krumme Str 64, D-10627 Berlin, Germany



## Motivation:

### What?


Unlock new frontier by allowing the **quantification of lithium** in complex materials such as lithium bearing minerals and materials used for **lithium batteries** with a technique that is both accessible and cost-effective.

### The challenges:

- Low fluorescence yield
  - Low energy emission and
  - Strong self-absorption
- Make the detection with a **non destructive technique** an unresolved issue.

### How?

- Probing the fine electronic structure of the valence band using soft x-ray emission spectroscopy technique (SXES) integrated in EPMA
- Tracking the changes in the chemical environment and emission spectrum (Intensity, shape, shift, band width...)

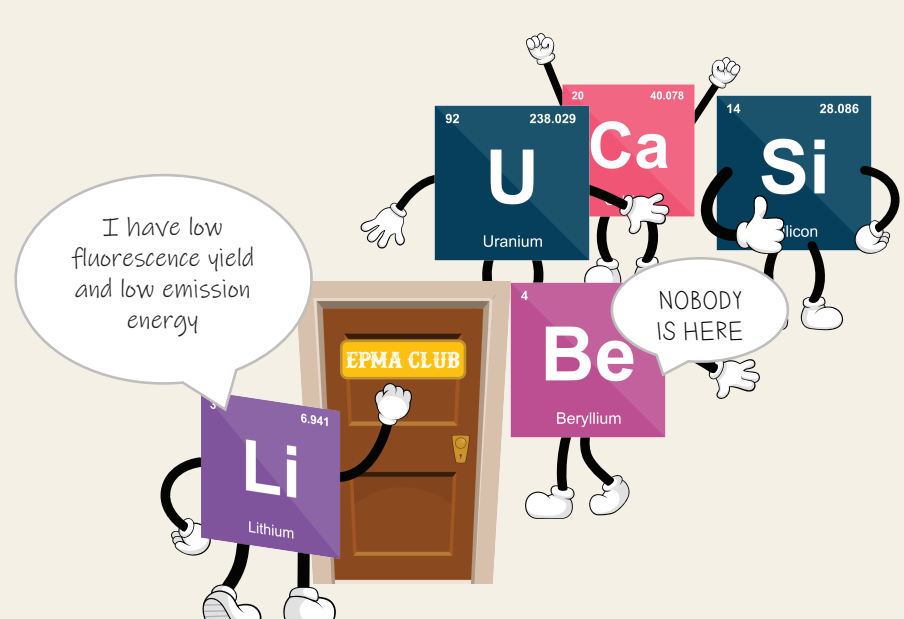


Can we detect lithium and study its electronic structure in complex materials?

Yes!

## Results:

The Li K emission occurs when a **2p** valence states transition to **1s** core states in lithium atoms and therefore the spectrum reflects the distribution of the p states within the valence band [2]. The Li K emission of lithium in LiF, Lepidolite  $K(Li,Al)_3(Si,Al)_4O_{10}(OH)_2$ , Bikitaite  $LiAlSi_2O_6(OH)_2$ , and Amblygonite  $(Li,Na)AlPO_4(F,OH)$  is presented below.

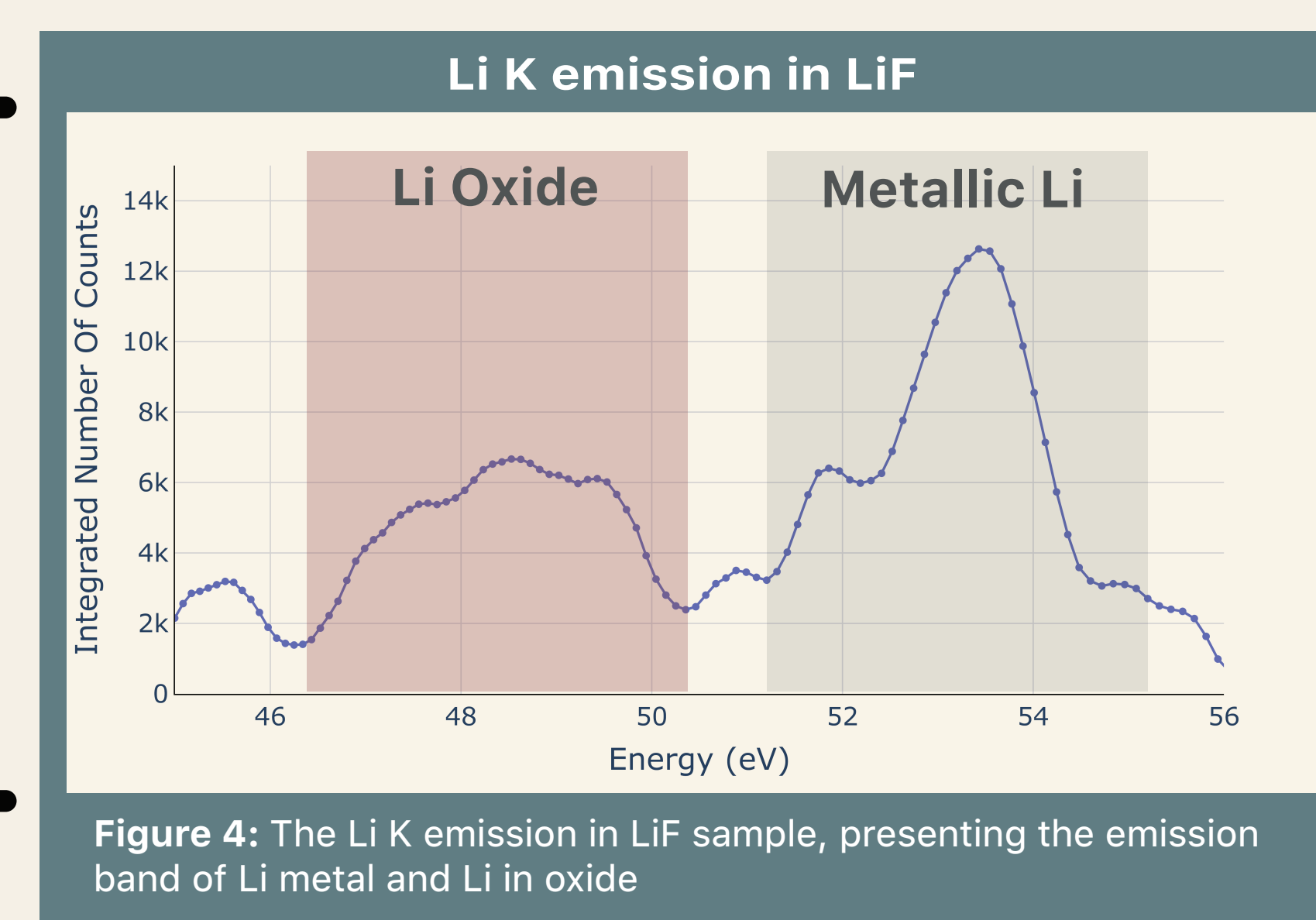

Accumulation of Li on the surface after F being driven off the sample resulting to an emission identical to Li metal.

Figure 4: The Li K emission in LiF sample, presenting the emission band of Li metal and Li in oxide

## X-ray emission spectra of minerals


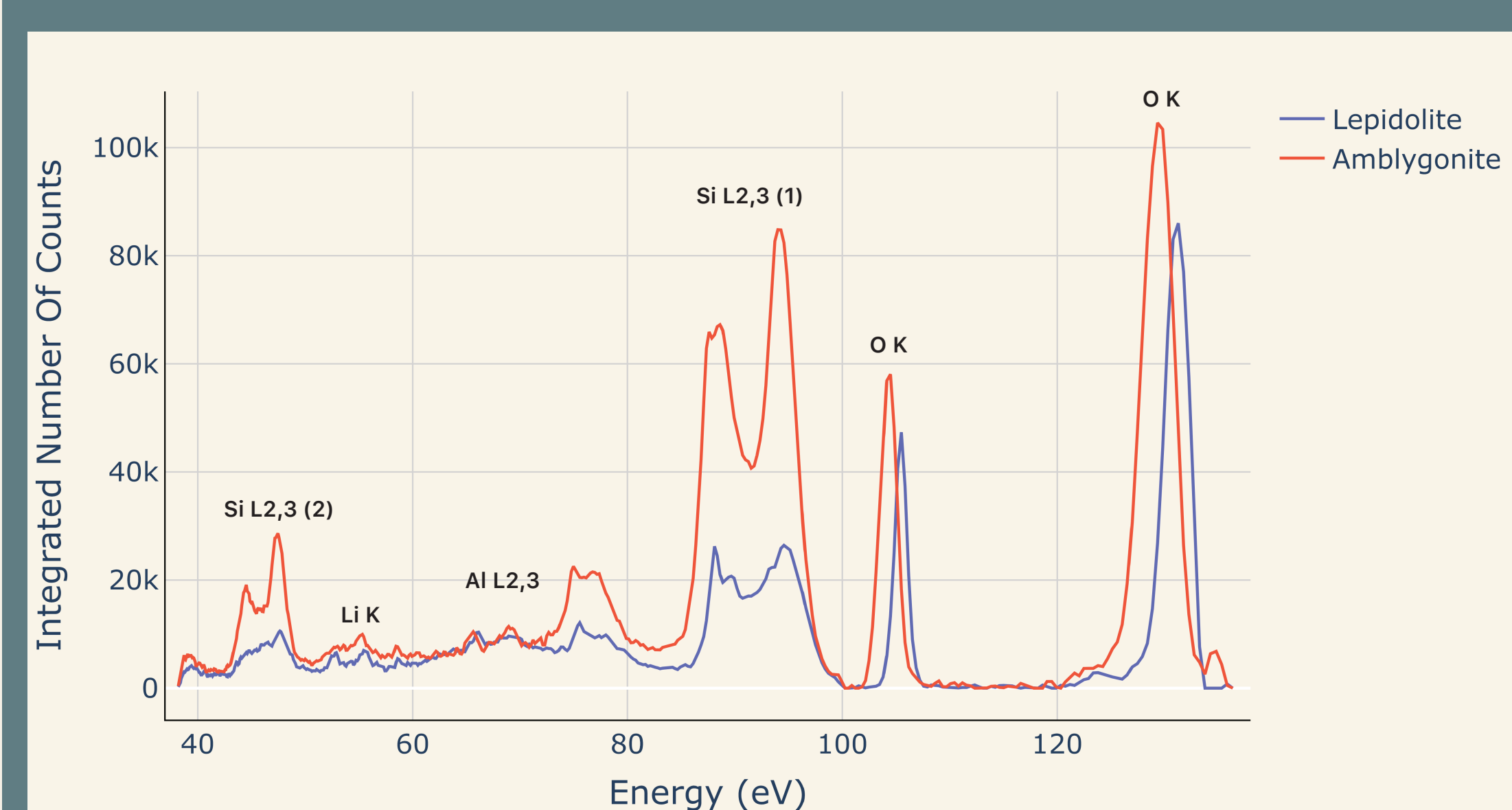


Figure 5: The x-ray emission spectra of Lepidolite and Amblygonite in the full range of the spectrometer

## Li K emission reported for the first time in minerals

Figure 6: The Li K emission in the minerals: Lepidolite, Amblygonite, and Bikitaite.

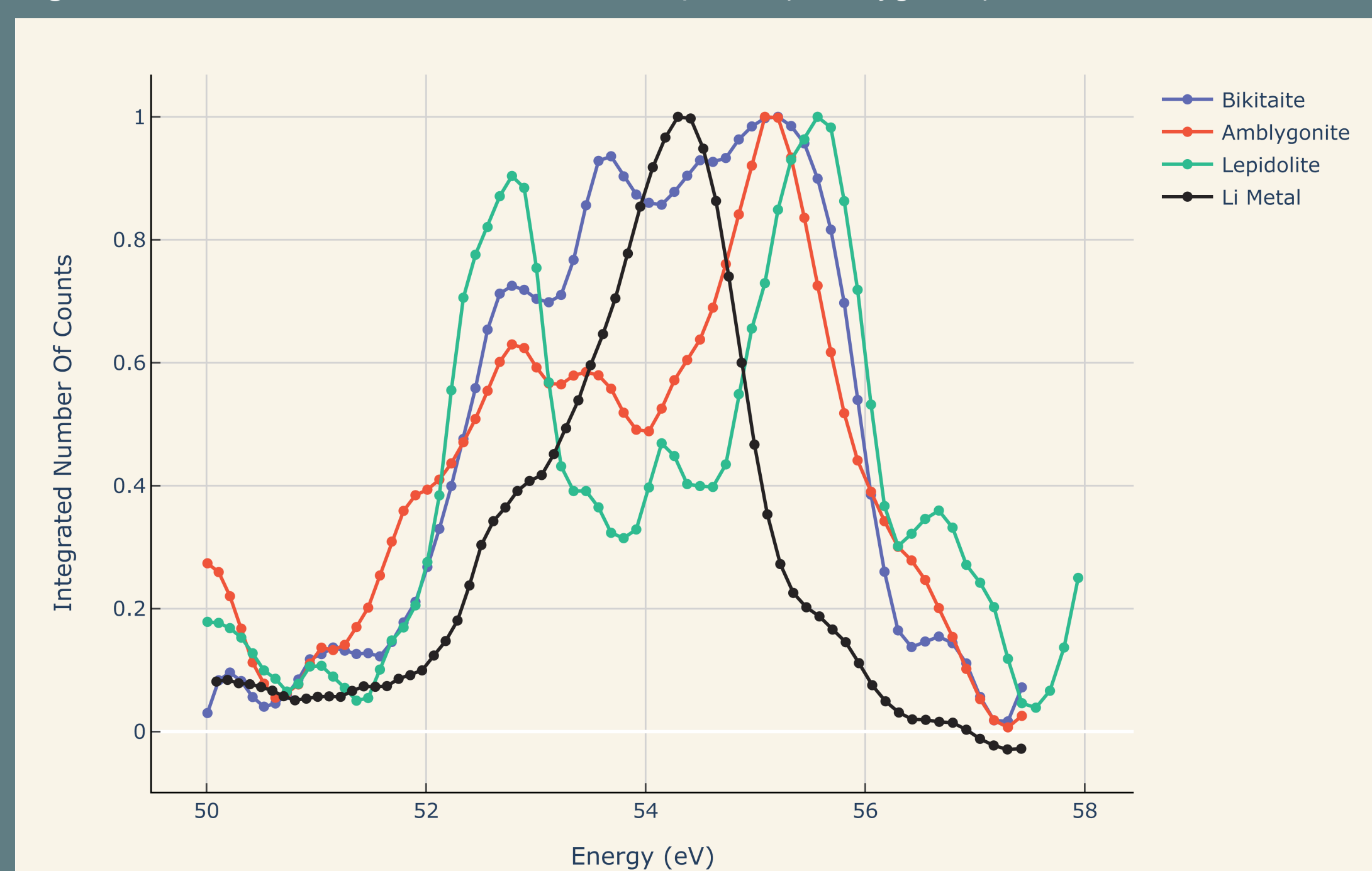

## Methods:

### The x-ray emission process:

1. **Excitation:** The sample is being bombarded with and accelerated beam of charged particles (electrons).
2. **Ionization:** The interaction of the beam of electrons with matter leads to the ionization of the core states.
3. **Emission:** the vacancy is filled with an electron from an outer shell, resulting in the emission of an x-ray photon that correspond to the energy difference between the two states.

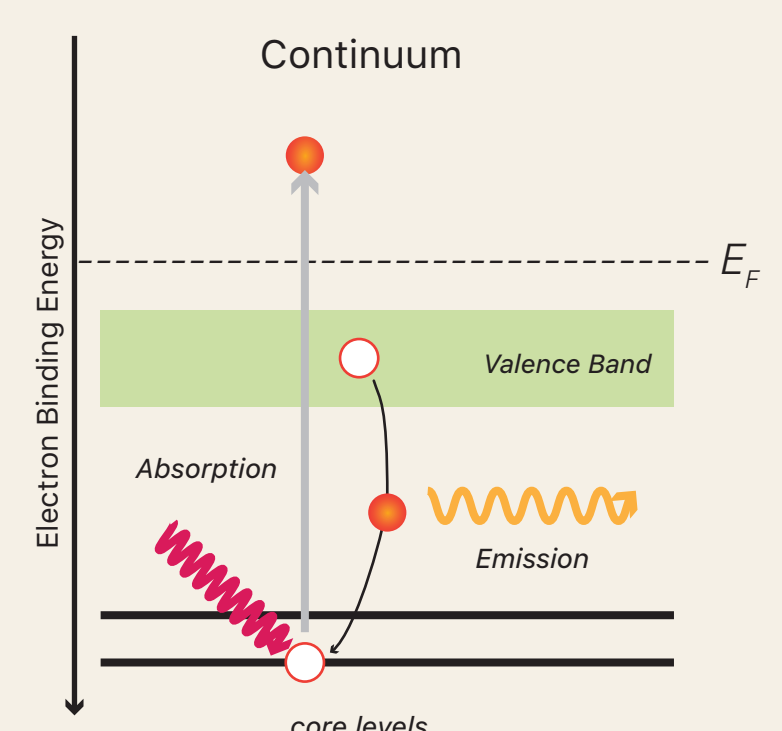


Figure 1: Illustration of the X-ray emission process

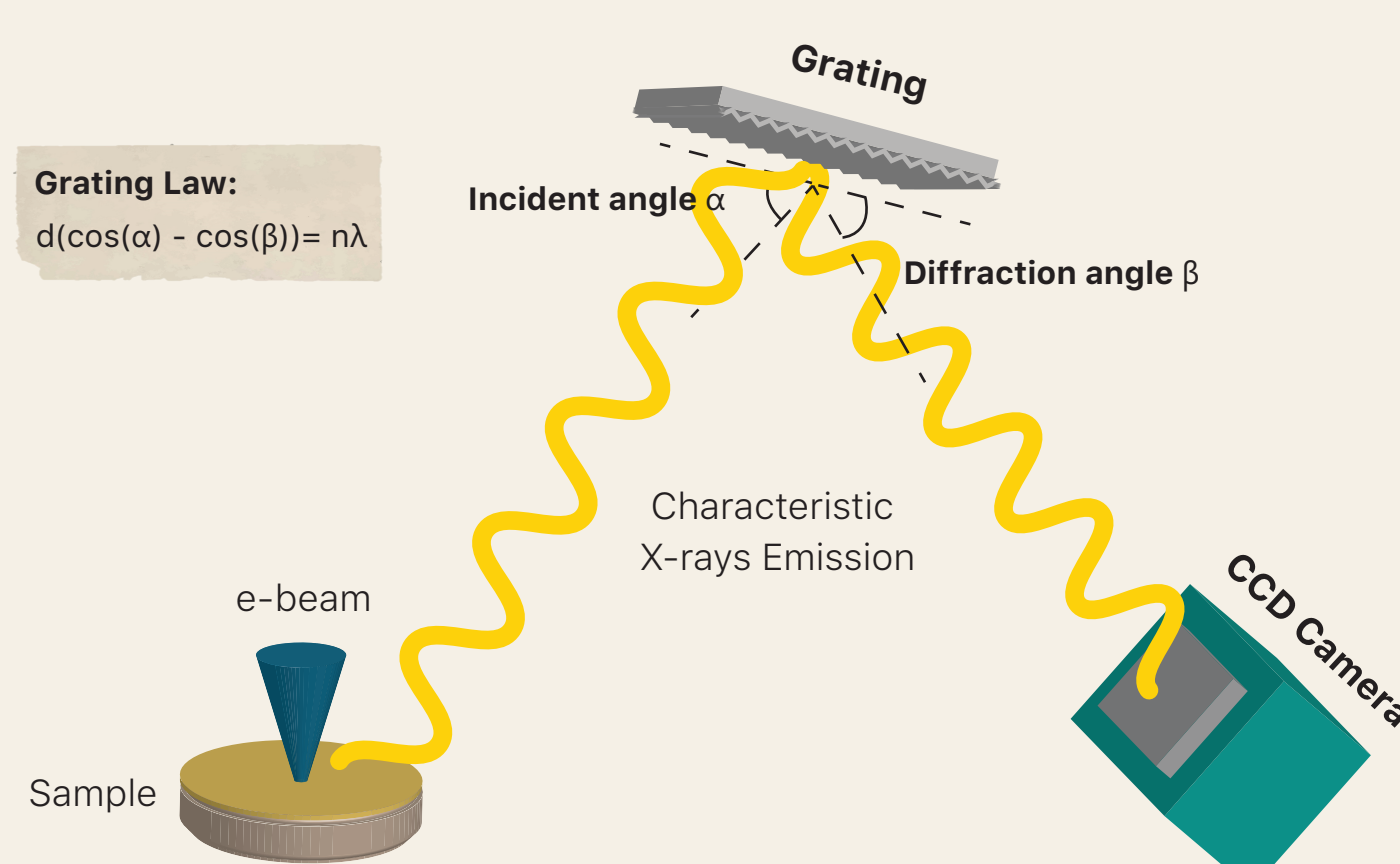


Figure 2: Wavelength dispersive spectroscopy (WDS) technique implemented inside the EPMA to acquire emission spectra

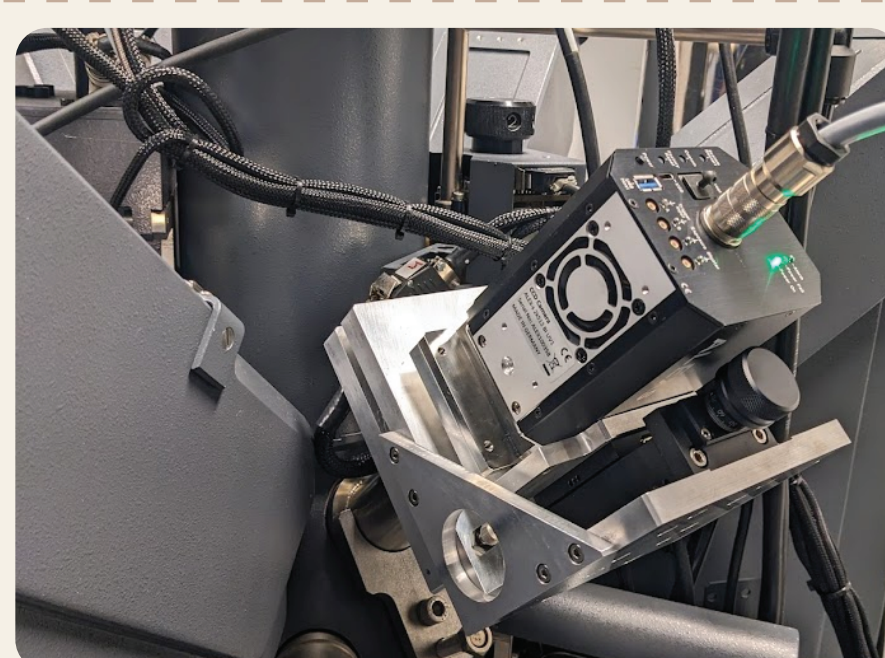


Figure 3: Image of reflection zone plate spectrometer mounted on an electron probe microanalyser (Cameca SX100).

### Why RZP instead of Multilayers?

- No scanning mechanism
- Higher resolution (~0.5 at Al L2,3 edge and ~0.25 eV at Mg L2,3 edge)
- Cover wider energy range

### Detection:

Reflection zone plates (RZP) are used in XES as a wavelength-dispersive element [1]. They consist of a special lamellar groove profile, each with a slightly different thickness. When x-rays interact with the zone plate, the different zones act as diffraction gratings, resulting in the x-rays diffracting and interfering with each other while focusing them at the same time.

### Sample analysis in three easy steps:

1. **Polishing**  
To remove the surface contaminations
2. **Coating**  
Sputter coating for non conducting samples
3. **Analyzing**  
Introducing the sample into the microprobe

## Conclusion:

We demonstrate the possibility to spectroscopically analyze Li K emissions in different compounds using a non destructive and accessible technique. **We also report for the first time the Li K emission bands in minerals.** Our results unlocked potential for various applications to study in details the fine structure of valence band of lithium.

## References:

- [1] A. Erko et al., Opt. Express, vol. 22, no. 14, p. 16897, Jul. 2014, doi: 10.1364/OE.22.016897.  
 [2] K. Hassebi et al., X-Ray Spectrometry, doi: 10.1002/xrs.3329



You want to keep in touch!  
 Scan the QR code or contact me at:  
 khalil.hassebi@sorbonne-universiter.fr