

4D electroactive and microporous polyHIPE-PEDOT scaffolds as a dynamic in vitro cell culture platform

Franziska Hahn, Ana Ferrández Montero, Mélodie Culot, Cedric Vancaeyzeele, Cedric Plesse, Rémy Agniel, Johanne Leroy-Dudal

▶ To cite this version:

Franziska Hahn, Ana Ferrández Montero, Mélodie Culot, Cedric Vancaeyzeele, Cedric Plesse, et al.. 4D electroactive and microporous polyHIPE-PEDOT scaffolds as a dynamic in vitro cell culture platform. ESB2023 (33rd Annual Conference of the European Society for Biomaterials), Sep 2023, Davos, Switzerland. hal-04400290

HAL Id: hal-04400290

https://hal.science/hal-04400290

Submitted on 17 Jan 2024

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.











4D electroactive and microporous polyHIPE-PEDOT scaffolds as a dynamic in vitro cell culture platform

Franziska Hahn ^{1,2}, Ana Ferrández Montero ^{1, 2}, Cédric Vancaeyzeele ¹, Cédric Plesse ¹, Rémy Agniel ², Johanne Leroy-Dudal ²

1) CY Cergy Paris Université, Laboratoire de Physicochimie des Polymères et des Interfaces (LPPI), Neuville sur Oise, France 2) CY Cergy Paris Université, Equipe de Recherche sur les Relations Matrice Extracellulaire-Cellule (ERRMECe), Mec'up group, Neuville sur Oise, France

Introduction

- In vivo, cells are surrounded by an extracellular matrix (ECM)
- Signals from the ECM such as biochemical or biophysical cues regulate cell behaviour
- 3D cell cultures offer biologically superior structures, but don't transmit dynamic mechanical signals that exist in the cell microenvironment in vivo
- To better recapitulate the cell microenvironment in vitro the aim of this work is to develop a 4D microporous and electroactive scaffold as an innovative cell culture platform:
- In situ electromechanical stimulation of cells
- Monitoring of the cell behavior

"4D cell culture platform" pore 10 -100 μm nterconnection up to 10 μm


Synthesis of polyHIPE-PEDOT scaffolds






Functionalization with electroactive polymer


PEDOT

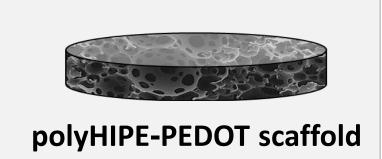



Electroactive polymeric scaffold:


Preserved polyHIPE morphology

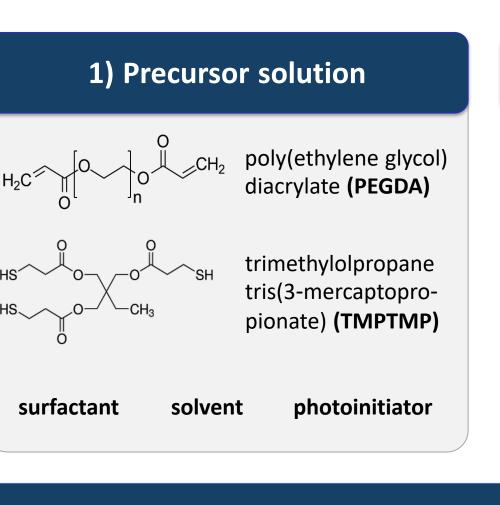


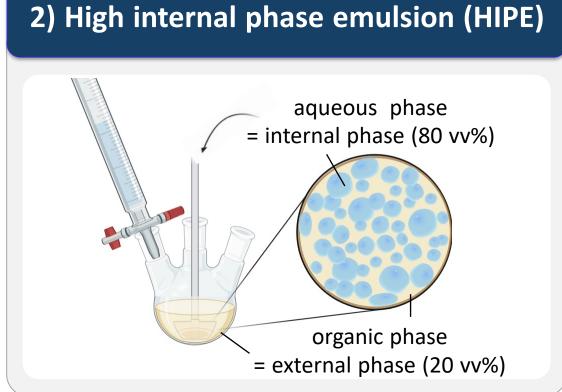

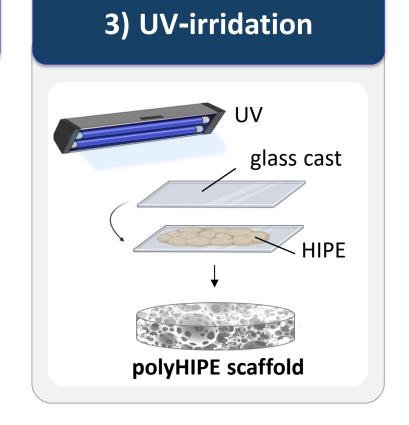


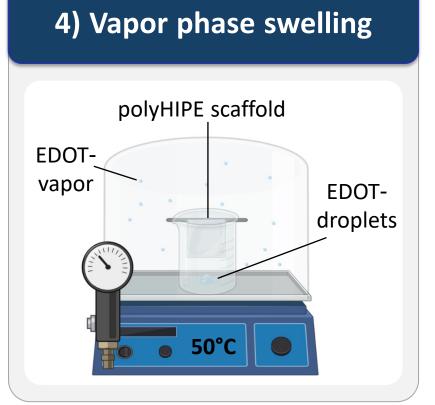


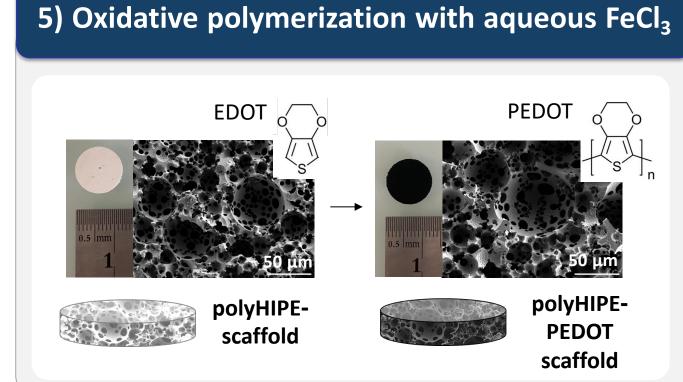


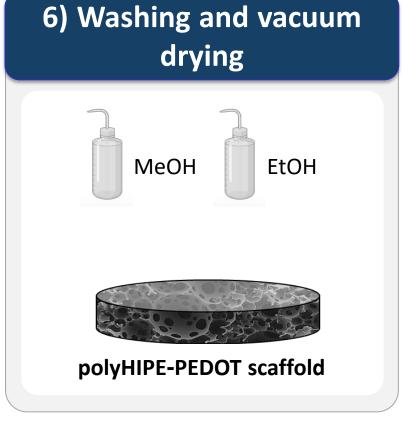


Actuation behavior of polyHIPE-PEDOT scaffolds

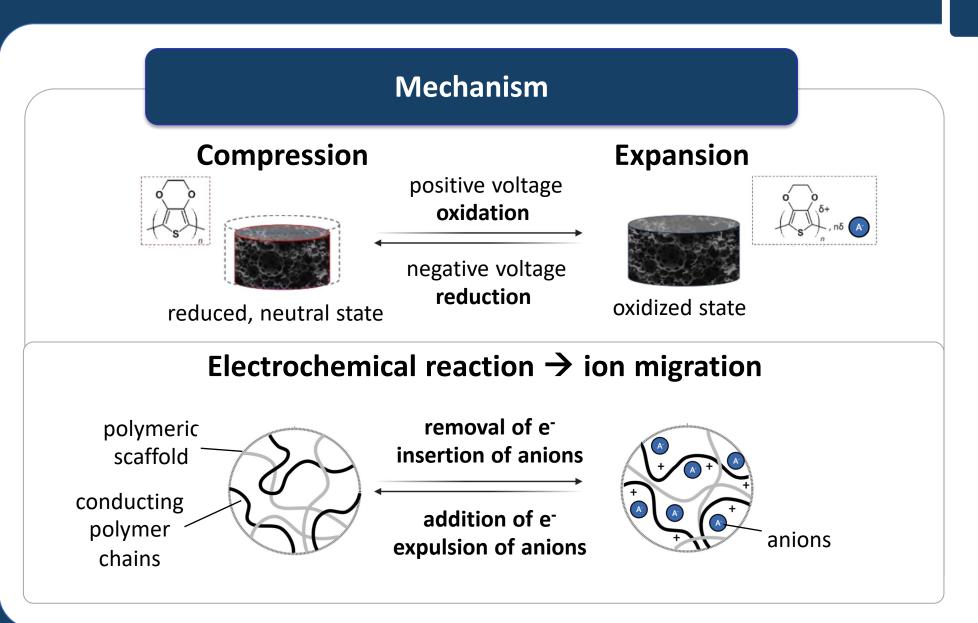

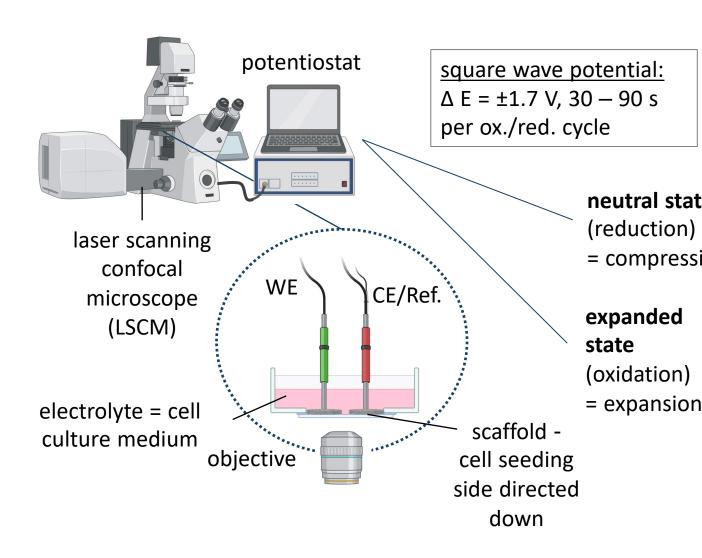

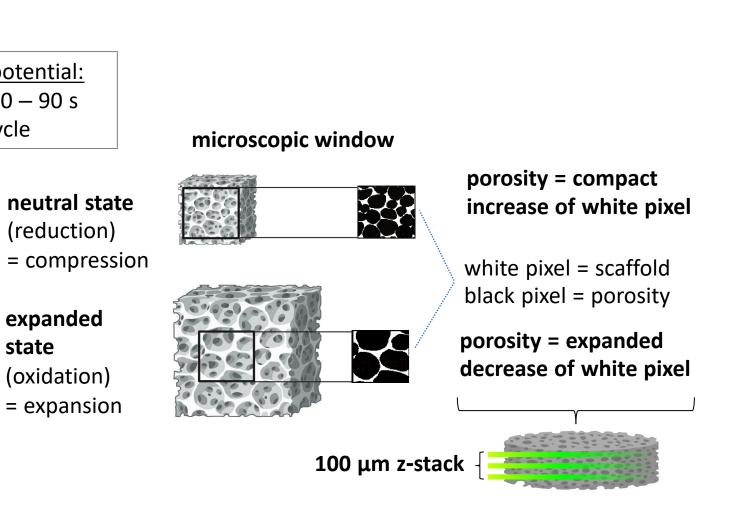


Actuation response of polyHIPE-PEDOT scaffolds to electrochemical stimulation

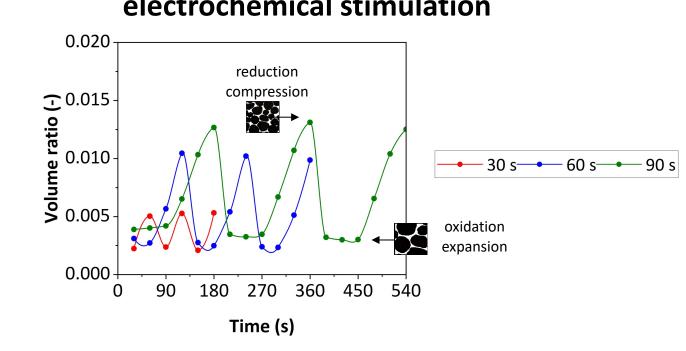

- PolyHIPE-PEDOT scaffolds showed direct response to electrochemical stimulation
- Reduction/oxidation time determines actuation frequency and volumetric expansion


Cytocompatibility

Cytotoxicity via LDH-activity

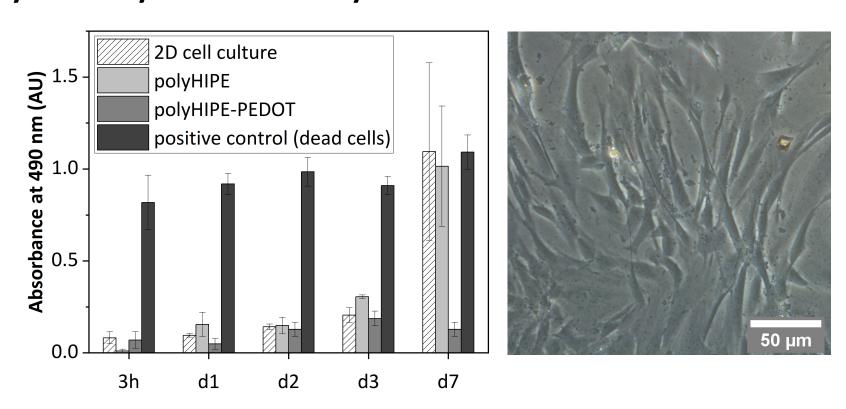


PolyHIPE-PEDOT scaffolds showed no cytotoxic effect over time for human dermal fibroblasts

Cell colonization

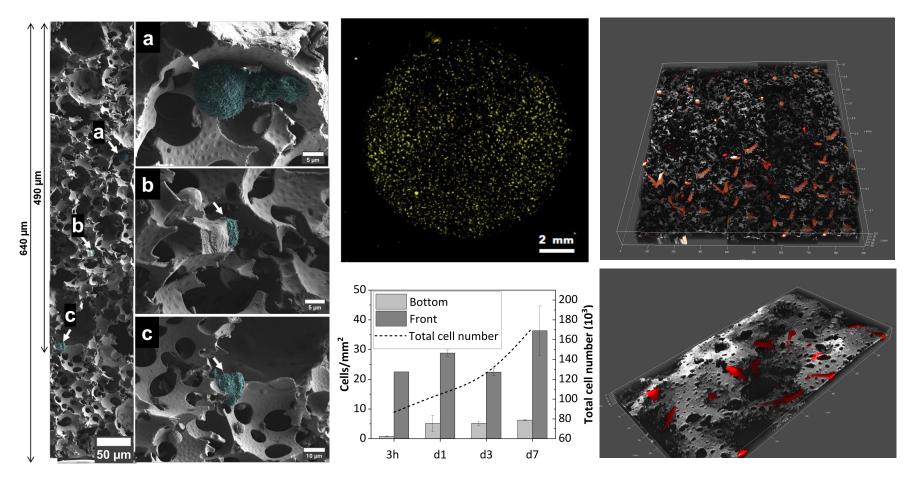

PolyHIPE-PEDOT scaffolds allow cells to adhere, spread and migrate into the microporous structure


Cell morphology and ECM

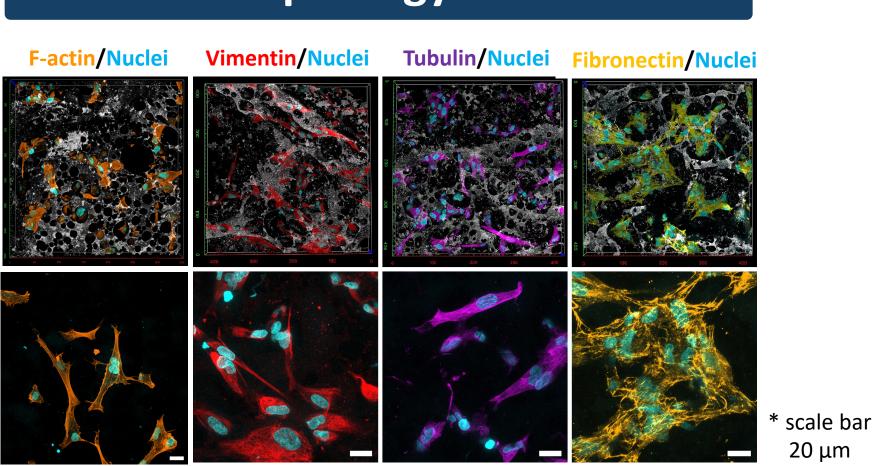


Fibroblast cells show typical spindle shaped cell morphology on polyHIPE-PEDOT scaffolds and secrete and organize fibronectin

Electromechanical stimulation in vitro


Electromechanical stimulation setup

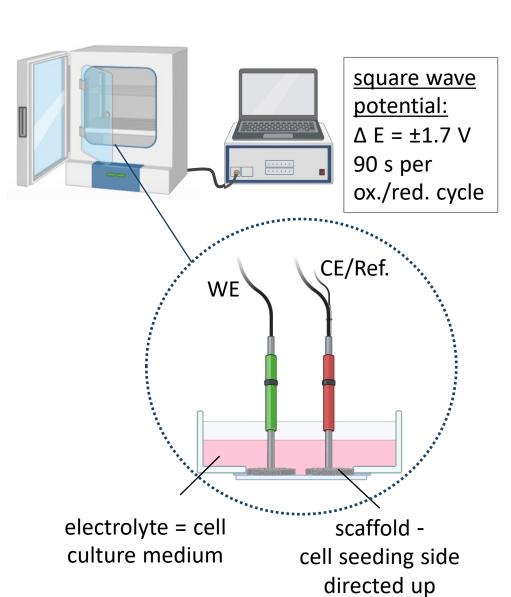


In situ and real-time monitoring

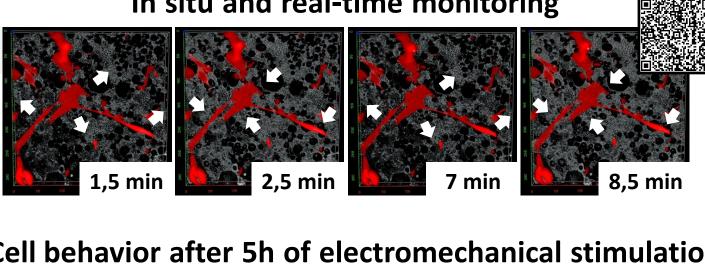

Cell behavior after 5h of electromechanical stimulation

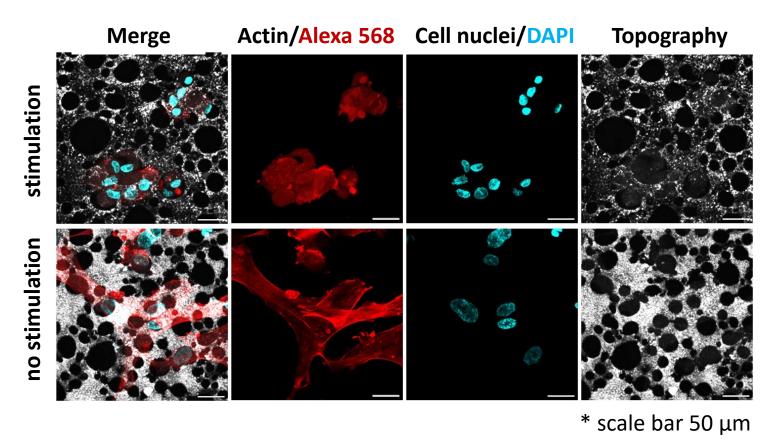

Cytotoxicity after 5h of electromechanical

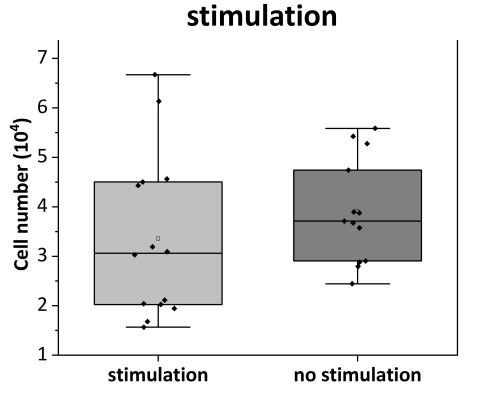

- 5h of electromechanical stimulation had no cytotoxic effect on cells
- Stimulated human dermal fibroblasts showed membrane deformation and changes in their morphology

Conclusion

- 4D microporous and electroactive scaffolds were successfully processed by the functionalization of a polyHIPE-scaffold with the conducting polymer PEDOT
- PolyHIPE-PEDOT scaffolds comply with the requirements of cell culture and allow cell adhesion, spreading and migration inside the microporous structure
- Typical 3D cell morphology and ECM production on polyHIPE-PEDOT scaffolds were found
- Electromechanical stimulation setup was designed and implemented for a CLSM and allowed the *in situ* and real-time monitoring of cells during stimulation
- 5h of electromechanical stimulation had no cytotoxic effect on human dermal fibroblasts
- Electromechanical stimulation leads to changed cell phenotype of human dermal fibroblasts, which remains to be explored

Acknowledgments