

## The High-Precision Liquid Chromatography with Electrochemical Detection (HPLC-ECD) for Monoamines Neurotransmitters and Their Metabolites: A Review

Bruno P Guiard, Guillaume Gotti

### ▶ To cite this version:

Bruno P Guiard, Guillaume Gotti. The High-Precision Liquid Chromatography with Electrochemical Detection (HPLC-ECD) for Monoamines Neurotransmitters and Their Metabolites: A Review. Molecules, 2024, 29 (2), pp.496. 10.3390/molecules29020496 . hal-04410209

## HAL Id: hal-04410209 https://hal.science/hal-04410209

Submitted on 22 Jan 2024

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.






# **The High-Precision Liquid Chromatography with Electrochemical Detection (HPLC-ECD) for Monoamines Neurotransmitters and Their Metabolites: A Review**

Bruno P. Guiard <sup>1,2</sup> and Guillaume Gotti <sup>1,2,\*</sup>

- <sup>1</sup> Centre de Recherches sur la Cognition Animale (CRCA), CNRS UMR5169, 31062 Toulouse, France; bruno.guiard@univ-tlse3.fr
- <sup>2</sup> Centre de Biologie Intégrative (CBI), Faculté Sciences Ingénierie (FSI), Université de Toulouse III, 31062 Toulouse, France
- \* Correspondence: guillaume.gotti@univ-tlse3.fr

**Abstract**: This review highlights the advantages of high-precision liquid chromatography with an electrochemical detector (HPLC-ECD) in detecting and quantifying biological samples obtained through intracerebral microdialysis, specifically the serotonergic and dopaminergic systems: Serotonin (5-HT), 5-hydroxyindolacetic acid (5-HIAA), 3,4-dihydroxyphenylacetic acid (DOPAC), dopamine (DA), 3-metoxytryptamin (3-MT) and homovanillic acid (HVA). Recognized for its speed and selectivity, HPLC enables direct analysis of intracerebral microdialysis samples without complex derivatization. Various chromatographic methods, including reverse phase (RP), are explored for neurotransmitters (NTs) and metabolites separation. Electrochemical detector (ECD), particularly with glassy carbon (GC) electrodes, is emphasized for its simplicity and sensitivity, aimed at enhancing reproducibility through optimization strategies such as modified electrode materials. This paper underscores the determination of limits of detection (LOD) and quantification (LOQ) and the linear range (L.R.) showcasing the potential for real-time monitoring of compounds concentrations. A non-exhaustive compilation of literature values for LOD, LOQ, and L.R. from recent publications is included.

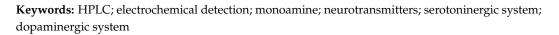

#### 1. Introduction

Monoamine neurotransmitters (NTs) are a class of NTs that include brain chemicals involved in the transmission of signals between nerve cells, or neurons. The two main types of monoamine NTs are catecholamines and indolamines: catecholamines include dopamine (DA), norepinephrine (NE), and epinephrine (E) and indolamines include serotonin (5-hydroxytryptamine or 5-HT).

These monoamine NTs play a crucial role in the regulation of various physiological and behavioral processes, such as the regulation of mood, sleep, appetite, reward, attention, and stress. Imbalances in the functioning of monoamine NTs are associated with various neurological and psychiatric disorders.

#### 1.1. The Serotonergic System

The serotonergic system utilizes 5-HT as a neurotransmitter and plays a pivotal role in regulating diverse physiological and psychological processes, including mood, sleep, appetite, thermoregulation, and heart rate [1].

5-HT is synthesized from the essential amino acid tryptophan (TRP). Within neurons, TRP undergoes conversion to 5-hydroxytryptophan (5-HTP), catalyzed by the tryptophan hydroxylase enzyme. Subsequently, an enzymatic reaction involving aromatic amino acid decarboxylase (AADC) transforms 5-HTP into 5-HT (Figure 1). Following synthesis, the



Citation: Guiard, B.P.; Gotti, G. The High-Precision Liquid Chromatography with Electrochemical Detection (HPLC-ECD) for Monoamines Neurotransmitters and Their Metabolites: A Review. *Molecules* 2024, 29, 496. https://doi.org/10.3390/ molecules29020496

Academic Editor: Ritin Sharma

Received: 15 December 2023 Revised: 15 January 2024 Accepted: 17 January 2024 Published: 19 January 2024



**Copyright:** © 2024 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (https:// creativecommons.org/licenses/by/ 4.0/). 5-HT is stored in vesicles via the vesicular monoamine transporter (VMAT) within the nerve endings of serotonergic neurons. Upon the arrival of an action potential at the nerve ending and the membrane depolarization, it triggers the release of 5-HT from vesicles into the synaptic cleft. Here, it can interact with receptors located on the postsynaptic membrane, known as serotonin receptors or 5-HT receptors. Extracellular 5-HT may also bind and activate inhibitory receptors expressed on 5-HT neurons (5-HT1A/1B/1D) thereby exerting negative feedback on its own release. Intracellular excess of 5-HT can be sequentially broken down by monoamine oxidase (MAO) and aldehyde dehydrogenase, resulting in the accumulation of 5-hydroxyindoleacetic acid (5-HIAA) in the intracellular space. All these receptors are present in the central nervous system and various tissues throughout the body. The effects of 5-HT vary based on the specific receptors it binds to.

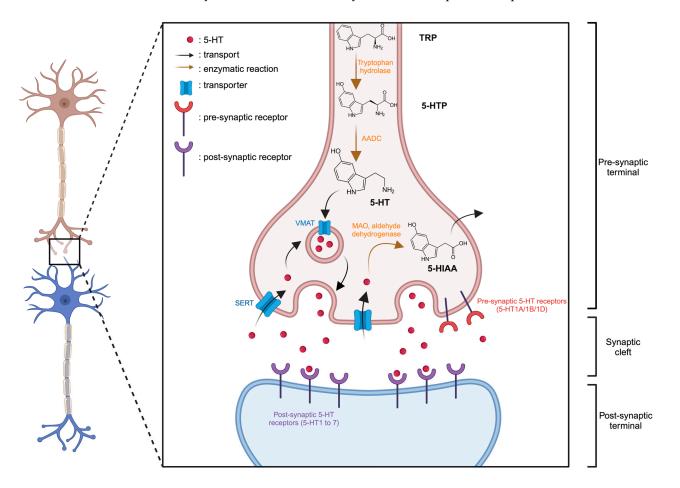


Figure 1. Representation of serotonergic synapse. Created with BioRender.com accessed on 15 January 2024.


Dysfunctions in the serotonergic system may contribute to several disorders, including neurodegenerative disease like Parkinson's [2] or Alzheimer's [3,4], and also psychiatric disorders such as anxiety [5] and depression [6].

#### 1.2. The Dopaminergic System

The dopaminergic system is a neurotransmission system in the brain that uses DA as a neurotransmitter. DA is a chemical compound that plays a crucial role in various brain functions, including the regulation of motivation [7], pleasure, reward, and movement [8].


The dopaminergic system is composed of several neural pathways, with the most well-known being the mesolimbic, mesocortical, and nigrostriatal pathways. It plays a key role in various physiological and behavioral processes, and its dysfunction is associated with several neurological and psychiatric disorders. For example, overactivity of the dopaminergic system is linked to schizophrenia [9,10], while deficiency is associated with Parkinson's disease [11]. Medications that modulate dopaminergic activity are often used to treat these disorders.

DA is synthesized in the neuron from tyrosine (TYR) through a two-step enzymatic reaction involving first, the tyrosine hydroxylase, synthesizing the levodopa (L-DOPA) and then, the DOPA decarboxylase. The DA can then be introduced into vesicles by the VMAT and released into the synaptic cleft. Degradation of DA can be catalyzed by MAO into DOPAC within the neuron. Two pathways are possible for its transformation into homovanillic acid (HVA), either through catecholamine-O-methyl-transferase (COMT) and MAO sequentially, with 3-metoxytryptamine (3-MT) as an intermediate, or by directly degrading DOPAC via COMT (Figure 2).

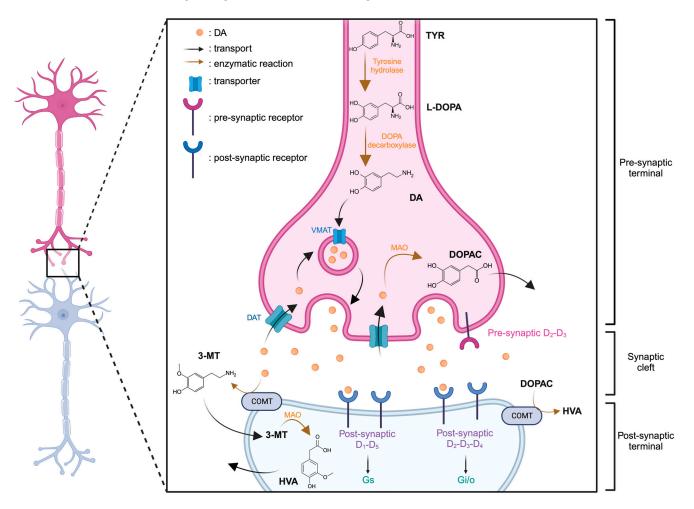

**Figure 2.** Representation of dopaminergic synapse. Created with BioRender.com accessed on 15 January 2024.

Two additional catecholamines, derived from dopamine and integral to the neurotransmitter family, are norepinephrine (NE) and epinephrine (E). E plays a dual role, functioning as a neurotransmitter within the central nervous system and as a hormone within the adrenal medulla. NE, conversely, acts as a metabolic precursor to epinephrine and operates as a neurotransmitter in the central nervous system. Both NE and E are crucial in regulating fundamental processes such as the sleep-wake cycle, vigilance, stress response [12], emotions [13], learning [14], and memory [15]. Their multifaceted roles underscore their significance in orchestrating key physiological and cognitive functions within the intricate framework of the human body. Microdialysis, employed for the past forty years [16,17], is a real-time sampling technique facilitating the analysis of extracellular NTs [18,19]. It functions on the principle of passive diffusion of small molecular-weight compounds through a porous membrane, transitioning from the region with the highest NTs concentration, namely the synaptic extracellular space, to the less concentrated area, the dialysis probe perfused with a neurotransmitter-free buffer solution at physiological pH. This technique allows the collection of dialysates at variable intervals, with flow rates ranging from 0.5 to  $1.5 \,\mu L \,min^{-1}$ , dependent on the experimental setup and brain area studied (Figure 3). These samples, which include monoamines and their metabolites, play a crucial role in comprehending the functioning of the serotonin and dopamine systems and associated disorders.

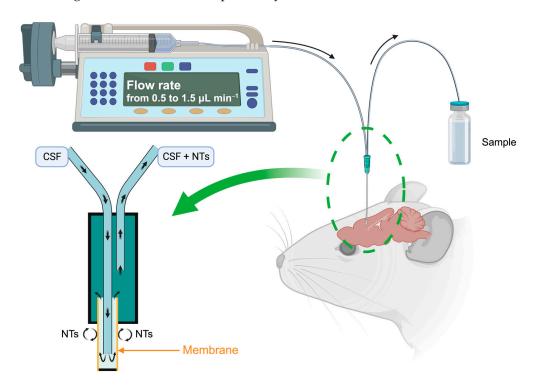

**Figure 3.** Schematic principle of the microdialysis on a brain mouse. Created with BioRender.com accessed on 15 January 2024.

The reported sensitivity limit for 5-HT in various scientific publications is approximately 0.5 fmol per sample with a signal-to-noise ratio = 2 [20]. To comprehensively grasp the sero-tonergic and dopaminergic systems functionalities and its associated disorders and diseases, a combination of microdialysis sampling methods and high-performance analytical techniques is crucial for quantifying variations in extracellular NTs and metabolite concentrations.

#### 3. Analytical Methods for Detection and Quantification of NTs and Metabolites


The analysis methods have been developed for the detection and quantification of NTs in vivo to understand specific neurotransmission dysfunctions in neurodegenerative diseases or in psychiatric disorders.

#### 3.1. Gas Chromatography and Capillary Electrophoresis

Gas chromatography (GC) [21,22] and capillary electrophoresis (CE) [23,24] are analytical techniques used in these research areas. However, they have several drawbacks: GC requires the use of compound derivatization to make them more volatile [25], and CE is a limited technique for the analysis of complex mixtures. HPLC, on the other hand, allows for the analysis of biological samples from microdialysis without the need for chemical derivatization or extraction steps.

#### 3.2. HPLC

Since the late 1980s, research teams led by Donzanti et al. [26], Peinado et al. [27], and Rogers et al. [28] have focused their studies on the detection of amino acid neurotransmitters, particularly  $\gamma$ -aminobutyric acid (GABA), achieving a sensitivity of about 0.1 nM. HPLC analysis of NTs is characterized by its speed and selectivity. Typically, a reverse-phase (RP) system is used to separate NTs due to their hydrophobic nature. For instance, Bidel et al. developed an HPLC method for the separation of monoamines and their metabolites (Figure 4a). They achieved a total elution time of about 12 min, a minimum resolution between each peak of 2.2, and a number of theoretical plates (N) ranging from 2900 to 3800 [29]. These values are largely sufficient to validate the use of HPLC.

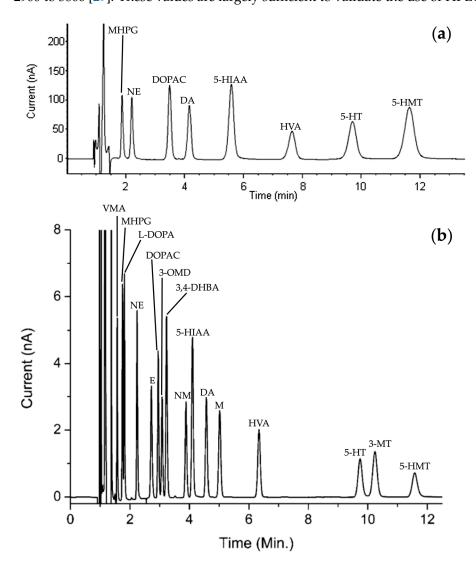

**Figure 4.** HPLC chromatogram reprinted and adapted with permission from [29] Copyright 2016 John Wiley and Sons (**a**) and UHPLC chromatogram reprinted and adapted with permission from [30] Copyright 2013 American Chemical Society (**b**) of a standards solution containing of each NTs and metabolites: DOPAC, DA, 5-HIAA, HVA, 3-MT and 5-HT. For information, other NTs and metabolites have been presented: 3-methoxy-4-hydroxyphenylglycol (MHPG), norepinephrine (NE), 5-hydroxymethyl tolterodine (5-HMT), D,L-4-hydroxy-3-methoxy-mandelic acid (VMA), L-3,4-dihydroxyphenylalanine (L-DOPA), epinephrine (E), 3,4-dihydroxybenzylamine hydrobromide (3,4-DHBA), 3-methoxy-L-tyrosine (3-OMD) and D,L-metanephirine (M).

By way of comparison, Reinhoud et al. developed a UHPLC method similar conditions (see Section 3.3 for more details) under reverse phase conditions, enabling the separation of

NTs in a 12 min run time elution (Figure 4b). The columns of UHPLC have smaller particle sizes. Consequently, chromatographic peaks will be finer and better defined, allowing for improved separation of compounds. Moreover, the number of theoretical plates in this study has been estimated to be close to 200,000 [30].

More recently, hydrophilic interaction liquid chromatography (HILIC) columns have been employed [31,32]. These columns offer reduced system pressure, allowing for enhanced flow rates and shorter analysis times. For example, in the study by Zhou et al., the use of a low-granularity column (1.7  $\mu$ m) combined with a flow rate of 0.4 mL min<sup>-1</sup> is employed. These experimental conditions result in very short analysis times, on the order of 4 min. It is noteworthy that the Multiple Reaction Monitoring (MRM) mode of the mass spectrometer enhances the selectivity of the measurements [33].

The utilization of HILIC column necessitates coupling with mass spectrometry because of its incompatibility with the salts used with the electrochemical detection.

#### 3.3. Detection of NTs and Their Metabolites

Various types of detectors can be employed, including fluorescence detectors (FLD) [34–36], mass spectrometry (MS) [37–39], sensors and biosensors [40–42] or electrochemical detector (ECD). ECD holds an advantage over MS detectors due to its lower operating cost and easy maintenance. For instance, weekly polishing of a glassy carbon (GC) working electrode with alumina powder is adequate to maintaining a reflective surface, thereby enhancing result repeatability and reproducibility.

ECD has also been extensively utilized for the detection and quantification of NTs and their metabolites [43–46]. It proves to be a straightforward and sensitive method for quantifying small amounts of NTs. One of its notable benefits is its simplicity and minimal maintenance.

Several types of electrochemical detectors exist, including coulometric and amperometric detectors. The principle of the coulometric detector lies in the measurement of the current charge produced during electrolysis at a constant potential. Compounds are then oxidized or reduced almost entirely. Coulometric detectors are less sensitive to variations in flow rate and temperature than amperometric detectors. However, electrodes are easily contaminated due to their large electrolysis surface area. The measurement of the charge during the electrochemical reaction is based on Faraday's law (Equation (1)).

$$Q = n F N_O \tag{1}$$

with Q the charge (C or A s), n the number of electrons exchanges during the electrochemical detection, F the Faraday's constant (96,485 C mol<sup>-1</sup>) and N<sub>O</sub> the number of moles of the O compound.

Amperometric detection generally provides high sensitivity and the response remains linear over a wide range of concentrations. For these reasons, they are more widespread than coulometric detectors. Compared to coulometric detection, electrolysis consumes only about 10% of the analyte. The value of the measured current in an amperometric detector depends on the geometry of the electrode used (Table 1).

Regardless of the detector used, the number of electrons exchanged determines the measured current during the analysis. In case of detecting 5-HT, the number of electrons exchanges is 2 during the oxidation, as shown in Figure 5.

Electrode Geometry	Limiting Current Equation			
Tubular	$i = 1.61 nFC_O \left(\frac{DA}{r}\right)^{2/3} \upsilon^{1/3}$			
Planar, parallel flow in channel	$i = 1.47 nFC_O \left(\frac{DA}{b}\right)^{2/3} v^{1/3}$			
Planar, perpendicular flow	$i = 0.903 n F C_O D^{2/3} \upsilon^{-1/6} A^{3/4} U^{1/2}$			
Wall jet	$i=0.898nFC_{O}D^{2/3}\upsilon^{-5/12}a^{-1/2}A^{3/8}\upsilon^{3/4}$			

**Table 1.** Limiting current equation as a function of geometry electrode used during the amperometric detection [47].

With: i the intensity (A),  $C_0$  concentration of the O compound (mol cm<sup>-3</sup>), D the diffusion coefficient (cm<sup>2</sup> s<sup>-1</sup>), A the electrode area (cm<sup>2</sup>), r radius of the tubular electrode (cm), v average volume flow rate (cm<sup>3</sup> s<sup>-1</sup>), b channel height (cm), U flow velocity (cm s<sup>-1</sup>) and a diameter of jet intel (cm).

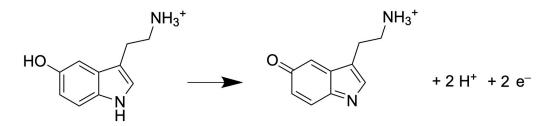



Figure 5. Half reaction of 5-HT oxidation [48].

For electroactive compounds, charge can be correlated with the moles of the considered compound, extending to concentration. It is crucial to emphasize that the number of electrons exchanged in the electrochemical reaction is contingent on the applied potential. If the overvoltage is insufficient, then the electrochemical reaction remains incomplete.


ECD has some limitations: this detector is not universal and compounds must be electroactive to be detectable. In the same domain, GABA, a neurotransmitter involved in neurodegenerative diseases [49,50], can be analyzed by HPLC-ECD. However, GABA is not electroactive and requires chemical derivatization. In the Panrod et al. study, two derivatization pathway involving aldehydes moieties have been employed to be able to detect GABA [51]. Therefore, the development of a single analytical method for all NTs with an electrochemical detector appears to be challenging in light of the mentioned limitations.

The two most frequently employed working electrode materials are GC and borondoped diamond (BDD). These materials showcase diverse physicochemical characteristics: the BDD electrode boasts the largest electrochemical window, facilitating the detection of numerous electroactive compounds. However, as indicated in the study by Zhang et al. in 2016, a high overpotential seems to reduce the peak area of the analytes (Figure 6). The sensitivity of electrode passivation is higher for a BDD electrode compared to a GC electrode [52]. Consequently, it seems that reproducibility is superior with a GC electrode.

The GC electrode exhibits a narrower electrochemical window compared to the BDD electrode. Consequently, the analysis range is restricted, and the oxidation of water impedes the detection of compounds with high oxidation potentials [54]. Analyzing compounds with potentials involving reduction is, to our knowledge, impossible due to the reduction of dissolved oxygen [55]. In Figure 7a, the detection current is depicted as a function of the potential applied for monoamines NTs with a GC working electrode. Ultimately, the choice between the two electrodes depends on the application when the compounds have relatively low oxidation potentials. Then, it is preferable to use a GC working electrode.

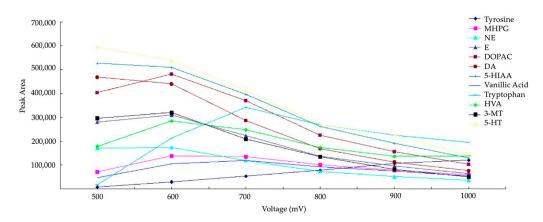

**Figure 6.** Hydrodynamic voltammograms of monoamines NTs on BDD working electrode reprinted and modified with permission from [53] Copyright 2016 Elsevier.

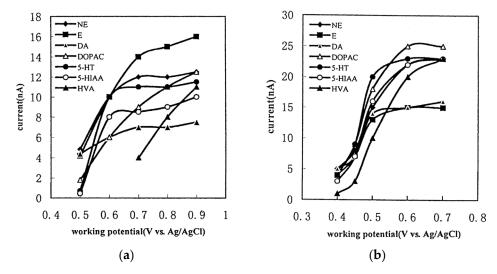

**Figure 7.** Hydrodynamic voltammograms of monoamines NTs and their main metabolites on GC working electrode (**a**) and P-pABA modified working electrode (**b**) reprinted with permission from [56] Copyright 2001 Elsevier.

To enhance the sensitivity of NTs detection, two optimizations are considered: firstly, increasing the surface-to-volume ratio of the electrochemical cell, and secondly, using a modified electrode material to enhance its physicochemical characteristics. In the study by Xu et al., the working electrode on GC is modified through an electrochemical reaction with poly(para-aminobenzoic acid) (P-pABA) (Figure 7b) [56]. With a 0.7 V vs. Ag/AgCl potential detection, the sensitivity of measurement improved from 8.5 nA to 23 nA, 11 nA to 23 nA, 7 nA to 16 nA, 9 nA to 25 nA and 4 nA to 23 nA for 5-HIAA, 5-HT, DA, DOPAC and HVA, respectively. Augmenting sensitivity enhances the quantification of NTs.

In conclusion, the ECD is widely employed for the detection of monoamines due to its user-friendly nature and low maintenance cost. Its sensitivity enables the quantification of very low concentrations of NTs. To expand the capabilities of this detector, it would be possible to develop chemical derivatizations either directly on the compounds or on the electrode surface. These chemical modifications could induce a broader spectrum of applications.

#### 3.4. Limit of Detection and Quantification

By convention, the limit of detection (LOD) and limit of quantification (LOQ) are established at signal-to-noise ratios of 3 and 10, respectively. The linear range (L.R.) denotes the measurement zone where the device-recorded signal is directly proportional to the quantity or concentration of the substance being measured.

Table 2 compiles the various values of LOD, LOQ, and the linear range for the NTs and their metabolites. Additionally, the different experimental conditions have been reported for informational purposes. The studies are ranked based on the particle size of the chromatographic  $C_{18}$  column used, from smallest to largest. The list of presented publications is not exhaustive. Other articles have been published on the use of HPLC-ECD but they did not detail the analytical development and LOD/LOQ of NTs. Those articles were not included in this summary table. For instance, we can cite several recent works by Guo et al. (2023) [57], Lei et al. (2023) [58], and Wang et al. (2022) [59].

According to these studies, the sensitivity of the targeted molecules can approach values close to  $10^{-9}$  mol L<sup>-1</sup> for HPLC system and less, close to  $10^{-11}$  mol L<sup>-1</sup>, in UHPLC cases. As mentioned earlier, the use of columns with small particle sizes refines elution peaks and consequently lowers quantification limits. However, using such columns induces high back pressure. These low values facilitate the collection of microdialysis samples in the range of microliter, enabling real-time monitoring of neurotransmitter concentrations.

In electrochemical detection processes, the use of salts becomes crucial to increase conductivity and consequently achieve optimal performance. Typically, the aqueous phase of the mobile phase incorporates essential components, such as chelating compounds (e.g., sodium citrate, EDTA), which play a role in preventing coagulation [60], an anionic surfactant (sodium dodecyl sulfate, octane sulfonic acid) to prevent interference with lipids and enhance the retention of anionic compounds [61], and finally, a buffer (phosphate, citrate).

Most of the time, the pH of the mobile phase is buffered in an acidic environment. The more acidic the eluent, the faster the compounds will elute. Additionally, some compounds are sensitive and degrade more rapidly at basic pH. For instance, this is the case with DA, which oxidizes at a significant rate under basic pH conditions: the half-life at pH 5.6, 7.1, and 7.4 is 13.5 days, 19.7 min, and 4.95 min, respectively [62].

Moreover, the mobile phase, predominantly consisting of the aqueous phase with salts, is often augmented with methanol and/or acetonitrile. This strategic addition enhances the eluent's apolar characteristics, ensuring a comprehensive and effective electrochemical detection setup.

Reference		HPLC-ECD Conditions				Application	
Reinhoud et al. (2013) [30]	Methods and Ap- plication	Mobile phase	Flow r. Column Injecti : 100 mM phosphoric acid, 1 ACN,	100 mm, 1.7 $\mu$ m particles size) ate = 50 $\mu$ L min <sup>-1</sup> temperature = 37 °C on volume = 5 $\mu$ L 00 mM citric acid, 8 mM KCl, 0 mM OSA. /H <sub>2</sub> O (8/92; $v/v$ ) pH = 3 gelectrode (+0.65 V vs. Ag/Ag6	Development of analytical method for the determination of monoamines concentrations and application to a microdialysis in a prefrontal cortex.		
	Detection (mol L <sup>-1</sup> )		5-HT	5-HIAA	DA	DOPAC	HVA
		LOD	$8.3 imes10^{-11}$	$3.5 imes10^{-11}$	$4.2  imes 10^{-11}$	$5.0 imes10^{-11}$	$4.7 imes10^{-11}$
Ferry et al. (2014) [63]	Methods and Ap- plication	$\begin{array}{c} C_{18} \mbox{ column } (0.32 \times 100 \mbox{ mm}, 1.9 \mbox{ \mum} \mbox{ particles size}) \\ Flow \mbox{ rate } = 8.5 \mbox{ \muL} \mbox{ min}^{-1} \\ Column \mbox{ temperature } = 40 \ ^{\circ}\text{C} \\ Injection \mbox{ volume } = 1 \mbox{ \muL} \\ Mobile \mbox{ phase: } 140 \mbox{ mM} \mbox{ potassium phosphate, 8 \mbox{ mM} \mbox{ KCl, } 0.1 \mbox{ mM} \mbox{ EDTA, 6 \mbox{ mM} \mbox{ OSA.} \\ MeOH/H_2O \mbox{ (} 6/94; v/v) \\ \mbox{ pH} = 5 \mbox{ (adjusted with 10 \mbox{ mM} \mbox{ NaOH}) \\ Detection: \mbox{ GC working electrode (+0.45 \mbox{ V vs.} \mbox{ Ag/AgCl}) \end{array}$			Development of analytical method for the determination of monoamines concentrations and application to microdialysis in the dorsal hippocampus		
	Detection - (mol L <sup>-1</sup> )	5-HT DA		3-MT			
		LOD $1.5 \times 10^{-9}$ $7.5 \times 10^{-10}$		$1.5 imes 10^{-9}$			
	. ,	LOQ	LOQ $5.0 \times 10^{-9}$ $2.5 \times 10^{-9}$		$5.0 imes10^{-9}$		

**Table 2.** Recent values of separation, detection and quantification of compounds involved in serotonergic system and their application (non-exhaustive list). All values provided have been standardized to the same unit (mol  $L^{-1}$ ) for the purpose of comparability.

Reference			Application					
Schou- Pedersen et al. (2016) [64]	Methods and Ap- plication	$\label{eq:heat} \begin{array}{l} \mbox{HPLC-ECD Conditions} \\ C_{18} \mbox{ column } (4.6 \times 100 \mbox{ mm}, 2.6 \mbox{ µm particles size}) \\ \mbox{Flow rate } = 0.8 \mbox{ mL min}^{-1} \\ \mbox{Column temperature } = 30 \ ^{\circ}\text{C} \\ \mbox{Injection volume } = 20 \mbox{ µLMobile phase: 70 \mbox{ mM potassium dihydrogen phosphate, 2 \mbox{ mM OSA,}} \\ \mbox{0.1 \mbox{ mM EDTA.}} \\ \mbox{MeOH}/H_2 O \ (10/90; v/v) \\ \mbox{pH = 3.12 (adjusted with 1 \mbox{ mctric acid})} \\ \mbox{Detection: Porous graphite working electrode} \\ \mbox{(+0.40 \mbox{ V vs. Pd)} \end{array}$				Development of chromatographic method for the quantification of monoamine NTs from sub-regions of guinea pig brain (intracellular and extracellular).		
	Detection (mol L <sup>-1</sup> )		5-HT	5-HIAA	DOPAC	DA	HVA	
		$LOQ \qquad 8.8\times10^{-9}$		$3.8 imes10^{-9}$	$3.6 imes10^{-9}$	$1.0 imes10^{-8}$	$1.2  imes 10^{-8}$	
		L.R.	$8.8\times 10^{-9}1.0\times 10^{-6}$	$3.8  imes 10^{-9}$ - $1.0  imes 10^{-6}$	$\begin{array}{c} 3.6 \times 10^{-9}  1.0 \\ \times 10^{-6} \end{array}$	$\begin{array}{c} 1.0 \times 10^{-8} 7.5 \\ \times 10^{-7} \end{array}$	$1.2 imes 10^{-8}$ -5.0 $ imes 10^{-7}$	
Van Dam et al. (2014) [65]	Methods and Ap- plication		$\begin{array}{c} C_{18} \mbox{ column } (1.0 \times 250 \mbox{ mm}, 3.0 \mbox{ \mum} \mbox{ particles size}) \\ Flow rate = 40  \mu L \mbox{ min}^{-1} \\ Column \mbox{ temperature } = \mbox{ from } 30 ^{\circ} C \mbox{ to } 36 ^{\circ} C \\ Injection \mbox{ volume } = 50  \mu L \\ Mobile \mbox{ phase: } 8 \mbox{ mM } KCl, 50 \mbox{ mM} \mbox{ phosphoric } acid, 50 \mbox{ mM} \mbox{ citric } acid, 0.1 \mbox{ mM } EDTA, \mbox{ from } 1.8 \mbox{ to } 2.2 \mbox{ mM } OSA. \\ MeOH/H_2O \mbox{ (from } 13/87 \mbox{ to } 17/83;  v/v) \\ \mbox{ pH } = \mbox{ from } 3.0 \mbox{ to } 3.6 \mbox{ (adjusted with } NaOH) \\ Detection: \mbox{ GC working electrode } (from +0.63 \mbox{ V to } +0.67 \mbox{ V vs. } Ag/AgCl) \end{array}$				Development of analytical method for the quantification of biogenic amines and metabolites in human brain tissue.	
			5-HT	5-HIAA	DOPAC	DA	HVA	
	Detection (mol L <sup>-1</sup> )	L.R.	$5.3 \times 10^{-10}$ - $1.8 \times 10^{-7}$	$2.2  imes 10^{-10}$ -7.7 $ imes 10^{-8}$	$\begin{array}{c} \textbf{2.8}\times10^{-10} \textbf{-1.6} \\ \times10^{-7} \end{array}$	$\begin{array}{r} 2.7 \times 10^{-10}  1.8 \\ \times 10^{-7} \end{array}$	$5.8  imes 10^{-10}$ -3.2 $ imes 10^{-7}$	
	Methods and Ap- plication		C <sub>18</sub> column (3.0 × 5 Flow ra Column f n volume = 25 μLMobile phase horic acid, 2 mM sodium dodec pH = 3.2 (adjus Detection: GC working	Development of analytical method for the quantification of NTs and metabolites in brain mice microdialysates.				
	Detection (mol L <sup>-1</sup> )	5-HT 5-HIAA DA		DA	HVA			
		LOD	LOD $3.1 \times 10^{-10}$ $5.4 \times 10^{-10}$ $3.7 \times 10^{-10}$		$4.1 imes 10^{-10}$			
		LOQ	LOQ $3.9 \times 10^{-10}$ $5.8 \times 10^{-10}$ $4.2 \times 10^{-10}$				$4.3 imes10^{-10}$	
Allen et al.	Methods and Ap- plication	Injection v	C <sub>18</sub> column (3.2 × Flow ra Column rolume = 40 μLMobile phase: sodium octyl sodium octyl ACN, pH = 3.3 (adjust Detection: Dual workin	Development of analytical method for the quantification of monoamines and metabolites in brain tissue of mice.				
(2017) [67]			5-HT	5-HIAA	DOPAC	DA	HVA	
	Detection	LOD	$1.7 imes10^{-9}$	$6.5 imes10^{-10}$	$7.4 imes10^{-10}$	$8.2  imes 10^{-10}$	$4.1  imes 10^{-10}$	
	$(mol L^{-1})$	LOQ	$3.6 imes10^{-9}$	$1.3  imes 10^{-9}$	$1.5  imes 10^{-9}$	$1.6 imes10^{-9}$	$1.4  imes 10^{-9}$	
		L.R.	$3.6  imes 10^{-9}$ - $1.7  imes 10^{-4}$	$1.3  imes 10^{-9}$ -7.9 $ imes 10^{-5}$	$\begin{array}{c} 1.5 \times 10^{-9}  8.9 \\ \times 10^{-5} \end{array}$	$\begin{array}{c} 1.6 \times 10^{-9}  2.0 \\ \times 10^{-4} \end{array}$	$1.4 imes 10^{-9}$ -8.2 $ imes 10^{-5}$	
Yardimci et al. (2023) [68] –	Methods and Ap- plication	Mobile pha	$\begin{array}{c} C_{18} \mbox{ column } (4.6 \times 10^{-3} \mbox{ Flow r} \\ \mbox{ Flow r} \\ \mbox{ Column i} \\ \mbox{ Injectio} \\ \mbox{ se: 35 mM citric acid, 19 mM s} \\ \mbox{ acid. Glacial acetic aci} \\ \mbox{ (0.11/0.3/ \mbox{ pH} = 4.9 (aci) \mbox{ pH} \\ \mbox{ Detection: GC working} \end{array}$	Analysis in hypothalamic and subcortical nuclei				
	Detection	5-HT 5-HIAA DA			DOPAC			
	Detection (mol L <sup>-1</sup> )		5-H1	5-HIAA	DA	DO	PAC	

#### Table 2. Cont.

Reference		HPLC-ECD Conditions					Application		
Du et al. (2018) [69] -	Methods and Ap- plication	$C_{18} \text{ column } (4.6 \times 250 \text{ mm, } 5.0  \mu\text{m particles size})$ $Flow \text{ rate } = 1 \text{ mL min}^{-1}$ $Column \text{ temperature } = 25 ^{\circ}\text{C}$ $Injection \text{ volume } = 20  \mu\text{L}$ $Mobile \text{ phase: } 25 \text{ mM sodium citrate, } 0.01 \text{ mM EDTA.}$ $ACN/H_2O (5/95; v/v)$ $pH = 4.5 \text{ (adjusted with } 1 \text{ M acetic acid})$ $Detection: BDD \text{ working electrode } (+0.70 \text{ V vs. } \text{Ag/AgCl})$				Analytical method development			
	Detection $\_$ (mol L <sup>-1</sup> )	5-HT					5-HIAA		
		LOD $2.1 \times 10^{-8}$					$1.6 imes 10^{-8}$		
		LOQ $2.8 \times 10^{-8}$				$4.2 imes10^{-8}$			
	_	L.R.		$\textbf{2.8}\times10^{-8}\textbf{-1.1}$	$\times 10^{-6}$		$2.6 imes 10^{-8}$ – $2.6 imes 10^{-6}$		
Zhang et al. (2016) [53]	Methods and Ap- plication	C <sub>18</sub> column (4.6 × 250 mm, 5.0 µm particles size) Flow rate = 1 mL min <sup>-1</sup> Column temperature = 25 °C Injection volume = 20 µL Mobile phase: 25 mM sodium citrate, 0.01 mM EDTA. ACN/H <sub>2</sub> O (5/95; $v/v$ ) pH = 4.5 (adjusted with 1 M acetic acid) Detection: BDD working electrode (+0.70 V vs. Ag/AgCl)			Determination of the concentrations of monoamines NTs in rat cortex and hippocampus tissues.				
	Detection – (mol L <sup>-1</sup> )		5-HT	5-HIAA	DA	DOPAC	3-MT	HVA	
		LOD	$2.3 imes10^{-8}$	$1.1 imes 10^{-8}$	$2.6 imes10^{-8}$	$1.2  imes 10^{-8}$	$3.6 imes10^{-8}$	$2.7 imes10^{-8}$	
		LOQ	$8.5 imes10^{-8}$	$3.1 imes10^{-8}$	$9.8 imes10^{-8}$	$4.8 imes10^{-8}$	$1.2  imes 10^{-7}$	$8.2 imes10^{-8}$	
		L.R.	$8.5 \times 10^{-8}  1.4 \times 10^{-6}$	$\begin{array}{c} 3.1 \times 10^{-8} 7.9 \\ \times 10^{-7} \end{array}$	$\begin{array}{c} 9.8 \times 10^{-8}  2.3 \\ \times 10^{-6} \end{array}$	${\begin{array}{r}{6.0 imes10^{-8}}{ m -3.0}\  imes10^{-6}\end{array}}$	$\begin{array}{c} 1.2 \times 10^{-7}  1.8 \\ \times 10^{-6} \end{array}$	$8.2  imes 10^{-8}$ -1.4 $ imes 10^{-6}$	
Jiang et al. (2015) [70]	Methods and Ap- plication	$\begin{array}{c} C_{18} \mbox{ column } (4.6 \times 250 \mbox{ mm}, 5.0 \mbox{ \mum particles size}) \\ \mbox{ Flow rate = 1 mL min^{-1}} \\ \mbox{ Column temperature = 30 °C} \\ \mbox{ Injection volume = 10  \mu L} \\ \mbox{ Mobile phase: 50 mM potassium dihydrogen phosphate, 0.1 mM octane sulfonic acid.} \\ \mbox{ MeOH/H}_2O (5/95; v/v) \\ \mbox{ pH = N/A} \\ \mbox{ Detection: GC working electrode (+0.70 V vs. Ag/AgCl)} \end{array}$					Development of the analytical method for monoamines NTs in human urine.		
	Detection	5-HT				DA			
	$(mol L^{-1})$		LOD	$6.1 imes10^{-8}$			$3.9 imes10^{-8}$		
Lokhande et al. (2022) [71]	Methods and Ap- plication	$\begin{array}{c} C_{18} \mbox{ column } (4.6 \times 250 \mbox{ mm}, 5.0 \mbox{ \mum} \mbox{ particles size}) \\ Flow rate = 1.3 \mbox{ mL} \mbox{ min}^{-1} \\ Column \mbox{ temperature} = 35 \mbox{ °C} \\ Injection \mbox{ volume} = 20 \mbox{ \muL} \\ Mobile \mbox{ phase: } 50 \mbox{ mM} \mbox{ potassium dihydrogen phosphate, } 0.99 \mbox{ mM} \mbox{ SOS and } 53 \mbox{ \muM} \mbox{ EDTA}. \\ MeOH/H_2O \mbox{ (12/88; $v'v)} \\ pH = 2.5 \mbox{ (adjusted with } 85\% \mbox{ phosphoric acid}) \\ Detection: \mbox{ BDD working electrode (+0.70 \mbox{ V vs. } Ag/AgCl)} \end{array}$					Development of the analytical method for metabolites quantification in human CSF.		
	Detection $-$ (mol L <sup>-1</sup> )	5-HIAA			HVA				
			LOQ $6.5 \times 10^{-8}$			$6.9 imes10^{-8}$			
	· / _	L.R. $6.5 \times 10^{-8}$ -2.6 × 10 <sup>-6</sup>				$6.9 imes 10^{-8}$ – $2.7 imes 10^{-6}$			

#### Table 2. Cont.

Abreviations: N/A: not available, LOD: limit of detection, LOQ: limit of quantification, L.R.: linear range, EDTA: ethylendiaminetetraacetic acid, OSA: octane sulfonic acid, SOS: sodium octane sulfonate, ACN: acetonitrile, MeOH: methanol.

#### 4. Conclusions

To study the brain serotonergic and dopaminergic systems, researchers employ intracerebral microdialysis coupled with high-performance liquid chromatography (HPLC) and electrochemical detector (ECD). HPLC, especially in reverse phase conditions, is a widely used separation technique for NT analysis, offering rapid and selective results. Electrochemical detection, particularly with glassy carbon electrodes, proves to be a sensitive and cost-effective method for quantifying NTs and their metabolites with values close to  $10^{-9}$  mol L<sup>-1</sup> for HPLC and  $10^{-11}$  mol L<sup>-1</sup> with UHPLC.

Optimizations, such as increasing the surface-to-volume ratio of the electrochemical cell and modifying electrode materials can enhance sensitivity. The determination of limits of detection and quantification, as well as the linear range, is crucial for evaluating the analytical performance of the HPLC-ECD system. Overall, the combination of microdialysis,

HPLC, and ECD provides a powerful approach for understanding the dynamics of NTs and their metabolites, contributing to advancements in neuroscientific research and potential clinical applications.

It should be noted, however, that the use of ECD for neurotransmitters is becoming less common. Indeed, mass spectrometry is becoming more affordable, and the preventive maintenance of the analyzer is becoming simpler. Consequently, due to its high sensitivity and selectivity, mass spectrometry is currently the most common detector in analysis laboratories.

Research perspectives in the field of high-performance liquid chromatography with an electrochemical detector offer a fertile ground for improving the sensitivity and specificity of neurotransmitter detection. One promising avenue involves exploring the modification of the electrochemical detector's surface. By altering surface properties, such as roughness or chemical composition, one can optimize the interaction with neurotransmitters, thereby enhancing detection efficiency.

Another research direction focuses on modifying the characteristics of the electrochemical cell itself, including its surface and dimensions. By adjusting these parameters, the signal-to-noise ratio can be increased, which is crucial for precise measurements. Electrochemical cells with an increased surface area can enhance the likelihood of interaction with analytes, thereby improving detection. Moreover, optimized dimensions can facilitate a clearer electrochemical response and reduce unwanted interferences.

By combining these approaches, it becomes conceivable to achieve quantification limits comparable to those obtained with mass spectrometry detectors, paving the way for significant advancements in understanding neurochemical mechanisms. These research endeavors hold the promise of bringing substantial improvements to neurotransmitter detection, opening new perspectives for understanding neurological disorders and mental illnesses.

**Author Contributions:** B.P.G.: review design and writing, G.G.: referencing and writing. All authors have read and agreed to the published version of the manuscript.

**Funding:** This research was funded by the Agence Nationale pour la recherche, grant number AIS-in-DEP (ANR-22-CE37-0018 Neuroscience Intégrative et Cognitive).

Institutional Review Board Statement: Not applicable.

Informed Consent Statement: Not applicable.

**Data Availability Statement:** No new data were created or analyzed in this study. Data sharing is not applicable to this article.

Conflicts of Interest: The authors declare no conflicts of interest.

#### References

- 1. Berger, M.; Gray, J.A.; Roth, B.L. The Expanded Biology of Serotonin. Annu. Rev. Med. 2009, 60, 355–366. [CrossRef]
- Factor, S.A.; McDonald, W.M.; Goldstein, F.C. The Role of Neurotransmitters in the Development of Parkinson's Disease-related Psychosis. *Eur. J. Neurol.* 2017, 24, 1244–1254. [CrossRef] [PubMed]
- Tajeddinn, W.; Fereshtehnejad, S.; Ahmed, M.; Yoshitake, T.; Kehr, J.; Shahnaz, T.; Milovanovic, M.; Behbahani, H.; Hoglund, K.; Winblad, B.; et al. Association of Platelet Serotonin Levels in Alzheimer's Disease with Clinical and Cerebrospinal Fluid Markers. J. Alzheimer Dis. 2016, 53, 621–630. [CrossRef]
- 4. Xu, Y.; Yan, J.; Zhou, P.; Li, J.; Gao, H.; Xia, Y.; Wang, Q. Neurotransmitter Receptors and Cognitive Dysfunction in Alzheimer's Disease and Parkinson's Disease. *Prog. Neurobiol.* **2012**, *97*, 1–13. [CrossRef] [PubMed]
- Bahi, A.; Dreyer, J.-L. Anxiety and Ethanol Consumption in Socially Defeated Mice; Effect of Hippocampal Serotonin Transporter Knockdown. *Behav. Brain Res.* 2023, 451, 114508. [CrossRef]
- Kraus, C.; Castrén, E.; Kasper, S.; Lanzenberger, R. Serotonin and Neuroplasticity—Links between Molecular, Functional and Structural Pathophysiology in Depression. *Neurosci. Biobehav. Rev.* 2017, 77, 317–326. [CrossRef]
- Yagishita, S. Transient and Sustained Effects of Dopamine and Serotonin Signaling in Motivation-related Behavior. *Psychiatry Clin. Neurosci.* 2020, 74, 91–98. [CrossRef] [PubMed]
- 8. Speranza, L.; Di Porzio, U.; Viggiano, D.; De Donato, A.; Volpicelli, F. Dopamine: The Neuromodulator of Long-Term Synaptic Plasticity, Reward and Movement Control. *Cells* **2021**, *10*, 735. [CrossRef]

- 9. Weinstein, J.J.; Chohan, M.O.; Slifstein, M.; Kegeles, L.S.; Moore, H.; Abi-Dargham, A. Pathway-Specific Dopamine Abnormalities in Schizophrenia. *Biol. Psychiatry* 2017, *81*, 31–42. [CrossRef]
- Di Forti, M.; Lappin, J.M.; Murray, R.M. Risk Factors for Schizophrenia—All Roads Lead to Dopamine. *Eur. Neuropsychopharmacol.* 2007, 17, S101–S107. [CrossRef]
- 11. Latif, S.; Jahangeer, M.; Maknoon Razia, D.; Ashiq, M.; Ghaffar, A.; Akram, M.; El Allam, A.; Bouyahya, A.; Garipova, L.; Ali Shariati, M.; et al. Dopamine in Parkinson's Disease. *Clin. Chim. Acta* **2021**, *522*, 114–126. [CrossRef] [PubMed]
- 12. Wong, D.L. Epinephrine Biosynthesis: Hormonal and Neural Control During Stress. *Cell. Mol. Neurobiol.* **2006**, *26*, 889–898. [CrossRef]
- 13. Harrison, N.A.; Morgan, R.; Critchley, H.D. From Facial Mimicry to Emotional Empathy: A Role for Norepinephrine? *Soc. Neurosci.* **2010**, *5*, 393–400. [CrossRef]
- 14. Alves, E.; Lukoyanov, N.; Serrão, P.; Moura, D.; Moreira-Rodrigues, M. Epinephrine Increases Contextual Learning through Activation of Peripheral B2-Adrenoceptors. *Psychopharmacology* **2016**, *233*, 2099–2108. [CrossRef] [PubMed]
- 15. Ferry, B.; Roozendaal, B.; McGaugh, J.L. Role of Norepinephrine in Mediating Stress Hormone Regulation of Long-Term Memory Storage: A Critical Involvement of the Amygdala. *Biol. Psychiatry* **1999**, *46*, 1140–1152. [CrossRef] [PubMed]
- 16. Tossman, U.; Ungerstedt, U. Microdialysis in the Study of Extracellular Levels of Amino Acids in the Rat Brain. *Acta Physiol. Scand.* **1986**, *128*, 9–14. [CrossRef]
- 17. Ungerstedt, U.; Hallström, Å. In Vivo Microdialysis-a New Approach to the Analysis of Neurotransmitters in the Brain. *Life Sci.* **1987**, *41*, 861–864. [CrossRef] [PubMed]
- Zhou, S.Y.; Zuo, H.; Stobaugh, J.F.; Lunte, C.E.; Lunte, S.M. Continuous in Vivo Monitoring of Amino Acid Neurotransmitters by Microdialysis Sampling with Online Derivatization and Capillary Electrophoresis Separation. *Anal. Chem.* 1995, 67, 594–599. [CrossRef]
- Van Schoors, J.; Lens, C.; Maes, K.; Michotte, Y.; Smolders, I.; Van Eeckhaut, A. Reassessment of the Antioxidative Mixture for the Challenging Electrochemical Determination of Dopamine, Noradrenaline and Serotonin in Microdialysis Samples. *J. Chromatogr.* B 2015, 998–999, 63–71. [CrossRef]
- Guiard, B.; Lanfumey, L.; Gardier, A. Microdialysis Approach to Study Serotonin Outflow in Mice Following Selective Serotonin Reuptake Inhibitors and Substance P (Neurokinin 1) Receptor Antagonist Administration: A Review. *Curr. Drug Targets* 2006, 7, 187–201. [CrossRef]
- Zhang, P.; Wang, B.; Sun, Y.; Gao, J.; Lian, K. Analysis of 5-Hydroxytryptamine and Its Related Indoles in Cerebrospinal Fluid of Leukemic Children by Gas Chromatography-Mass Spectrometry. J. Lab. Med. 2020, 44, 41–45. [CrossRef]
- Shi, H.; Wang, B.; Niu, L.; Cao, M.; Kang, W.; Lian, K.; Zhang, P. Trace Level Determination of 5-Hydroxytryptamine and Its Related Indoles in Amniotic Fluid by Gas Chromatography–Mass Spectrometry. *J. Pharm. Biomed. Anal.* 2017, 143, 176–182. [CrossRef] [PubMed]
- Tsai, H.; Whang, C. Capillary Electrophoresis of Monoamines and Catechol with Indirect Chemiluminescence Detection. *Electrophoresis* 1999, 20, 2533–2538. [CrossRef]
- Benturquia, N.; Couderc, F.; Sauvinet, V.; Orset, C.; Parrot, S.; Bayle, C.; Renaud, B.; Denoroy, L. Analysis of Serotonin in Brain Microdialysates Using Capillary Electrophoresis and Native Laser-Induced Fluorescence Detection. *Electrophoresis* 2005, 26, 1071–1079. [CrossRef] [PubMed]
- 25. Hattox, S.E.; Murphy, R.C. Mass Spectrometry and Gas Chromatography of Trimethylsilyl Derivatives of Catecholamine Related Molecules. *Biol. Mass Spectrom.* **1978**, *5*, 338–345. [CrossRef] [PubMed]
- 26. Donzanti, B.A.; Yamamoto, B.K. An Improved and Rapid HPLC-EC Method for the Isocratic Separation of Amino Acid Neurotransmitters from Brain Tissue and Microdialysis Perfusates. *Life Sci.* **1988**, *43*, 913–922. [CrossRef]
- 27. Peinado, J.M.; McManus, K.T.; Myers, R.D. Rapid Method for Micro-Analysis of Endogenous Amino Acid Neurotransmitters in Brain Perfusates in the Rat by Isocratic HPLC-EC. J. Neurosci. Methods **1986**, *18*, 269–276. [CrossRef]
- Rogers, K.L.; Philibert, R.A.; Allen, A.J.; Molitor, J.; Wilson, E.J.; Dutton, G.R. HPLC Analysis of Putative Amino Acid Neurotransmitters Released from Primary Cerebellar Cultures. J. Neurosci. Methods 1987, 22, 173–179. [CrossRef]
- Bidel, F.; Corvaisier, S.; Jozet-Alves, C.; Pottier, I.; Dauphin, F.; Naud, N.; Bellanger, C. An HPLC-ECD Method for Monoamines and Metabolites Quantification in Cuttlefish (Cephalopod) Brain Tissue: Biogenic Monoamines and Metabolites in Cuttlefish (Cephalopod) Brain. *Biomed. Chromatogr.* 2016, *30*, 1175–1183. [CrossRef]
- Reinhoud, N.J.; Brouwer, H.-J.; Van Heerwaarden, L.M.; Korte-Bouws, G.A.H. Analysis of Glutamate, GABA, Noradrenaline, Dopamine, Serotonin, and Metabolites Using Microbore UHPLC with Electrochemical Detection. ACS Chem. Neurosci. 2013, 4, 888–894. [CrossRef]
- Tufi, S.; Lamoree, M.; De Boer, J.; Leonards, P. Simultaneous Analysis of Multiple Neurotransmitters by Hydrophilic Interaction Liquid Chromatography Coupled to Tandem Mass Spectrometry. J. Chromatogr. A 2015, 1395, 79–87. [CrossRef] [PubMed]
- 32. Chirita, R.-I.; West, C.; Finaru, A.-L.; Elfakir, C. Approach to Hydrophilic Interaction Chromatography Column Selection: Application to Neurotransmitters Analysis. *J. Chromatogr. A* 2010, 1217, 3091–3104. [CrossRef]
- Zhou, G.-S.; Yuan, Y.-C.; Yin, Y.; Tang, Y.-P.; Xu, R.-J.; Liu, Y.; Chen, P.-D.; Yin, L.; Duan, J.-A. Hydrophilic Interaction Chromatography Combined with Ultrasound-Assisted Ionic Liquid Dispersive Liquid–Liquid Microextraction for Determination of Underivatized Neurotransmitters in Dementia Patients' Urine Samples. *Anal. Chim. Acta* 2020, 1107, 74–84. [CrossRef]

- Boulghobra, A.; Bonose, M.; Billault, I.; Pallandre, A. A Rapid and Sensitive Method for the Quantification of Dopamine and Serotonin Metabolites in Cerebrospinal Fluid Based on UHPLC with Fluorescence Detection. J. Chromatogr. B 2022, 1200, 123264. [CrossRef] [PubMed]
- 35. De Benedetto, G.E.; Fico, D.; Pennetta, A.; Malitesta, C.; Nicolardi, G.; Lofrumento, D.D.; De Nuccio, F.; La Pesa, V. A Rapid and Simple Method for the Determination of 3,4-Dihydroxyphenylacetic Acid, Norepinephrine, Dopamine, and Serotonin in Mouse Brain Homogenate by HPLC with Fluorimetric Detection. J. Pharm. Biomed. Anal. 2014, 98, 266–270. [CrossRef]
- 36. Zhao, H.-X.; Mu, H.; Bai, Y.-H.; Yu, H.; Hu, Y.-M. A Rapid Method for the Determination of Dopamine in Porcine Muscle by Pre-Column Derivatization and HPLC with Fluorescence Detection. *J. Pharm. Anal.* **2011**, *1*, 208–212. [CrossRef] [PubMed]
- 37. Carrera, V.; Sabater, E.; Vilanova, E.; Sogorb, M.A. A Simple and Rapid HPLC–MS Method for the Simultaneous Determination of Epinephrine, Norepinephrine, Dopamine and 5-Hydroxytryptamine: Application to the Secretion of Bovine Chromaffin Cell Cultures. *J. Chromatogr. B* 2007, *8*47, 88–94. [CrossRef]
- Xu, N.; Qiu, C.; Wang, W.; Wang, Y.; Chai, C.; Yan, Y.; Zhu, D. HPLC/MS/MS for Quantification of Two Types of Neurotransmitters in Rat Brain and Application: Myocardial Ischemia and Protection of Sheng-Mai-San. J. Pharm. Biomed. Anal. 2011, 55, 101–108. [CrossRef]
- Kovac, A.; Somikova, Z.; Zilka, N.; Novak, M. Liquid Chromatography–Tandem Mass Spectrometry Method for Determination of Panel of Neurotransmitters in Cerebrospinal Fluid from the Rat Model for Tauopathy. *Talanta* 2014, 119, 284–290. [CrossRef]
- 40. Shao, Z.; Chang, Y.; Venton, B.J. Carbon Microelectrodes with Customized Shapes for Neurotransmitter Detection: A Review. *Anal. Chim. Acta* 2022, 1223, 340165. [CrossRef]
- 41. Polo, E.; Kruss, S. Nanosensors for Neurotransmitters. Anal. Bioanal. Chem. 2016, 408, 2727–2741. [CrossRef] [PubMed]
- 42. Banerjee, S.; McCracken, S.; Hossain, M.F.; Slaughter, G. Electrochemical Detection of Neurotransmitters. *Biosensors* 2020, 10, 101. [CrossRef]
- Dicgory, G.L.; Buckett, W.R. An Automated Method to Measure Monoamines and Metabolites Using Elevated Temperature Reversed Phase HPLC with Electrochemical Detection Application to Striatal Dopamine and Hippocampal Serotonin Turnover. J. Pharmacol. Methods 1984, 11, 207–217. [CrossRef] [PubMed]
- Cannazza, G.; Di Stefano, A.; Mosciatti, B.; Braghiroli, D.; Baraldi, M.; Pinnen, F.; Sozio, P.; Benatti, C.; Parenti, C. Detection of Levodopa, Dopamine and Its Metabolites in Rat Striatum Dialysates Following Peripheral Administration of I-DOPA Prodrugs by Mean of HPLC–EC. J. Pharm. Biomed. Anal. 2005, 36, 1079–1084. [CrossRef] [PubMed]
- 45. Birch, P.J.; Fillenz, M. Measurement of Noradrenaline Synthesis in Rat Hippocampal Synaptosomes Using HPLC with ECD. J. *Neurosci. Methods* **1985**, *13*, 231–238. [CrossRef]
- Saito, H.; Murai, S.; Abe, E.; Masuda, Y.; Itoh, T. Rapid and Simultaneous Assay of Monoamine Neurotransmitters and Their Metabolites in Discrete Brain Areas of Mice by HPLC with Coulometric Detection. *Pharmacol. Biochem. Behav.* 1992, 42, 351–356. [CrossRef]
- 47. Bard, A.J.; Faulkner, L.R. *Electrochemical Methods, Fundamentals and Applications*; John Wiley & Sons Inc.: New York, NY, USA, 2001; ISBN 0-471-04372-9.
- 48. Sharma, S.; Singh, N.; Tomar, V.; Chandra, R. A Review on Electrochemical Detection of Serotonin Based on Surface Modified Electrodes. *Biosens. Bioelectron.* 2018, 107, 76–93. [CrossRef]
- 49. Błaszczyk, J.W. Parkinson's Disease and Neurodegeneration: GABA-Collapse Hypothesis. *Front. Neurosci.* **2016**, *10*, 239. [CrossRef]
- 50. Nägga, K.; Bogdanovic, N.; Marcusson, J. GABA Transporters (GAT-1) in Alzheimer's Disease. J. Neural Transm. 1999, 106, 1141–1149. [CrossRef]
- 51. Panrod, K.; Tansirikongkol, A.; Panapisal, V. Comparison of Validated High-Performance Liquid Chromatography Methods Using Two Derivatizing Agents for Gamma-Aminobutyric Acid Quantification. *Thai J. Pharm. Sci.* **2016**, *40*, 203–208.
- Narmadha, M.; Noel, M.; Suryanarayanan, V. Relative Deactivation of Boron-Doped Diamond (BDD) and Glassy Carbon (GC) Electrodes in Different Electrolyte Media Containing Substituted Phenols—Voltammetric and Surface Morphologic Studies. J. Electroanal. Chem. 2011, 655, 103–110. [CrossRef]
- Zhang, L.; Yang, J.; Luo, Y.; Shang, J.; Jiang, X. Simultaneous Determination of Eleven Compounds Related to Metabolism of Bioamines in Rat Cortex and Hippocampus by HPLC-ECD with Boron-Doped Diamond Working Electrode. *J. Pharm. Biomed. Anal.* 2016, 118, 41–51. [CrossRef] [PubMed]
- 54. Gotti, G.; Evrard, D.; Gros, P. Simultaneous Electrochemical Detection of Oxygen (O<sub>2</sub>) and Hydrogen Peroxide (H<sub>2</sub>O<sub>2</sub>) in Neutral Media. *Int. J. Electrochem. Sci.* 2023, *18*, 100262. [CrossRef]
- 55. Gotti, G.; Fajerwerg, K.; Evrard, D.; Gros, P. Kinetics of Dioxygen Reduction on Gold and Glassy Carbon Electrodes in Neutral Media. *Int. J. Electrochem. Sci.* 2013, *8*, 12643–12657. [CrossRef]
- Xu, F.; Gao, M.; Shi, G.; Wang, L.; Zhang, W.; Xue, J.; Jin, L.; Jin, J. Simultaneous Detection of Monoamines in Rat Striatal Microdialysate at Poly(Para-Aminobenzoic Acid) Modified Electrode by High-Performance Liquid Chromatography. *Anal. Chim. Acta* 2001, 439, 239–246. [CrossRef]
- Guo, Y.-X.; Xia, C.-Y.; Yan, Y.; Han, Y.; Shi, R.; He, J.; Wang, Y.-M.; Wang, Z.-X.; Zhang, W.-K.; Xu, J.-K. Loganin Improves Chronic Unpredictable Mild Stress-Induced Depressive-like Behaviors and Neurochemical Dysfunction. *J. Ethnopharmacol.* 2023, 308, 116288. [CrossRef] [PubMed]

- 58. Lei, C.; Chen, H.; Chen, K. Effects of Bone Marrow Mesenchymal Stem Cells on Myelin Repair and Emotional Changes of a Cuprizone-Induced Demyelination Model. *J. Integr. Neurosci.* **2023**, *22*, 40. [CrossRef]
- Wang, Z.-X.; Lian, W.-W.; He, J.; He, X.-L.; Wang, Y.-M.; Pan, C.-H.; Li, M.; Zhang, W.-K.; Liu, L.-Q.; Xu, J.-K. Cornuside Ameliorates Cognitive Impairments in Scopolamine Induced AD Mice: Involvement of Neurotransmitter and Oxidative Stress. J. Ethnopharmacol. 2022, 293, 115252. [CrossRef]
- 60. Kalantary, R.R.; Moradi, M.; Pirsaheb, M.; Esrafili, A.; Jafari, A.J.; Gholami, M.; Vasseghian, Y.; Antolini, E.; Dragoi, E.-N. Enhanced Photocatalytic Inactivation of E. Coli by Natural Pyrite in Presence of Citrate and EDTA as Effective Chelating Agents: Experimental Evaluation and Kinetic and ANN Models. *J. Environ. Chem. Eng.* **2019**, *7*, 102906. [CrossRef]
- 61. Broch, S.C.; Celia García Alvarez-Coque, M.; Broch, S.C.; Esteve-Romero, J.S. Liquid Chromatographic Determination of Some Thiazide Diuretics in Pharmaceuticals with a Sodium Dodecyl Sulfate Mobile Phase. *Analyst* **1998**, *123*, 301–306. [CrossRef]
- Umek, N.; Geršak, B.; Vintar, N.; Šoštarič, M.; Mavri, J. Dopamine Autoxidation Is Controlled by Acidic pH. Front. Mol. Neurosci. 2018, 11, 467. [CrossRef] [PubMed]
- 63. Ferry, B.; Gifu, E.-P.; Sandu, I.; Denoroy, L.; Parrot, S. Analysis of Microdialysate Monoamines, Including Noradrenaline, Dopamine and Serotonin, Using Capillary Ultra-High Performance Liquid Chromatography and Electrochemical Detection. *J. Chromatogr. B* **2014**, *951–952*, 52–57. [CrossRef] [PubMed]
- 64. Schou-Pedersen, A.M.V.; Hansen, S.N.; Tveden-Nyborg, P.; Lykkesfeldt, J. Simultaneous Quantification of Monoamine Neurotransmitters and Their Biogenic Metabolites Intracellularly and Extracellularly in Primary Neuronal Cell Cultures and in Sub-Regions of Guinea Pig Brain. J. Chromatogr. B 2016, 1028, 222–230. [CrossRef] [PubMed]
- Van Dam, D.; Vermeiren, Y.; Aerts, T.; De Deyn, P.P. Novel and Sensitive Reversed-Phase High-Pressure Liquid Chromatography Method with Electrochemical Detection for the Simultaneous and Fast Determination of Eight Biogenic Amines and Metabolites in Human Brain Tissue. J. Chromatogr. A 2014, 1353, 28–39. [CrossRef] [PubMed]
- 66. Pantiya, P.; Guiard, B.P.; Gotti, G. Sensitive and Fast Detection of Monoamines and Their Metabolites by High-Performance Liquid Chromatography Coupled with an Electrochemical Detector (HPLC-ECD) under Isocratic Conditions. Application to Intracerebral Microdialysis in Mice Treated by Fluoxetine and Atomoxetine. *Chromatographia* **2024**, *in press*.
- Allen, S.A.; Rednour, S.; Shepard, S.; Pond, B.B. A Simple and Sensitive High-performance Liquid Chromatography– Electrochemical Detection Assay for the Quantitative Determination of Monoamines and Respective Metabolites in Six Discrete Brain Regions of Mice. *Biomed. Chromatogr.* 2017, *31*, e3998. [CrossRef]
- 68. Yardimci, A.; Ertugrul, N.U.; Ozgen, A.; Ozbeg, G.; Ozdede, M.R.; Ercan, E.C.; Canpolat, S. Effects of Chronic Irisin Treatment on Brain Monoamine Levels in the Hypothalamic and Subcortical Nuclei of Adult Male and Female Rats: An HPLC-ECD Study. *Neurosci. Lett.* **2023**, *806*, 137245. [CrossRef]
- 69. Du, T.; Cui, T.; Qiu, H.; Wang, N.; Huang, D.; Jiang, X. Simultaneous Determination of Tryptophan, Kynurenine, Kynurenic Acid and Two Monoamines in Rat Plasma by HPLC-ECD/DAD. *J. Pharm. Biomed. Anal.* **2018**, *158*, 8–14. [CrossRef]
- Jiang, L.; Chen, Y.; Chen, Y.; Ma, M.; Tan, Y.; Tang, H.; Chen, B. Determination of Monoamine Neurotransmitters in Human Urine by Carrier-Mediated Liquid-Phase Microextraction Based on Solidification of Stripping Phase. *Talanta* 2015, 144, 356–362. [CrossRef]
- Lokhande, R.V.; Bhagure, G.R.; Dherai, A.J.; Naik, P.R.; Udani, V.P.; Desai, N.A.; Ashavaid, T.F. Analytical Method Validation for Estimation of Neurotransmitters (Biogenic Monoamines) from Cerebrospinal Fluid Using High Performance Liquid Chromatography. *Ind. J. Clin. Biochem.* 2022, *37*, 85–92. [CrossRef]

**Disclaimer/Publisher's Note:** The statements, opinions and data contained in all publications are solely those of the individual author(s) and contributor(s) and not of MDPI and/or the editor(s). MDPI and/or the editor(s) disclaim responsibility for any injury to people or property resulting from any ideas, methods, instructions or products referred to in the content.